


ETDFL MAT + BELT

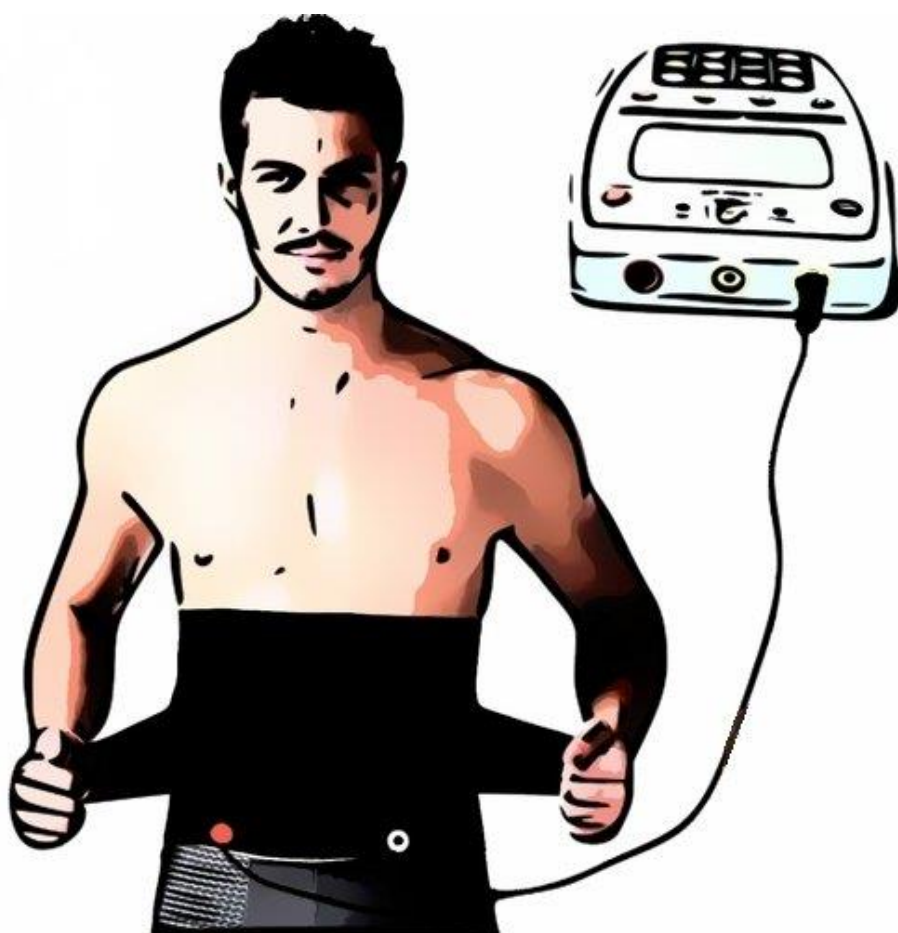

INSTRUCTION BOOK

English

© August / 2017 - © 2018

Index	page 2
Introduction	page 3
. Introduction to Electrotherapy	page 4
Technical Description: ETDFL MAT & BELT	Page 5
. HOOKING UP ETDFL MAT:	page 6
. ETDFL MAT: Bed Mat Method:	page 7
. ETDFL MAT: Body Wrap Method, & Extension	page 8
. Body Wrap, Best Voltage Settings	page 9
.	
ETDFL BELT: Wrap Methods	page 10
. ETDFL Head-Belt Forehead Method	page 10
. ETDFL CHIN / MOUTH / UPPER NECK Method	page 11
. ETDFL FULL NECK Belt Method	page 12
. ETDFL Leg-Belt Wrap Method	page 13
. ETDFL Arm-Belt Wrap Method	page 14
FAQ	Page 15-16
Warranty & Disclaimer	Page 17
ABOUT ETDFL	Page 18
World Sales and Service Centers	Page 19

Introduction:

Congratulations on your purchase of your ETDFL Series Healing MAT or ETDFL Healing Belt, also known as Dr Rife Healing Mat, and Dr Rife Healing Belt.

Since 2016, the ETDFL MAT and ETDFL Belt have been used by thousands of people with Electro Therapy Devices, to transmit the Healing Frequencies to specific parts of the body.

The New Design 2017-2018 ETDFL Mats use advanced 3M conductive rubber technologies built into a body wrap mat, which can now allow your mat to be affixed to any area of the body. The 3M Rubber is a safe and hygienic property, and can be easily cleaned for multiple users.

ETDFL Recommend you only use the Rife Digital Professional V2 (19v) * and Zapper Digital Professional V2 (19v) * with the ETDFL Mat/Belt

However, The ETDFL Mat & Belt are also compatible with:

Rife Digital Professional Series 1 (16v)
Rife Digital Ultimate Series 1 (16v)
Rife Digital Ultimate TRIO (16v)
Zapper Digital Professional Series 1 (16v)
Zapper Digital Ultra Series 1 (16v)
Zapper Digital Ultra TRIO (16v)
Dr Bob Beck Professional Machine

*** Note: The ETDFL Mat and ETDFL Belt are designed for direct use with**

*** Rife Digital Professional V2 (19v)**
*** Zapper Digital Professional V2 (19v)**

The V2 Series have voltage adjustment which allows you to use the ETDFL Mat/Belt directly against the skin/body. Direct contact provides the most efficient use of the ETDFL Mat technology.

Introduction to Electrotherapy with ETDFL Mat, and Belt

The following guide provides an introduction to Electrotherapy and recommendations for the use of your Mat or Belt with your Electro Therapy Device (ETD).

CAUTION: Read this brochure carefully and use your Mat or Belt exactly as directed!

Basics:

What is electrotherapy Mat / Belt treatment with the ETD for?

- Use the ETDFL Healing Mat on any larger area of the body (e.g. Stomach, Chest, Back to focus the electrotherapy at those precise areas for more effective transmission of frequencies
- Use the thinner, compact ETDFL Healing Belt around the Head, Neck, Arm or Leg, for treatment of smaller areas.
- Use electrotherapy to fight certain viruses, bacteria, parasites, etc. The electrotherapy works in a similar way as Homeopathy: using natural principles to heal the body. Everything around us already has an electrical charge and thus a certain electrical frequency - we are always and everywhere exposed to it. If you have a specific virus, it lives in a given frequency field. Let us assume for simplicity it lived only at one single electrical frequency. When you feed the same frequency generated by an ETD to your body (with very particular waves), the virus gets over-stimulated and dies.
- Use electrotherapy to restore vitality and health: The programs are producing a harmonic resonance in the body. Each symptom brings electric dissonance in the organism. With the more general programs such as back pain, detoxification programs for lung, liver, etc., the dissonance in the body will be transformed into harmonic resonance. Such programs are particularly useful for meridian treatments connected with 4 strap-combinations.
- Use electrotherapy for Detoxification. Detoxify Heavy Metals, Chemicals and other unwanted toxins in the organs
- Use electrotherapy for the treatment of all disease Conditions

Technical Description: ETDFL MAT & BELT

Frequency range

Tested 0:01 Hz - 50 MHz

Frequency Wave and Sideband

Transmits all ETDFL frequencies with +99%† conductivity to the body on the RDP-V2* and ZDP-V2* machines. All Dr Rife MOR Carrier Wave and RF frequencies are reproduced through the Positive and Grounded areas of the Mat and Belt with +99%† efficiency.

Voltage input

Set your RFP-V2* / ZDP-V2* onto Voltage Level 12.5v DC. Only at these lower voltages can you use the mat, or the belt directly against the skin. If the frequency output is too intense, then you can use a t-shirt between the Mat and your body. Should you feel to increase voltage, do this one step at a time, move to 14.5v and then 17.5v. If any sensitivity is noticed, then switch voltage down.

- Note: All other models (16v) will require a T-shirt between the mat and body as there are no voltage adjustments on First Series Rife Machine, and TRIO series

Clinic Use, and Cleaning

If using the Mat or Belt for multiple people, or for a clinic, wipe down the rubber area after use with a mild detergent. Do not use alcohol directly on the rubber. You cannot wash the mat in a washing machine, as this will affect the 3M Conductive Rubber surface, You can wash by hand in warm water with a little laundry liquid, however it is suggested to only use a mild detergent "wipe down" method directly on the dry mat.

**** (RDP-V2) Rife Digital Professional V2 2018 Series.***

**** (ZDP-V2) Zapper Digital Professional V2 2018 Series.***

† 99% Conductivity Relative to Copper/Gold at 100%


HOOKING UP ETDFL MAT:

Locate your "Leads" in your Rife Machine Bag, They are white in color with 2 Buttons, Red and Black



Plug your power into the Red socket of machine, and RUN/SWEEP your Program. Now Plug the Leads into One of the Yellow Sockets. If using 2 mats then insert 2 sets of leads.

ETDFL MAT: Bed Mat Method:

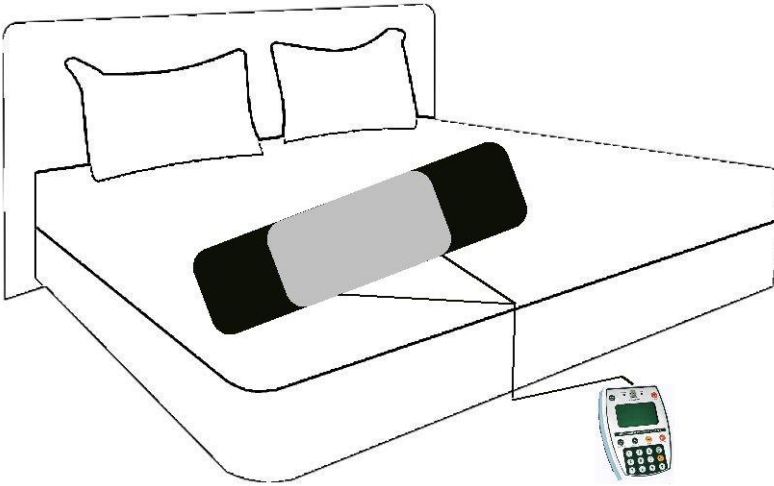

The Bed Mat Method is the traditional method for overnight treatments with the Healing Mat and the RDP.V2

- Plug the Red and Black Leads into the ETDFL Mat. You will notice that there is a RED button and a BLACK button on the top of the mat, indicating where you connect each of the leads.
- Lay Mat on the bed. Make sure the Connection Point is UNDER the mat, and the top side only shows the Rubber with the Red and Black indicators (no wires on top of mat)
- Set your Machine to 12.5v and see if you can lay on the mat directly. If there is too much sensitivity from the frequencies then place a cotton sheet over the mat, and lay on a sheet.

If Sensitivity is not an issue, then switch the voltage up slowly from 12.5v to 14.5v, and finally 17.5v, providing there are no sensitivity issues

ETDFL MAT: Body Wrap Method, & Extension

The New Series ETDFL Mat now have Velcro added to the sides, which permits a "Body Wrap". The Body Wrap Method, lets you put the mat directly at the point of the body where the frequencies can be focused most effectively.

Use the ETDFL MAT, Body wrap method, for treatment of Breasts, Chest, Stomach, Groin area, Upper Back, Lower Back, Buttocks. Basically any larger surface area, where the mat can be attached around you.



An Extension Belt is provided with your ETDFL Mat, to suit XL, XXL, 3XL sizes. Just wrap the mat around your stomach, or back, the rubber side facing the body. If you perceive that the Velcro on the existing ETDFL mat does not "meet", then add the Extension Belt to the end of the mat, and wrap this over the existing MAT. This method will allow all Extra Large Sizes to be able to use the Body Wrap method



- **Make Sure the Mat is a comfortable fit, not too tight.**
 - **Make Sure the Connection and wires are on the outside**
 - **Ensure the Rubber mat side is facing the body**
 - **Secure the outer Velcro wings down for a good fit**
- Set your Machine to 12.5v and see if you can use the mat against the skin directly. If there is too much sensitivity from the frequencies then wear a thin cotton T-Shirt, and put the Mat over the T-Shirt. Adjust Voltage up (14v, 17.5v) until a comfortable level of frequencies can be received, without sensitivity issues.

ETDFL BELT: Wrap Methods

The New Series ETDFL Belt now have extended Velcro added to the sides, which allows greater versatility and placement possibilities. The Belt Wrap Method, lets you put the mat directly at the point of the body where the frequencies can be focused most effectively.




Head-Belt Forehead Method

- **Make Sure the Belt is a comfortable fit, not too tight.**
- **Make Sure the Connection and wires are on the outside**
- **Ensure the Rubber mat side is facing the Head**

ETDFL CHIN / MOUTH / UPPER NECK Belt Method

- **Wrap Belt over head with Rubber under Chin**
- **Make Sure the Connection and wires are on the outside**
- **Ensure the Rubber mat side is facing upper neck / chin**

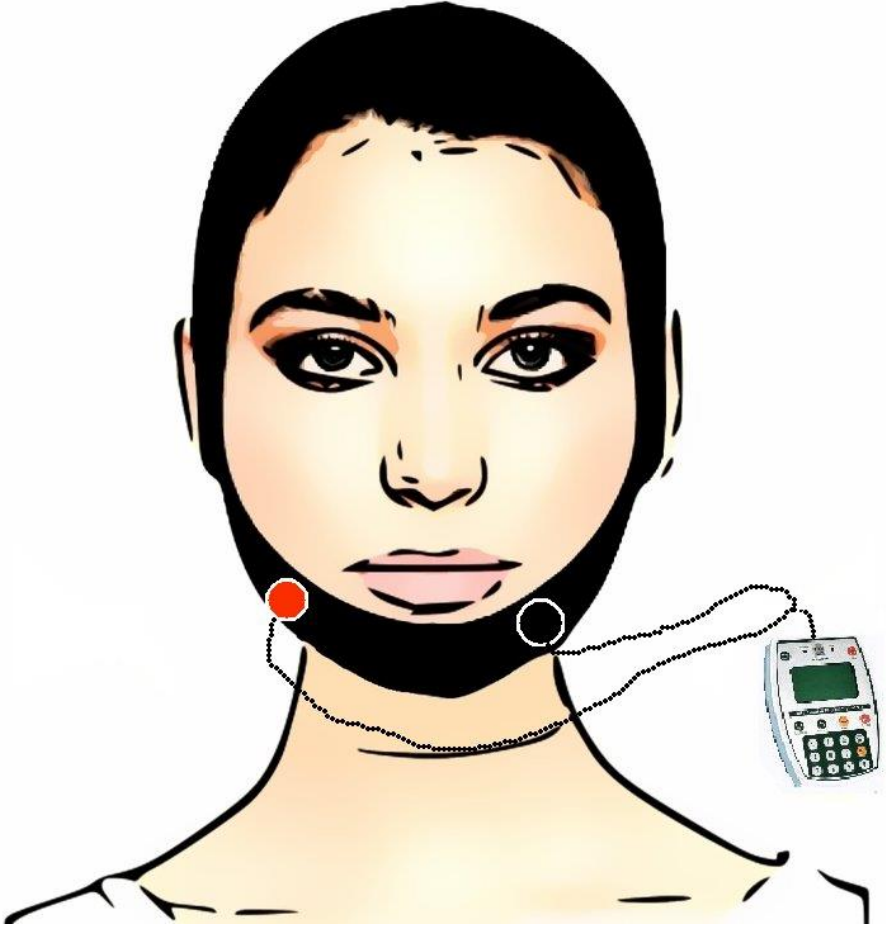

Effective Belt Positioning which allows multiple point treatment, Upper throat, Under Chin, Mouth, Jaw, Face Treatment points.

ETDFL FULL NECK Belt Method

- **Wrap Belt around neck with Rubber facing neck**
- **Make Sure the Connection and wires are on the outside**
- **Ensure the Rubber mat side is facing neck**




Effective Belt Positioning around the neck which allows treatment of the entire neck.

- **Treatment of front of neck, or back of neck**
- **Effective treatment of Throat Bacteria**
- **Treatment for Neck Injury**
- **Treatment of Neck diseases or disorders**

ETDFL Leg-Belt Wrap Method

- **Wrap Belt around Leg with Rubber facing Leg**
- **Make Sure the Connection and wires are on the outside**
- **Ensure the Rubber mat side is facing Leg**

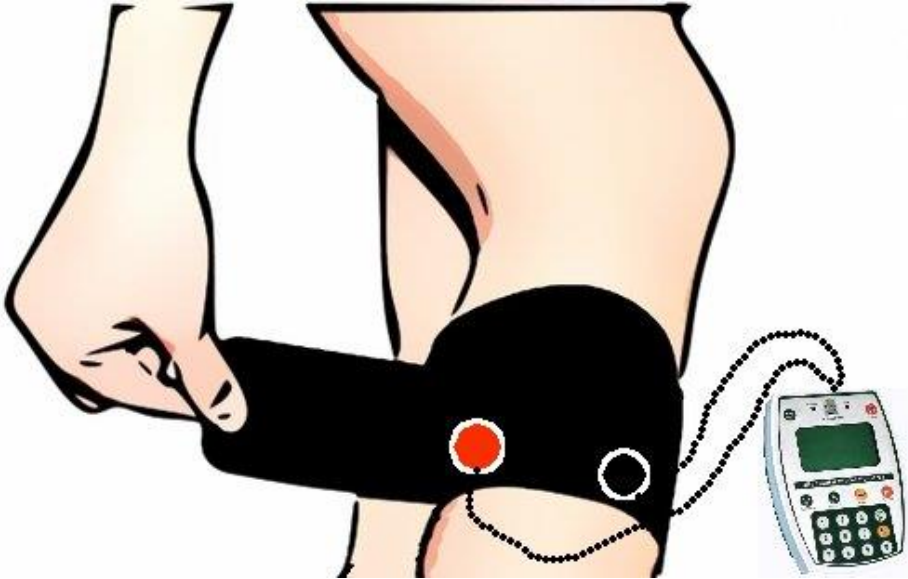


Effective Belt Positioning around the Leg which allows treatment of any part of the Leg, Upper or Lower

- **Treatment of Upper and Lower Leg**
- **Treatment for Leg Injury**
- **Treatment of Leg diseases or disorders**

ETDFL Arm-Belt Wrap Method

- **Wrap Belt around Arm with Rubber facing Arm**
- **Make Sure the Connection and wires are on the outside**
- **Ensure the Rubber mat side is facing the Arm**

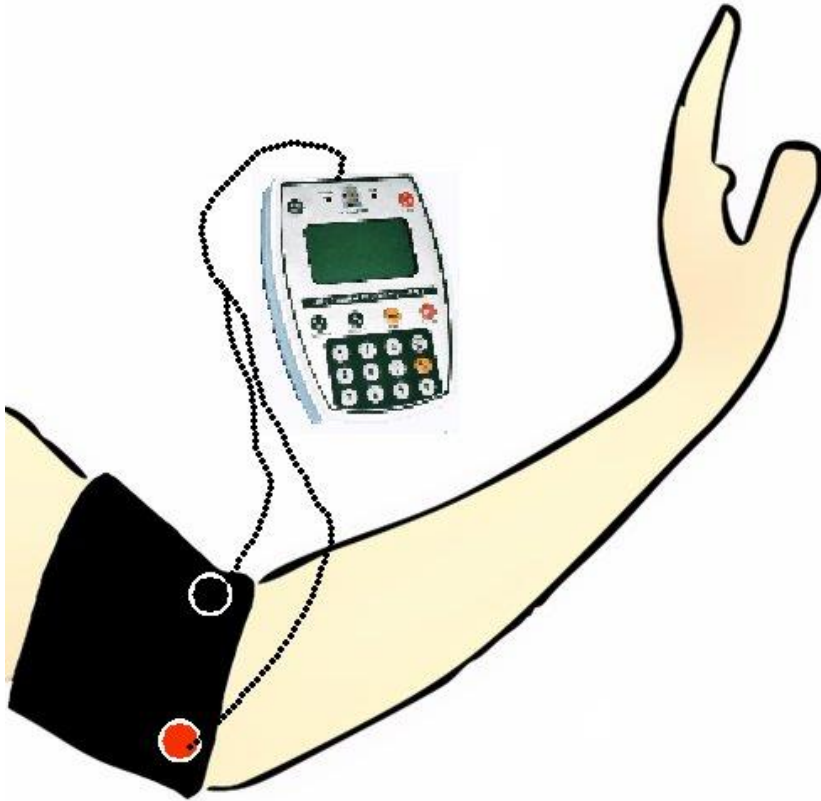

Effective Belt Positioning around the Arm which allows treatment of any part of the Arm, Upper or wrist. You may have to wrap the belt over twice due to smaller arm circumference.

- **Treatment of Upper and Middle arm, or wrist**
- **Treatment for all Arm Injuries**
- **Treatment of Arm diseases or disorders**

FAQ (Frequently Asked Questions)

What setting on the Rife Professional V2 do I use with the Mat or Belt?

- We suggest you look up the exact frequency of the disorder for the treatment. Often you will find the disorder listed within the Presets of the machine. if not use the ETDFL book (Free Download at ETDFL.COM) to look up the frequencies for the disease to be treated.

I have a GB-4000 can I use the Mat and Belt on this machine?

The GB-4000 uses old style banana plug fittings, you will have to find someone to modify the wire system to make it work.

Can I use my Mat or Belt with the RFP.V2 Foot Straps?

Yes you can combine a Mat or Belt treatment while also using the foot straps supplied with RDP.V2, use the second Yellow socket on RDP.V2 for the second set of straps

Where can I buy the most up to date Rife Machine to use my Mat/Belt accessory?

You can purchase the Rife Professional V2 directly from ETDFL.com via shopping cart.

I use the Mat against my skin, and I felt a prickling sensation, and high sensitivity:

Make sure your RDP.V2 machine is set to 12.5v if using directly against the skin. If 12.5v is still too sensitive, then wear a light thin T-shirt. Once the sensitivity issue is resolved, you can crank up the voltage a little at a time (14.5v, 17.5v) permitting the sensitivity issue does not return. The treatments are most effective at 17.5v setting.

I tried to wrap the belt around my arm but it is too long, Help!

Ok, the arm has a smaller circumference, than say the leg, so you can still use the belt effectively on the arm, but you must wrap the Velcro around twice, overlapping the Velcro until you have a good fit. Just make sure that the middle area where the rubber is split by material have both sides of the rubber facing the arm (contacting Positive & Ground)

I have a serious mouth/ tongue "Sarcoma" type issue, how would I use the belt?

Ok, Wrap Belt over head with Rubber under Chin, - Make Sure the Connection and wires are on the outside, Ensure the Rubber mat side is facing upper neck / chin, This will treat the whole mouth, tongue and part of the neck area. you can also opt for the neck belt method also. See page 11 and 12

Can I use the Belt on my 1 year old child.

Yes the Belt may be more suitable than the larger Mat. With small children, always use a t-shirt, or cotton top before fitting Belt. Fit the belt over the T-shirt. Small children must be under constant supervision when wearing Belt accessories, do not fit belt around neck for safety reasons, caution must be exercised when fitting belt around head. Do small treatment times, and keep checking there are no sensitivity issues or rash. Only use 12.5v setting.

Can I use these accessories for my Dog or Cat, Animals?

Yes the Mat can be placed in the bedding of a Dog or Cat for treatment, the belt can be used as a Dog/Car body wrap. You can use the larger mat on horses, but will need to fit the extension.

WARRANTY:

DOA Warranty: If your Mat or Belt arrives with any damage caused by transport/ courier, it will be replaced immediately

3 Month Replacement Warranty: If your Mat or Belt have broken with normal use, we will replace the item within 90 days. Cost of shipping will be charged for the return of a new item to your address.

Due to sanitary reasons a Mat or Belt will not be repaired by our service center. We will always issue a new product, under warranty.

Disclaimer

The E.T.D.F.L. (ETDFL.com) Labs holds no guarantee, that the ETDFL Mat, Or ETDFL Belt will cure any specific diseases. The Belt and Mat technologies have been tested for over 3 years in EU by researchers, who have had successful results using the MAT in conjunction with the Rife Digital Series only. The Mat has not been tested with any other brand Rife Machine, so there are no guarantees should you alter the equipment for other uses. The Mat itself is a Rubber Conductive Conduit, and in itself does not perform any process without a machine attached.

You will be shipped this same Instruction book whether you purchase ETDFL Healing Mat or the ETDFL Healing Belt

ETDFL Healing Mat, and ETDFL Healing Belt are two unique products which can be purchased separately. This book provides instructions for both products to build a greater awareness of our new product: "ETDFL Healing Belt".

ABOUT ETDFL

E.T.D.F.L is an ongoing project by 12 Clinics across the planet who use the Quantum SCIO Bioresonance machine to record Disease frequencies of actual clients, and then produce the 10 most accurate frequencies which are deemed suitable for treatment with the Royal Rife Machine. This project began in 2006 when the original CAFL frequency list was discontinued. ETDFL became the new list, now used by more than 100,000 Royal Rife Machine users across the planet. The accuracy of the frequencies is seen in the positive results of the clients who use this Frequency List in clinics and homes in every country.

ETDFL Berlin Germany, develop the latest accessories and equipment for Electro Frequency Therapy. The ETDFL Mat and ETDFL Belt are a unique product developed by ETDFL Labs Germany.

ETDFL is updated once a year. Latest PDF versions available at: ETDFL.com

PDF versions of ETDFL at **ETDFL.com** allow you to search the ETDFL Database for Disease frequencies which may be more up to date than the List provided within this book.

The Digital Professional V2 Device is licensed to use all ETDFL Royal Rife Frequency Sets © 2017 - 2018 ETDFL Berlin Germany

The ETDFL Mat and ETDFL Belt Devices are licensed to use all ETDFL Royal Rife Frequency Sets © 2017 - 2018 ETDFL Berlin Germany

INTERNATIONAL SALES AND INQUIRIES

USA / CANADA / AUSTRALIA / World
Sales and Service

(web) www.rifedigital.com

(E) healthproducts2@gmail.com

USA ☎ 206-708-6408

International ☎ 0011-1-206-708-6408

United Kingdom & European Union,
Eastern Europe & World Wide

London UK Store, Manager: Rosalie

(web) www.rifedigital.net (E): rifedigitaluk@gmail.com

(web) www.drclark.co.uk (E): clarkzappersuk@gmail.com

UK ☎ 07376056374

International ☎ 0044 7376056374

JAPAN: この店舗の日本人顧客のみを対象とした大規模割引
ウェブサイト rifedigital.jp

マネージャー ピーター・ウィリアムズ博士

Eメール rifedigital@gmail.com



© 2017-2018, *Electro Therapy Devices Berlin Germany*

© 2017-2018, *ETDFL.com, Berlin Germany*