


THE DR BOB BECK PROTOCOL

With The Dr Bob Beck Professional Machine

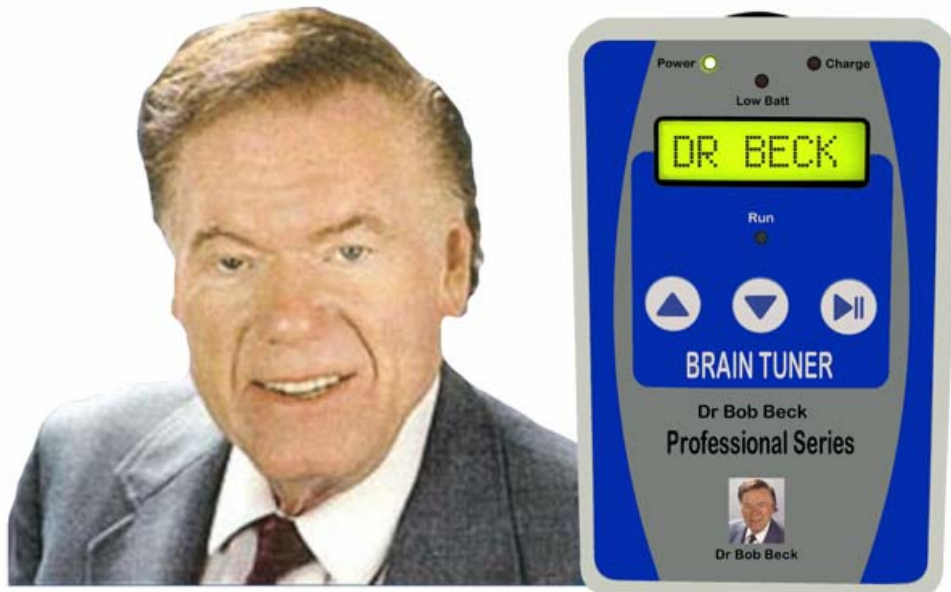

Over 22 Original Dr Bob Beck Frequencies including BT4, BT5, BT6, BT7, Brain and Disease Treatments and The Colloidal Silver Program, in a Portable Hand Held Machine.

© 2017 . The Bob Beck Organization LLC
www.drbeck.net

1. Index	page 2-3
2. Introduction	page 4
3. Technical information	page 5-6
4. Basics	page 7-8
5. Plugging in, Setting Up, Charging	page 9-10
6. Using the Machine Touchpad	page 11-12
7. Plugging in the Accessories	page 13-18
8. The Colloidal Silver Apparatus	page 19-20
9. The Conductive Breast Pads	page 21-23
10. Earclips for Brain Tuner Series (BT)	page 24-25
11. The 22 Bob Beck Professional Programs	page 26-27
-Program 1: Dr Beck Brain Tuner BT-4	page 28
-Program 2: Dr Beck Brain Tuner BT-2	page 29
-Program 3: Dr Beck Brain Tuner BT-6	page 30
-Program 4: Dr Beck Brain Tuner BT-7	page 31
-Program 5: Schumann Earth Frequency Program	page 32
-Program 6: Alpha Wave Program	page 33
-Program 7: Beta Wave Program	page 34
-Program 8: Gamma Wave Program	page 35
-Program 9: Delta Wave Program	page 36
-Program 10: Theta Wave Program	page 37
-Program 11: Lymph Treatment	page 38
-Program 12: Blood Treatment	page 39
-Program 13: Blood Pulser	page 40
-Program 14: Addictions Treatment	page 41
-Program 15: Bacterial Infections	page 42
-Program 16: Virus	page 43
-Program 17: HIV/ Aids	page 44
-Program 18: Herpes	page 45

-Program 19: Cold and Flu	page 46
-Program 20: Colloidal Silver Generator	page 47
-Program 21: Cancer Treatment	page 48
-Program 22: Allergy Treatment	page 49

12. FAQ. Frequently Asked Questions **page 50**

- Finding the Right Frequency for illness	page 50
- How much Silver Can you Drink ?	page 50
- Silver Atomizer: How Long to Use Machine	Page 50
- Program for ADHD Treatment	page 50
- Program for Anxiety and Depression	page 51
- Program for Learning and memory retention	page 51
- Program for General Meditation	page 51
- Deep Meditation, Transcendental Meditation	page 52
- Program for Schizophrenia and Bi-Polar?	page 52
- Charging the Dr Beck Professional Machine	Page 52
- Beck Machine vs Rife Digital Professional	Page 53
- Using Beck Machine at Workplace or School	Page 53
- Using Beck Machine for Chronic Conditions	Page 53
- Dr Beck Professional compared to SOTA	Page 54
- Dr Beck Protocol. Web Link to Shopping Cart	Page 54
- Dr Beck machine Language	Page 55
- Beck Machine in non US countries 110v - 240v	Page 55
- Beck Machine Charging and Using Device?	Page 55
- How Long Do I need to use the Machine?	Page 55
- Is Colloidal Silver a new invention?	Page 55
- How do you use Colloidal Silver?	Page 56
- How Does Colloidal Silver work?	Page 56
- Essential Oils & Silver with the Dr Beck Atomizer	Page 57
- Flex Straps and Velcro Straps on the wrists?	Page 57
- Warranty Period on the Dr Beck Machine	Page 57
- Notes	Page 58
- Contact, Telephone, Sales, Inquiries, Email	Page 59

2. Introduction:

Congratulations on your purchase of the Dr Bob Beck Professional Machine. The following Instruction book provides an introduction to brain treatments, disease treatments and electrotherapy, colloidal silver and guide for using your device properly.

Also provided is the official Brain Tuner Program list, which includes all the "BT" Machines 4, 5, 6, 7 created by Dr Beck, as well as Brain Wave sets, Alpha, Beta, Gamma, Delta, Theta, Schumman Earth Frequencies, Disease treatment programs, and a built in Colloidal Silver program running at Alpha Wave frequency for creating at home, the most bio-active Colloidal Silver available on the market. All features are integral to the Dr Bob Beck Protocol. Never before has any "Beck Machine" included the complete Beck Protocol until now: The Dr Bob Beck Professional machine brings the full 22 Dr Beck Programs on a single hand held portable device.

CAUTION: This machine is an electro-therapy device with a very high frequency output potential. Read this book carefully and use your Dr Bob Beck Professional machine exactly as directed!



3. Technical Information:

Program Type

22 Fixed Internal Frequency Programs

*** Dr Beck has designed all the programs within the machine to an exact standard of the "Dr Beck Protocol", so programs cannot be added or altered.**

Frequency range: 0.01 Hz - 10 kHz

Frequency waveform: Alters from Square to Sine according to the Selected Program

Battery Type: 4 x SONY AA 2000+ mAh High Power Rechargeable

Use only Battery type AA-1.25v-mAh2000 up to mAh3000

The Machine is supplied with 4 AA Batteries in the bag

Battery Charging

Internal, Batteries are charged while inside machine using charger

Charging Time:

7 Hours First Charge (or until charge light is off)

1 Hour Quick Charge (charge light is Red but holding charge)

4 Hours Normal Charge (charge light goes off after 3-4 hours)

Machine ON time while connected to the machine with Leads

48 - 72 Hours non-stop treatments while connected to machine

*** Try to Recharge Fully to ensure longer treatments**

*** Rechargeable Batteries loose Amp/Voltage over time this is normal**

*** Replace AA Batteries with a good brand (Sony/Kodak) with mAh +2000**

*** Your reseller will likely to be able to supply the correct type replacements**

Voltage Output: 5v DC at 2000+ mAh

Voltage input using Charger: 9v DC (110 V / 250 V SMPS) All world voltages accepted, external DC switching-mode power supply adapter. * Charger is only for recharging the device.

*** Do not use the machine while recharging as the machine produces exact frequency waves on a DC Cycle, This machine is not to be connected to a person during Recharge Cycle on the AC Mains Charger.**

CPU Operating Temperature (Safe) *

26° - 33° Celcius

78° - 87° Fahreneit (Normal Operating Temperature)

RCA Output socket

- Use With Leads to Connect all Dr Beck Accessories

Cleaning the Bob Beck Machine

Clean Machine with a soft cloth if required but do not allow water to come near the sockets. Do not clean with Alcohol or Abrasives

Accessories

The Bob Beck Machine is shipped with many different types of accessories, which plug into the standard "RCA" Plug Leads. Some Accessories such as GEL PADS, Flexible Straps, Colloidal Silver Cup Clip, are standard in the package. Other accessories can be purchased separately as they are specific for women, or specialized other purposes. (e.g. the conductive breast treatment pads, Dr Beck Magnetic Pulser, or The Dr Beck Atomizer

4. Basics:

What does Electrotherapy Treatment with the Dr Beck Machine achieve for Disease Conditions?

Use The Bob Beck Professional machine to fight certain viruses, bacteria, parasites, etc. The Beck Machine works in a similar way as Homeopathy: using natural principles to heal the body. Everything around us already has an electrical charge and thus a certain electrical frequency - we are always and everywhere exposed to it. If you have a specific virus, it lives in a given frequency field. Let us assume for simplicity it lived only at one single electrical frequency. When you feed the same frequency generated by an The Bob Beck machine to your body (with very particular waves), the virus gets over-stimulated and dies.

Use the Dr Beck Machine to restore vitality and health: The programs are producing a harmonic resonance in the body. Each symptom brings electric dissonance in the organism. With the more general programs such as back pain, detoxification programs for lung, liver, etc., the dissonance in the body will be transformed into harmonic resonance. Such programs are particularly useful for meridian treatments connected with 4 strap-combinations.

Specific Brain Treatments are achievable using the BT Program Series which have undergone years of research by Dr Bob Beck for treatment of many Brain Related conditions. Originally Dr Bob Beck released a different machine with each frequency set, the first machine which reached the public was the BT4, followed by BT5, BT6 and BT7.


As it was not practical to buy 4 different machines to receive these amazing treatment programs, we have incorporated all 4 BT Programs on this single machine, along with all known Brain Wave type programs, and the Shumman Earth Frequency

Specific Frequencies for Disease treatments are listed later in the book and follow the Dr Beck research into reversal of disease conditions using specific frequencies and microcurrents.

Remember that Bob Beck used unique methods of delivering those frequencies into the body which included the Pulser (a Program on the machine) which incorporates the use of conductive Gel Pads on the inside of the leg (or wrist) to pulse and electrify the blood which passed through the arteries and veins in this area. As the current is low this Pulse is not felt by the user, so do not expect any shock or zap at the location of the Gel Pads.

Use the Colloidal Silver Program as the second most important contribution to your therapy. Use the Colloidal Silver in conjunction with all Dr Beck Protocols and programs. The Colloidal Silver Apparatus is included with your Bob Beck Professional Machine, and is in the Bag.

You can also Purchase a separate device, The Dr Beck Atomizer (accessory) for Lung Specific treatments with the Silver you make on the Beck Machine.

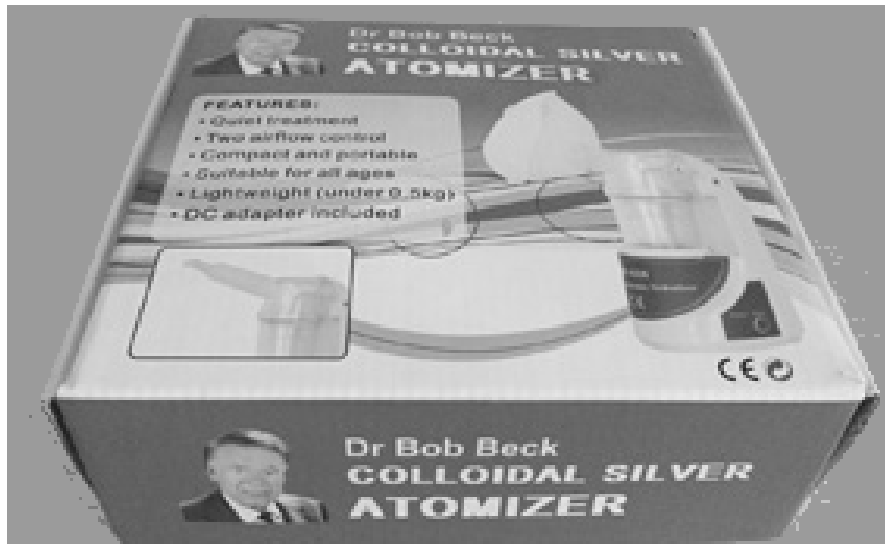

PLUGGING IN, SETTING UP, CHARGING

Insert the 4 Batteries into the Back Compartment of the Dr Beck Professional



Only use Rechargeable 1.25v AA batteries with a mAh Amperage rating of 2000 mAh - 3000 mAh. Make sure the Batteries are the correct way around or the machine will not work, or you could create a problem with the battery or machine. Generally speaking the FLAT (-) side of the AA battery should mount on top of the spring in the battery compartment. Once batteries are installed, place cover back over battery compartment. Make sure the Blue ON/OFF at top of machine is at OFF position.

BUY AN EXTRA SET OF BATTERIES IF YOU ARE A HEAVY USER

If you are a "Heavy" user of the machine, and require 7 day extensive treatments with little "down time" then purchase 4 extra AA 1.25v batteries at +2000 mAH and perhaps "4 Pack AA Battery Charger", so you can quickly dismount the internal batteries and load in another pack

CHARGING THE DR BECK PROFESSIONAL

Plug in Charger for first time and allow the machine to charge for 3-4 hours. usually the Red charging light will turn off when batteries are fully charged.

When Charge light turns off indicating a fully charged battery, then turn on the machine by sliding the Blue ON/OFF button to the left. Button is located at the top of the machine. Left is the ON Position




Using the Machine Touchpad

Selecting a Program:

In order to select a Program on the Dr Bob Beck Professional Machine, you must ensure you have fully charged batteries installed in the machine

(1) Slide the ON/OFF Button to the LEFT (see picture) The Power Light will then come on

- Machine LED will Light up with [Dr Beck] on the Screen, and then move to Program 1 (of 22) which is BT-4

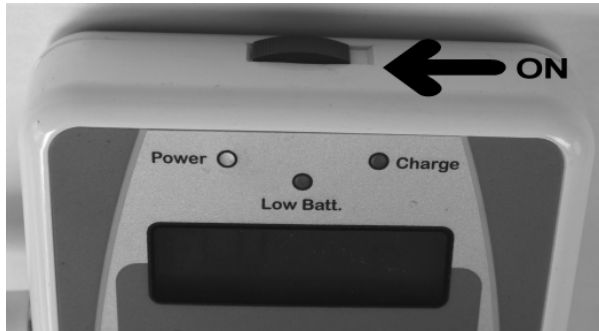

(2) The First Program shown is BT-4 in "i" IDLE mode (not operating)



- You may Press the UP (▲) or DOWN (▼) arrow to make another selection from the 22 Programs installed.



(3) Press (▶||) Play / Pause Button to Start a Program or Pause the program




- Once the ▶|| Play Button has been pressed, then an "r" will appear on the right side of the screen, indicating the program series for BT-4 in this example above is operating (running)

- Note: On the first series Dr Bob Beck Professional models, the RUN Bulb was de-activated to reduce battery consumption. So In order to check if your machine is in RUN mode, please look at the LCD screen for the letter "r" (Run Mode)


Plugging in the Accessories

(1) Please Locate the RCA Leads in your Dr Beck Bag

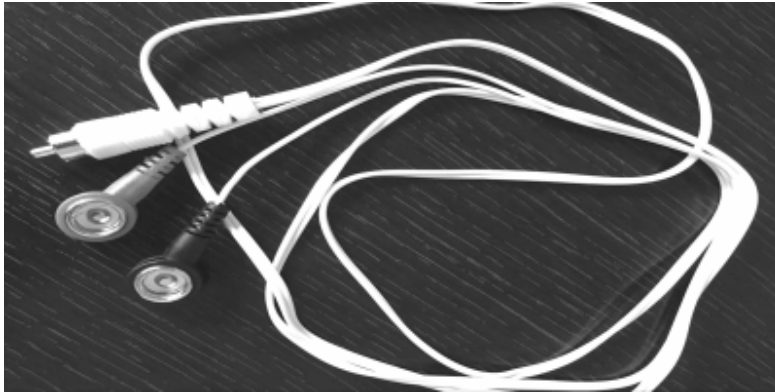


(2) Plug the White RCA Plug into the RCA socket on the left side of the Dr Bob Beck Professional Machine

(3) For General Body treatment: Clip the two Flexible White Straps into the Button Clip on the RCA Leads, and place over the hands or feet

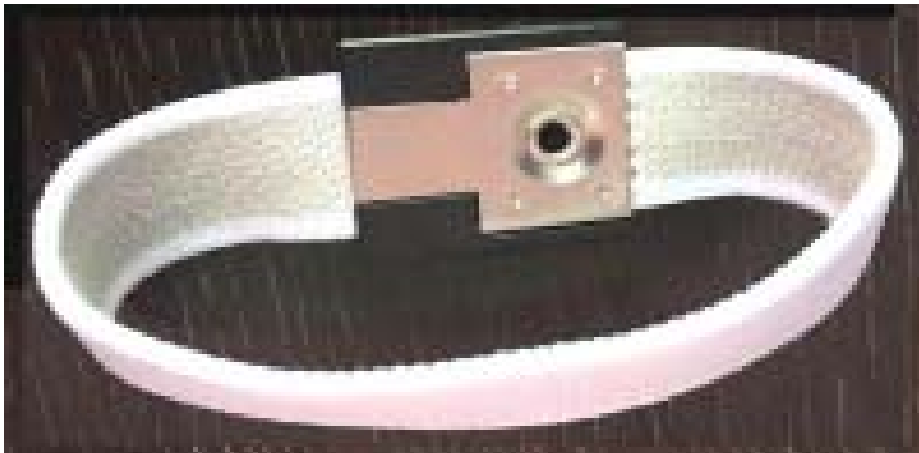

VELCRO Straps

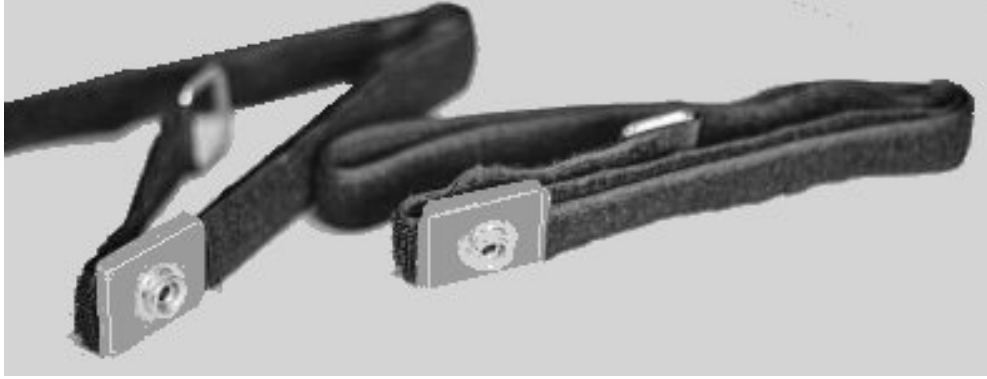

Example: Straps attached to soles of feet



Example: Straps attached to Palms of Hands




(4) For Blood Pulsing, Blood Cleansing, Lymph Treatment, or specific body area treatments; Connect the GEL Pads into the RCA Leads with PINS




Gel Pads Can Be used anywhere on the body safely




Gel Pads Can be re-used if you place them back on the plastic sheet after use.

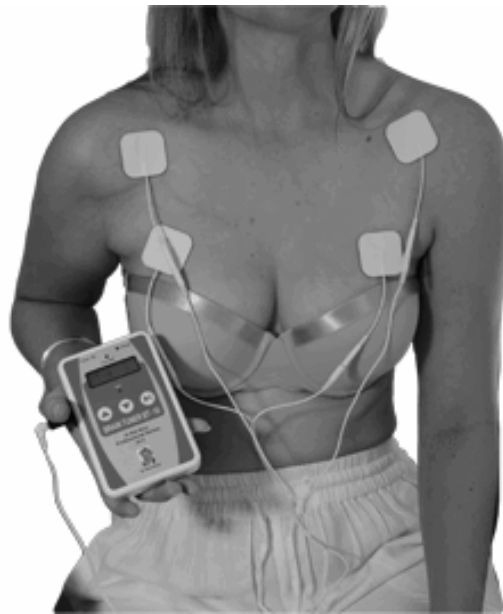


Gel Pads are the most versatile accessory: See above: A Comprehensive LYMPH Program using 4 Gel Pads. Gel pads are perfect for the Dr Beck Blood Cleanser and Blood Pulsar, by locating Gel Pads on the inner arterial position of the foot

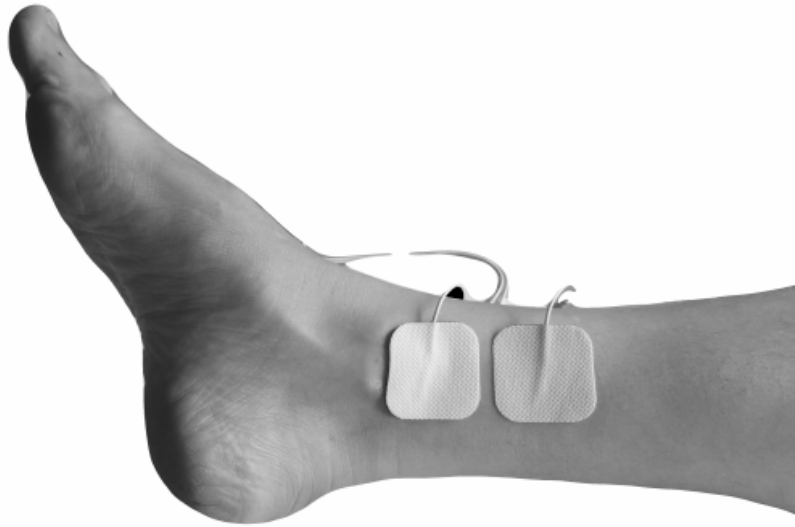

Gel pads are perfect for the Dr Beck Blood Cleanser and Blood Pulser, by locating Gel Pads on the wrists across the main arteries

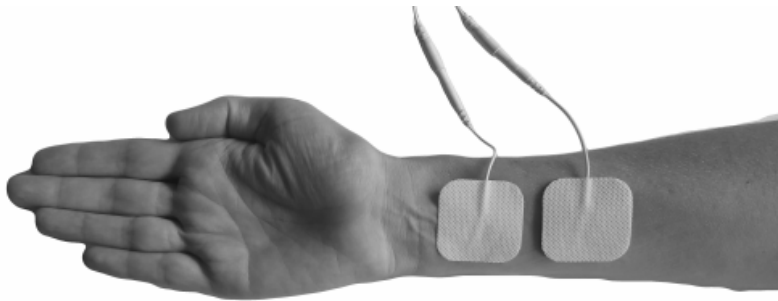

Gel Pads can also be used for any of the Brain Treatment Programs by locating Gel Pads on the temples.

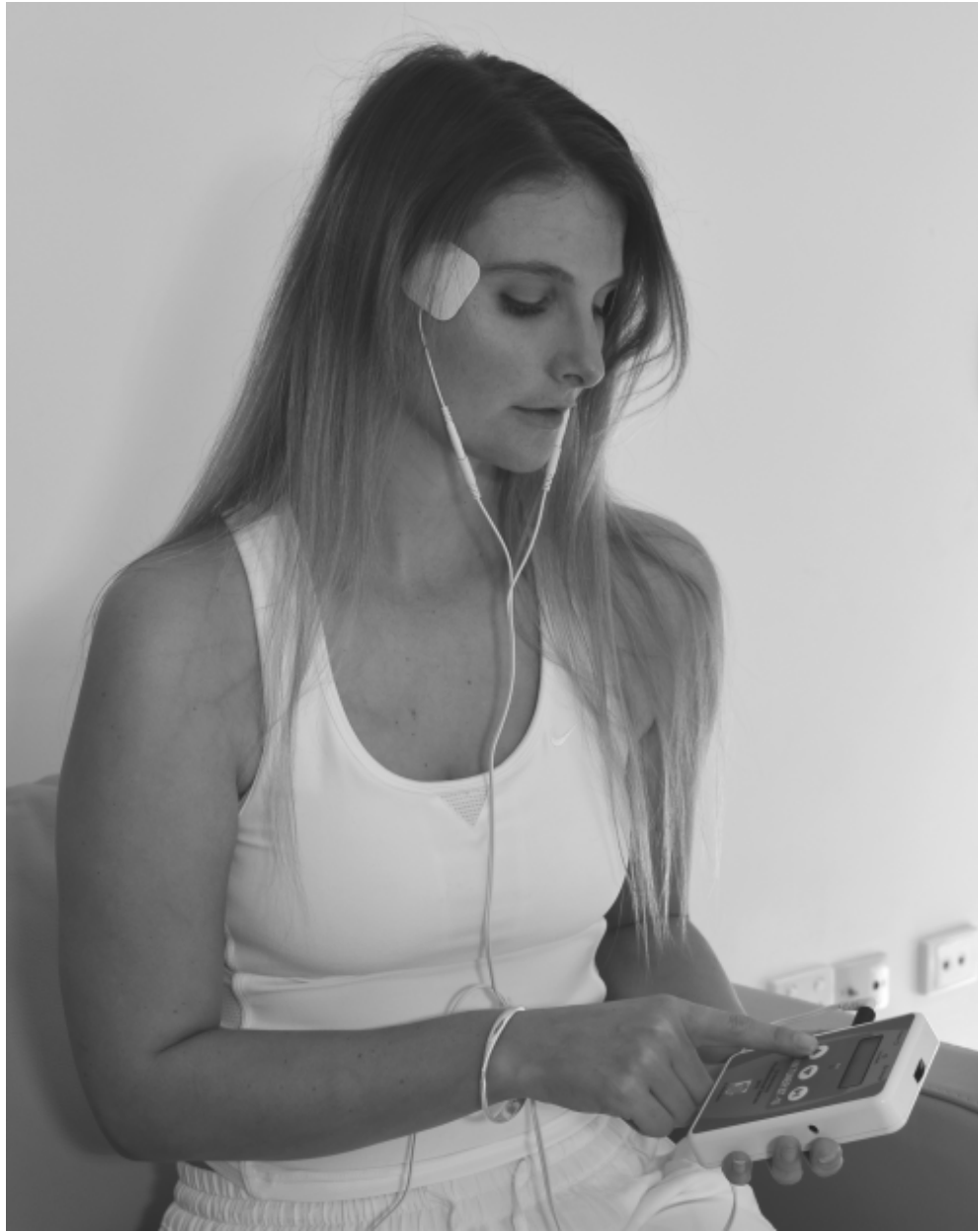

(4) Colloidal Silver Apparatus:

Move the UP (▲) or DOWN (▼) arrows until you reach the [SILVER] setting.


Plug in the Cup-Clip with Silver Rods installed, place the Rods in a cup of water, and the other cable into the RCA Socket of the Dr Beck Professional machine. Make sure the Rods are not touching each other, and properly immersed in the cup of water.

Press ►|| Play Button and the Program will begin making Colloidal Silver in the ALPHA Wave frequency, which provides the most Bio-Active Silver on the planet. The Program will stop in about 5 minutes, and when the Colloidal Silver has been dissolved throughout the water

Press ►|| Play Button again to double the strength of silver in the water.

Colloidal Silver is an important part of the Dr Beck Protocol and should be used in conjunction with all other Dr Beck Treatment Programs.


Purchase the Dr Beck Colloidal Silver Atomizer machine, if you wish to make further use for the Colloidal Silver you have just made on the Beck machine. This Liquid Silver can be put in the Atomizer machine to assist internally for Lung and Chest Infections, or any Lung related Disease

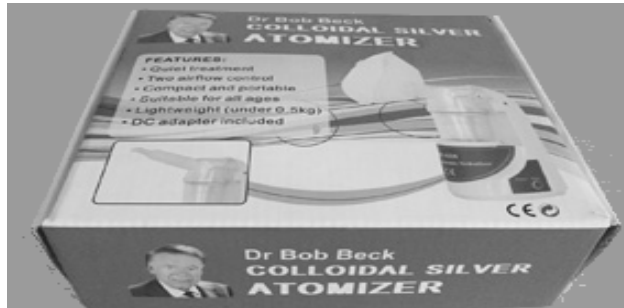



The Conductive Breast Pads


(4) The Conductive Breast Pad may be supplied as an Optional Accessory. The Conductive Breast Pads come in a packet of 2 and come with the RCA Pins Leads

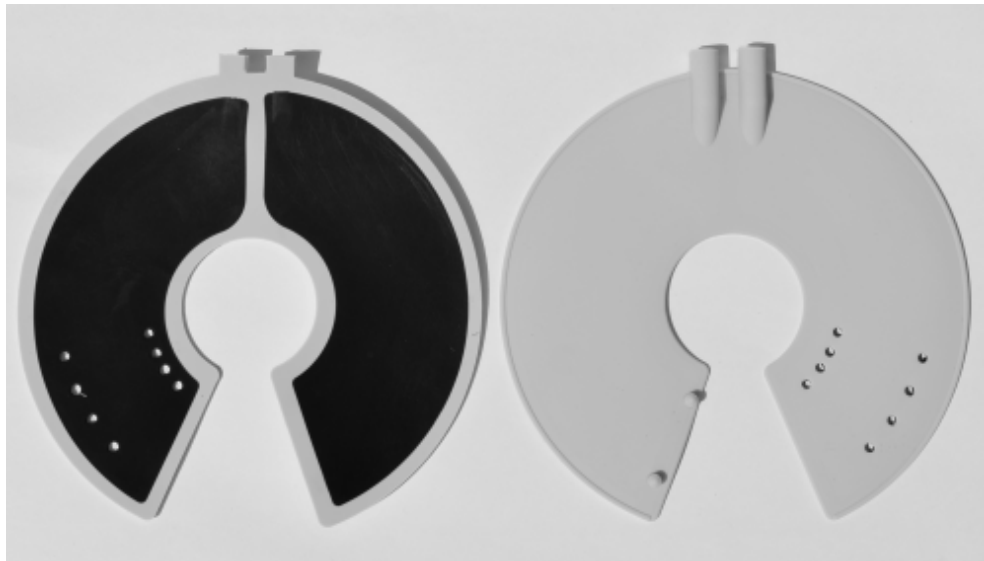

The Breast Pads, for Lymph and Cancer Treatments:

The materials of the Breast Pad item are made from rubber and do not have chemicals or solvents in them.

The Black Side of the Breast Pad is the Conductive Side and should face the skin. The Pink side is the outer rubber and is not conductive. Make sure the Pink Side faces outward and the Black side is against the skin. Move the Rubber plugs until they clip into the opposite holes and clamp down. Make sure you align the holes according to your breast size, and there are no air gaps under the rubber due to incorrect size adjustments.

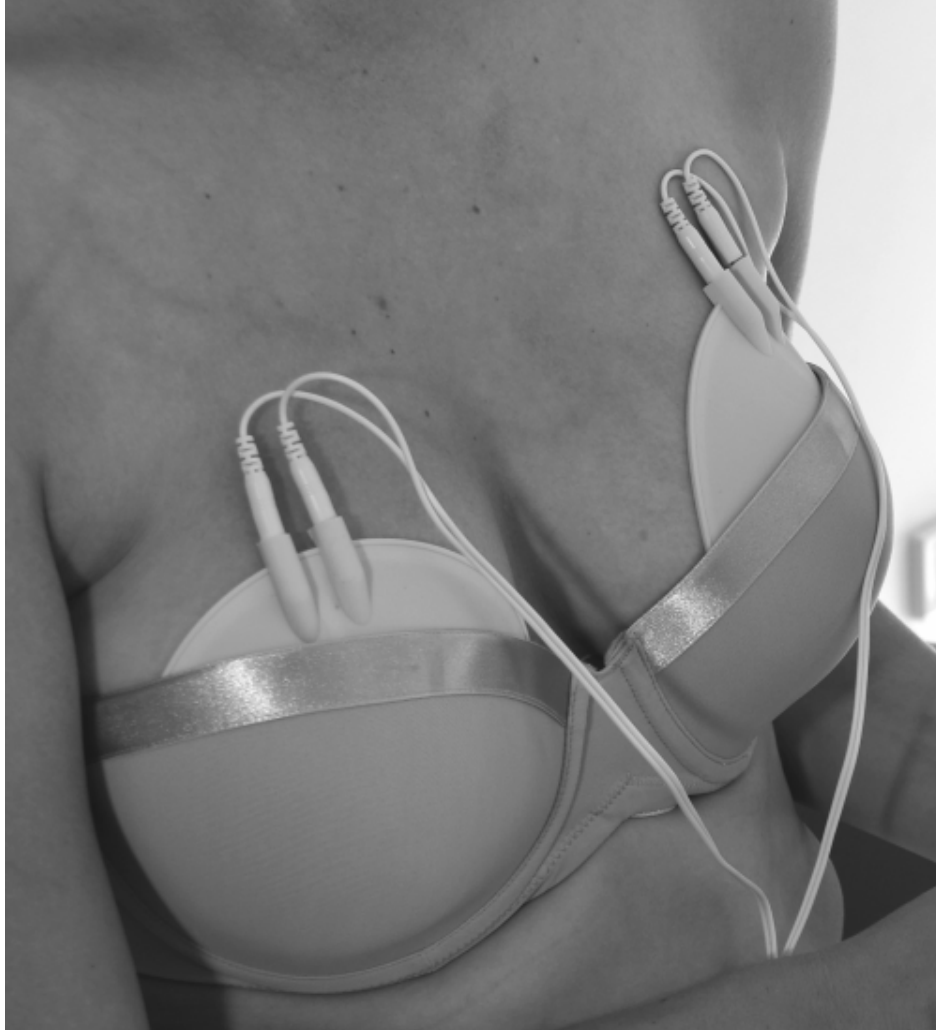


Plug the Breast Pads into the RCA Pin Leads, and then Plug the leads into the Dr Beck Professional machine.

Select the suitable frequency for treatment.

For example for Breast Cancer treatment, use Program 21 [CANCER]

For treatment of Lymph Glands in the areas above and around the breast, use the Program 11 [LYMPH]

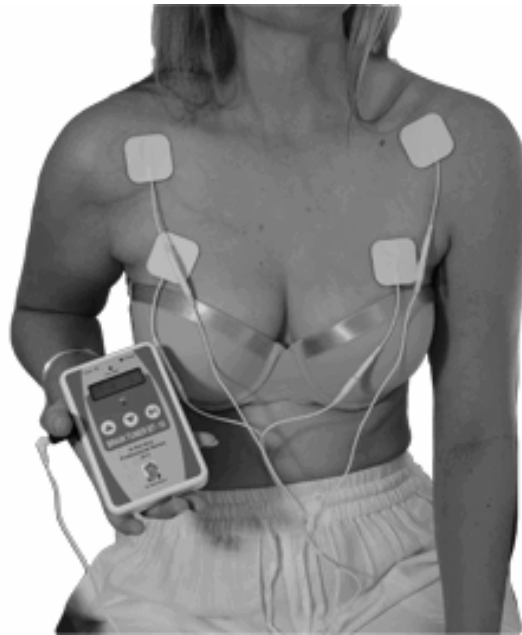


There are many other types of treatment Programs that can be used with the Breast Pads, and it is possible to do "Combination" Treatments by using both Gel Pads and Breast Pads at the same time for greater coverage during treatments.

Treatment times, and other FAQ can be found under the page for Program 21 [CANCER], and Program 11 [LYMPH]

Earclips:

For Dr Beck BT Series, Brain Wave, & Schumman Earth Frequencies

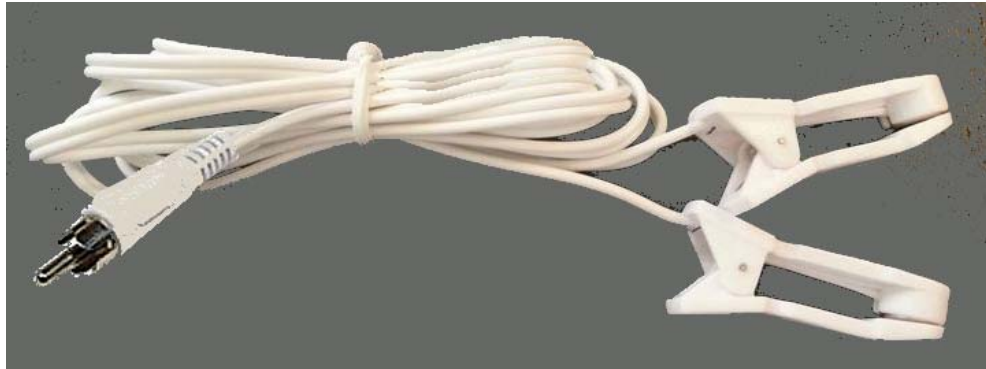

Earclips are specific to Brain Treatments as they permit delivery of frequencies directly to the brain through the earlobes



**Earclips are compatible with the Dr Beck Professional Machine
Program numbers:**

1	BT-4	Dr Becks Original Machine
2	BT-5	Dr Becks Second Machine also known as BT PLUS
3	BT-6	Dr Becks Third Machine
4	BT-7	Dr Becks Fourth Machine
5	SCHUM	Schumman Earth Freq.
6	Alpha	Alpha Wave
7	Beta	Beta Wave
8	Gamma	Gamma Wave
9	Delta	Delta Wave
10	Theta	Theta Wave
11	Lymph	Lymph Treatment
22	ALRGY	Allergy Treatment
19	FLU	Cold and Flu

The Earclips focus the frequency delivery to the brain by locating the frequency transmitting earclip terminals on the earlobes.

The Earclips are supplied in the Bag with a Conductive Soft-Pad covering a metal contact underneath. Your Package will be supplied with extra adhesive Soft-Pads should the pads wear out or be pulled off.

To avoid possible irritation, do not connect the ear-clips to your ear directly on the metal contact. Always stick another adhesive soft-pad to the metal part of the earclip should the pad fall off. Adhesive Replacement Pad Packs are available from your Dr Beck Store, or through the online website.

Please see individual Program listings later in the book to see how the Earclips can benefit you while you use each Brain Focus Program.

Note: Use the Earclips to assist with Some programs such as Cold & Flu where the frequencies need to be concentrated at the head area. Gel Pads can also be used on the neck for Cold/Flu treatments.

THE 22 DR BECK PROFESSIONAL INTERNAL PROGRAMS

There are 22 Dr Beck Professional Programs installed within this machine. The Programs cannot be modified as they are designed exactly according to the Dr Bob Beck Protocol, so you can only RUN the programs.

Unlike other "Beck machines" with 1 or 2 programs, we believed it important to captivate the essence of Dr Becks entire philosophy by incorporating all of Dr Becks programs in one device, the Dr Beck Professional machine.

Unlike other "Beck Machines" we have included the Colloidal Silver maker INSIDE THE MACHINE !

This Beck Machine, is without doubt the most state of the art device to be found on the market, not only is it portable, but uses a vast range of accessories, Gel Pads, Arm Bands, Leg Bands, Ear-Clips, Breast Pads to name a few, in order to deliver these frequencies into the Brain or the body in the most effective manner.

We have selected the most reliable and common battery, four of the AA Rechargeable 2000+ maH to power the Dr Beck Machine, and have customers using this machine non stop, off the Grid, on the move, for 3 days straight while connected

A Perfect machine for Work, School, Home, Driving, Flying. You can easily place the Beck Machine inside a pocket, and connect to the body using Gel Pads completely out of sight. Do your programs discretely while walking or riding a bicycle.

The Following Pages will dedicate ONE Page to each of the 22 Dr Beck Programs and give you more detailed information about each program, and the best accessory to use to complete that program

Remember: Serious Diseases require the longest treatment time, do not expect a "Silver Bullet" result using the machine for an hour or two. Best results are achieved with the longer use, and the overnight treatments which can last 7-9 hours per night.

This extended treatment time can also apply to the Brain Treatment programs for some Brain Conditions. The Body adjusts to the longer treatments and syncs with the frequency waves. The extended treatment time can make a significant effect on reducing symptoms and reversing all types of disease conditions.

Dr Bob Beck discovered the 1000hz Frequency which he called the "magic frequency" that was found to work where other frequencies were struggling with results, so remember to incorporate BT-6 1000hz Non Pulsing Program when other programs are not performing. It is best to use the BT-6 in conjunction with another program. For example use the [Lymph] Program for 6 hours on an overnight treatment, and switch to the BT-6 for the final 3 hours. Dr Beck would incorporate the BT-6 Machine with other machine programs for effective treatment of clients.

You May use the Bob Beck Professional in conjunction with other Bioresonance machines, such as the worlds most powerful Rife Machine (Dr Royal Rife) "The Rife Digital Professional", People have used both the Rife Digital Professional and the Dr Bob Beck Professional simultaneously and have said the results were amplified, the symptoms and disease were eliminated faster, or the Brain Treatments took on a more powerful effect.

In the Pages which follow, we will see the Programs in detail, followed by a suggestion of the type of accessory which will deliver the Program effectively to the body for best treatment results.

See FAQ, Frequently asked Questions, and answers about the Machine, treatments, and the Colloidal Silver program

PROGRAM # 1, Dr Beck (BT-4)

The Original Bob Beck Machine, the BT4 was a simple one frequency device, and was the first machine produced by Dr Beck, and so therefore is the first frequency provided on the Dr Beck Professional Machine.

The Bob Beck designed blood electrifier is a mutivibrator circuit that produces a square wave output at 3.92 Hz. creating a small current in the 50-100 microampere range to flow in both the Ulnar and Radial arteries located on underside of the forearm just behind the wrists and also the inner leg.

The information on using Beck's blood electrification device (sometimes called a "blood electrifier" or the Beck zapper) is provided in the 1998 Nexus magazine article titled: The Beck Zapper: Renewed Hope Against AIDS.


BT-4 gives a continuous output of 3.92Hz (Always On), there is no Pulse on this program. Alternatively use Program # 13, PULSE (Blood Pulser) should you require this same frequency set 3.92Hz in Pulse Mode. Use Machine for 7-9 Hours overnight for one week for a complete and thorough Blood Cleansing

ACCESSORIES TO USE WITH PROGRAM BT-4

[GELPADS] on the inner lower Leg on the arteries



[GELPADS] on the wrist over the arteries

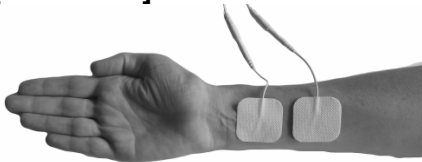

PROGRAM # 2, Dr Beck (BT-5) also known as the BT PLUS

The BT Plus uses the frequencies of 111 Hz pulse, and was referred to as the "The Black Box - Secret ... Treatment of Rock Superstars" in OMNI Magazine. The BT-PLUS black box was used by Dr. Margaret A. Patterson, MD., to cure British rock star Pete Townsend of "The Who" of his addiction to heroin.

Use the BT-5 frequency for Drug Addiction
Or use frequency #14 [HABIT] to assist with alcohol addiction issues.

BT-5 Program works well for Drug addiction specific problems, such as addiction to Nicotine, Heroin, Cocaine, Barbiturates and addictive Pharmaceutical drugs. BT-5 Program requires a minimum of one month overnight treatment 7-9 hours for positive results.

ACCESSORIES TO USE WITH PROGRAM BT-5

[Earclips] For Brain Treatment, Clip to Ear-Lobes



[Gelpads] For Brain Treatment, Attach to Temples



PROGRAM # 3, Dr Beck (BT-6)

1000hz. This is the "magic frequency" that was found to work where other's failed. Dr Beck often used the BT5 and BT6 frequencies together for treatments. BT5 - 111Hz & BT6-1000hz

BT-6, "The Magic Frequency" works to boost the effectiveness of other frequencies when certain frequencies are not working.

So an ideal program suggestion would be to select your usual frequency from the Dr Beck machine, and use for the suggested time; then follow this program with the BT-6 to boost the effectiveness of the original program

BT-6 works in unity with all other programs, and would be the second program of choice to your usual program

ACCESSORIES TO USE WITH PROGRAM BT-6

[Earclips] For Brain Treatment, Clip to Ear-Lobes



[Gelpads] For All Treatments



PROGRAM # 4, Dr Beck (BT-7)

Pulsed 0.5Hz. Voltage modulation corresponds to the Lower Band of the Delta brainwave frequency. This frequency, in terms of EEG brainwave, is associated with deep REM sleep. BT-7 can also be used in combination with the [DELTA] Program #9 which uses the full spectrum of 10 frequencies 0.01Hz - 1.5Hz

Dr Beck selected the 0.5Hz on the BT-7 specifically for Brain Treatments and assist with insomnia conditions. Slow Wave Sleep or SWS, is the deepest of sleep states and it plays a vital role in health and well being. During this phase of the sleep cycle, the brain begins producing very slow, large Delta waves. Even if your lifestyle doesn't allow for the luxury of a full eight hours of sleep, a few hours of on the BT-7, 0.5Hz Program will allow your brain to receive the restorative sleep it needs. Other Spectrums of the Delta Range of frequencies can be accessed using the [DELTA] Program # 9

ACCESSORIES TO USE WITH PROGRAM BT-6

[Gelpads] Use with BT-7 and attach Gelpads to Temples



[Earclips] For BT-7 Treatment, Clip to Ear-Lobes



PROGRAM # 5, SCHUMMAN EARTH FREQUENCY [SCHUM]

In the 1970s Dr Beck discovered the Schumann frequencies which were are a set of spectrum peaks in the extremely low frequency (ELF) portion of the Earth's electromagnetic field spectrum. Schumann resonances are global electromagnetic resonances, excited by lightning discharges in the cavity formed by the Earth's surface and the ionosphere. The Schumann Earth Frequency can clear negative energy and emotions within the body, The Frequency can ground (earth) the body and facilitate a better platform for healing.

By tuning into the Schumann Resonance, the body benefits from this healing frequency and can more effectively rebuild and repair itself. The [SCHUM] 7.83Hz Program has been shown in scientific studies to be an essential requirement for physical and psychological health. The Schumann Resonance in the body and helps to recover autonomic nerve balance, remove stress, and intensify concentration.

ACCESSORIES TO USE WITH PROGRAM SCHUMMAN EARTH

[Gelpads] Use with [SCHUM] and attach Gelpads to Temples



[Earclips] for Ear-Lobes

[Flex] For Hands or Feet



PROGRAM # 6, ALPHA WAVE FREQUENCY [ALPHA]

Alpha Waves Exist between 8.0Hz and 12.0Hz, on the Dr Bob Beck Professional, the [ALPHA] Program will run 10 Alpha Frequency programs continuously within the entire ALPHA Bandwidth.

Dr Bob Beck Used all 5 Brainwaves in his study of healing the Brain and the Body, Dr Beck discovered that Alpha Waves promotes a state of meditation and relaxation

Alpha is the bridge between Beta and Theta. In an Alpha state your brainwaves move at 7 and 14 cycles per second. The Alpha Program promotes a meditative state of mind. Alpha is the frequency between our conscious thinking and subconscious mind. Alpha Waves resolve Anxiety, High Stress, Insomnia, OCD. Alpha has been known to take away pain and is useful in healing. Alpha Program is not suitable for memory recall and learning, so use for Meditation and relaxation

ACCESSORIES TO USE WITH PROGRAM [ALPHA] WAVE

[Gelpads] Use with [ALPHA] and attach Gelpads to Temples



[Earclips] for Ear-Lobes

[Flex] For Hands or Feet



PROGRAM # 7, BETA WAVE FREQUENCY [BETA]

BETA Waves Exist between 12.0Hz and 30.0Hz, on the Dr Bob Beck Professional, the [BETA] Program will run 10 Beta Frequency programs continuously within the entire BETA Bandwidth.

Dr Bob Beck Used all 5 Brainwaves in his study of healing the Brain and the Body, Dr Beck discovered that Beta Waves promotes a state of mind activity.

Beta waves have a frequency of 14–28 cycles per second. Beta Waves are associated with normal waking consciousness and a heightened state of alertness, logic and critical reasoning. Having the right amount of Beta allows us to focus to our task whether at school, work or sports performance. [BETA] Program resolves ADHD, daydreaming, depression and poor cognition. Beta waves will permit Conscious focus, memory, problem solving

ACCESSORIES TO USE WITH PROGRAM [ALPHA] WAVE

[Gelpads] Use with [ALPHA] and attach Gelpads to Temples



[Earclips] for Ear-Lobes



[Flex] For Hands or Feet



PROGRAM # 8, GAMMA WAVE FREQUENCY [GAMMA]

GAMMA Waves Exist between 25.0Hz and 100.0Hz, on the Dr Bob Beck Professional, the [GAMMA] Program will run 10 GAMMA Frequency programs continuously within the entire GAMMA Bandwidth. Dr Bob Beck used all 5 Brainwaves in his study of healing the Brain and the Body, Dr Beck discovered that GAMMA Waves promotes a state of mind activity.

Gamma Brain Waves are the fastest brainwaves (high frequency) and the most recently discovered brain wave state. In this state, your brainwaves loop between 40 and 5,000 cycles a second. [GAMMA] Program promotes higher processing brain tasks as well as cognitive functioning; This is important for learning, memory and information processing. Gamma Waves resolve ADHD, Depression, Learning Disabilities. Gamma waves will allow greater Learning ability, Cognition, information processing, perception and night use will put the Brain into REM sleep

ACCESSORIES TO USE WITH PROGRAM [ALPHA] WAVE

[Gelpads] Use with [ALPHA] and attach Gelpads to Temples



[Earclips] for Ear-Lobes



[Flex] For Hands or Feet



PROGRAM # 9, DELTA WAVE FREQUENCY [DELTA]

DELTA Waves Exist between 0.01Hz to 1.5Hz, on the Dr Bob Beck Professional, the [DELTA] Program will run 10 DELTA Frequency programs continuously within the entire DELTA Bandwidth. Dr Bob Beck used all 5 Brainwaves in his study of healing the Brain and the Body, Dr Beck discovered that DELTA Waves promote The Deep Sleep State.

Delta are Low Frequency Brain Waves which exist at 0-4 cycles per second. Delta is experienced in a deep dreamless sleep or deep transcendental meditation. Delta is usual mostly in Infants. Deep sleep is important for the healing process – as it's linked with deep healing and regeneration. Delta Waves resolve an Inability to rejuvenate body, inability to revitalize the brain, poor sleep and insomnia. Delta waves will boost the Immune System, promote natural healing and a restorative sleep and deep sleep in infants and adults. Will also assist with infants who have sleeping problems.

ACCESSORIES TO USE WITH PROGRAM [ALPHA] WAVE

[Gelpads] Use with [ALPHA] and attach Gelpads to Temples



[Earclips] for Ear-Lobes



[Flex] For Hands or Feet



PROGRAM # 10, THETA WAVE FREQUENCY [THETA]

THETA Waves Exist between 4Hz to 7Hz, on the Dr Bob Beck Professional, the [THETA] Program will run 10 THETA Frequency programs continuously within the entire THETA Bandwidth. Dr Beck discovered that THETA Waves promotes a very deep state of relaxation

A Theta state is a very deep state of relaxation. This is the state used in hypnosis and dreaming. In it, the brainwaves are slowed to a frequency of 4 - 7 cycles per second. Theta brainwaves can be thought of as the subconscious; they govern the part of our mind that is layered between the conscious and the unconscious. They hold memories and sensations. They also govern our attitudes, beliefs and behavior. They are always creative, inspirational and characterized by very spiritual sensations.

Theta Brain Waves occur most often in sleep but are also dominant during deep meditation. In theta we are in a dream; vivid imagery, intuition and information beyond normal consciousness awareness. It helps us improve our intuition, creativity, and makes us feel more natural. This is also the Flow State Zone.

Theta Brain Waves resolve Anxiety issues, poor emotional awareness & stress. Theta brain waves promote creativity, emotional connection, intuition and relaxation.

ACCESSORIES TO USE WITH PROGRAM [ALPHA] WAVE

[Gelpads] Use with [ALPHA] and attach Gelpads to Temples




PROGRAM # 11, [Lymph] Lymph Treatment Program

The Lymph Treatment Program exists on the 3.54Hz - 3.99Hz Frequency bandwidth on the Dr Beck Machine. The [LYMPH] Program will run 10 LYMPH Frequency programs continuously. The Lymph is a complicated network of fluid-filled nodes, vessels, glands and organs, the lymphatic system touches almost every part of the body. The Lymph is one of the most important systems of the human body.

The lymphatic system's main function is to cleanse toxins and protect against harmful invaders. It tackles toxins that are introduced to the body from both external means (food, air, personal care products, water) as well as internal ones (damaged protein, Cell & Metabolic waste), making it a key detoxification pathway.

Through lymph nodes and the lymphatic network, your immune cells can travel around fighting pathogens, such as bacteria and mold, and preventing disease and infection. This is why keeping your lymphatic system functioning properly is directly related to the overall health of the body: a stronger lymphatic system means a more resilient and reactive immune response and defense. Use LYMPH Program overnight for 7-9 hours with GEL-PADS

USE GELPADS WITH THE PROGRAM [LYMPH]

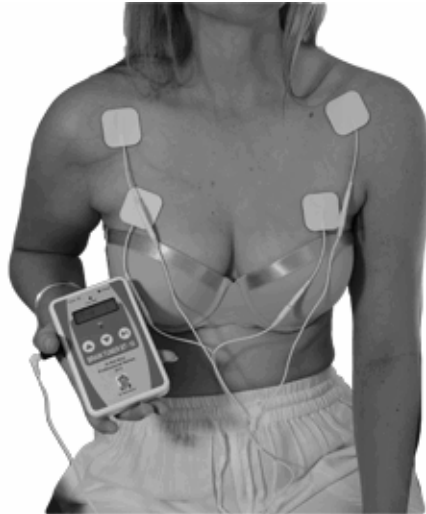

PROGRAM # 12, [Blood] Blood Electrifier Program

The Blood Treatment Program exists on the 100Hz - 150Hz Frequency bandwidth on the Dr Beck Machine. The [BLOOD] Program will run 10 BLOOD Programs continuously. Use for 7 - 9 Hours a night.

Dr Beck researched using electricity in the blood to disable microbes his research describes two blood electrification methods as being able to: "provide electric current flow through the blood sufficient to render the bacteria, virus, parasites and/or fungus ineffective to infect or affect normally healthy cells." Dr Beck then invented a simpler non-invasive means of getting electricity into the blood. He then did off-shore testing of the device to verify its effectiveness (mostly against viral disease since viruses are always circulating in the blood). Since then, Bob Beck has personally seen hundreds of people helped by his blood electrifier

Bob Beck has seen dozens of his friends use these devices and regain so much immune strength that they didn't lose a single day's work in over three years from getting sick! Bob Beck also knows of an MD who claims over two hundred documented cancer remissions by using only blood electrification via the Blood Electrifier.

ACCESSORIES TO USE WITH PROGRAM #12 [BLOOD]

[GELPADS] on the inner lower Leg on the arteries



[GELPADS] on the wrist over the arteries



PROGRAM # 13, [PULSE] Blood Pulser Program

The Blood Pulser Program is the 3.92Hz Frequency bandwidth on the Dr Beck Machine. The [BLOOD] Pulser will run this one BLOOD Pulse Frequency program continuously. Use Program 13 for 7 - 9 Hours overnight for one week for a complete and thorough Blood Pulsing


The Pulser aids transport of immune cells to the bloodstream to destroy the microbes deactivated by electricity in the blood, and can assist with all kinds of diseases including AIDS, prostate cancer, rheumatoid arthritis, Hepatitis A and B, Herpes B, Chronic Fatigue Immune Dysfunction Syndrome, cervical cancer, Lyme disease, and lupus. All of these are caused by viruses and bacteria which is the forte of the Blood Electrifier Blood Pulser, which means that the list of diseases these devices may reverse could be very long. Dr Bob Beck also knows of an MD who claims over two hundred documented cancer remissions by using only blood electrification via the Blood Electrifier. Bob Beck has personally helped more AIDS patients than any other type and feels that his devices are 100% effective if the patient uses them exactly as directed. In the early years of his electromedicine research he funded a study with two medical doctors and 8 AIDS patients

ACCESSORIES TO USE WITH PROGRAM #13 [PULSE]

[GELPADS] on the inner lower Leg on the arteries



[GELPADS] on the wrist over the arteries

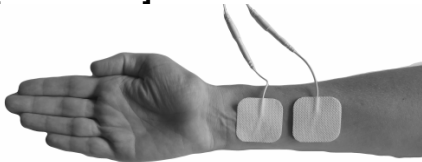

PROGRAM # 14, [HABIT] ADDICTIONS TREATMENT

The HABIT (Addictions) Program is the 100Hz Frequency bandwidth on the Dr Beck Machine. The [HABIT] Program will run a single 100Hz Frequency program continuously. Use Program #14 for 7 - 9 Hours overnight for one week to one month. Use for extended periods if required.

The 100hz Frequency works effectively for Alcohol specific addictions, however you may combine this frequency with the General Dr Beck BT-5 Addictions program

Concerning Alcohol addiction; the brain ceases making some of its own neurotransmitters, for example, its own pain-killers called endorphins. Putting a weak electromagnetic signal into the body, especially the head area, stimulates the brain into producing its own neurotransmitters again, including endorphins. Amazingly, the craving for alcohol and drugs disappear in three to ten days.

ACCESSORIES TO USE WITH PROGRAM 14, [HABIT]

[Gelpads] Use with [ALPHA] and attach Gelpads to Temples



[Earclips] for Ear-Lobes



[Flex] For Hands or Feet



PROGRAM # 15, [BACTERIA]

The [BACTE] BACTERIA Program is the 920Hz Frequency bandwidth on the Dr Beck Machine. The Bacteria Program will run a single 920Hz Frequency program continuously. Use Program #15 for 7 - 9 Hours overnight for one week to one month. Use for extended periods if required.

Use The Bob Beck Machine to fight Bacteria and common Bacterial Infections. The BACTERIA Program works in a similar way as Homeopathy: using natural principles to heal the body. Everything around us already has an electrical charge and thus a certain electrical frequency - we are always and everywhere exposed to it. If you have a specific Bacteria, it lives in a frequency field. When you transmit the Bacteria Frequency into your body, the bacteria gets over-stimulated and dies. Use Bacteria Program to restore vitality and health. The program is producing a harmonic resonance in the body.

ACCESSORIES TO USE WITH PROGRAM 15, [BACTERIA]

[Gelpads] Use with [ALPHA] and attach Gelpads to the Neck area



[Earclips] for Ear-Lobes



[Flex] For Hands or Feet



PROGRAM # 16, [VIRUS]

The [VIRUS] Program is the 250Hz Frequency bandwidth on the Dr Beck Machine. The Bacteria Program will run a single 220Hz Frequency program continuously. Use Program #16 for 7 - 9 Hours overnight for one week to one month. Use for extended periods if required

Viruses are susceptible to the same kind of mechanical excitation as seen in the Bacteria program. Pulses of Electronic Frequencies will induce destructive vibrations in the "outer shell" of a Virus

The specific Bandwidth of the frequency creates an impulsive push which gets the virus shaking, which actually induces a destructive vibrations in the virus shells. You can use the [VIRUS] program with the Bacteria Program

ACCESSORIES TO USE WITH PROGRAM 16, [VIRUS]

[Gelpads] Use with [VIRUS] and attach Gelpads to the Neck area



[Earclips] for Ear-Lobes



[Flex] For Hands or Feet



PROGRAM 17 [HIV] HIV/AIDS AND SYMPTOMS PROGRAM

The [HIV] Program is the 150Hz Frequency bandwidth on the Dr Beck Machine. The HIV Program will run a single 150Hz Frequency program continuously. Use Program #17 for 7 - 9 Hours overnight for one week to one month. Extended periods are required.

"Bioelectrification" or "blood electrification", a simple electronic therapy that was discovered to stop the replication of the virus that causes AIDS. Subsequently, Bob discovered that many other serious health conditions, including malaria, responded favorably to blood electrification. Many thousands, perhaps hundreds of thousands, of people, worldwide, have benefited enormously from the use of blood electrification due to the lecturing efforts of Robert C. Beck.


Use the Gel Pads and the Blood Cleanse methods for the treatment of HIV/AIDS. This is the exact method used by Dr Beck to eliminate this Virus in the blood. Longer sessions are imperative successful treatment

ACCESSORIES TO USE WITH PROGRAM #17 [HIV]

[GELPADS] on the inner lower Leg on the arteries



[GELPADS] on the wrist over the arteries

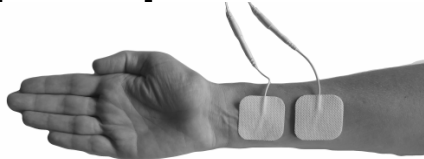

PROGRAM 18 [HERPES] HERPES AND SYMPTOMS PROGRAM

The [HIV] Program is the 950Hz Frequency bandwidth on the Dr Beck Machine. The HERPES Program will run a single 950Hz Frequency program continuously. Use Program #18 for 7 - 9 Hours overnight for one week to one month. Extended periods are required.

"Bioelectrification" or "blood electrification", a simple electronic therapy that was discovered to eliminate Herpes Virus. Dr Bob Beck discovered that many other serious health conditions, including malaria, responded favorably to blood electrification. Many thousands, perhaps hundreds of thousands, of people, worldwide, have benefited enormously from the use of blood electrification due to the lecturing efforts of Robert C. Beck.


Use the Gel Pads and the Blood Cleanse methods for the treatment of HERPES (All types). This is the exact method used by Dr Beck to eliminate this Virus in the blood. Longer sessions are imperative successful treatment

ACCESSORIES TO USE WITH PROGRAM #17 [HIV]

[GELPADS] on the inner lower Leg on the arteries



[GELPADS] on the wrist over the arteries

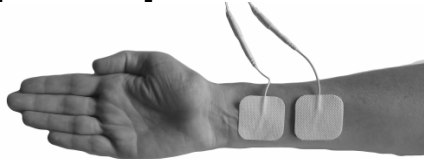

PROGRAM 19 [FLU] COLD AND FLU PROGRAM


The [FLU] Program is the 850Hz Frequency bandwidth on the Dr Beck Machine. The FLU Program will run a single 850Hz Frequency program continuously. The influenza FLU virus, as well as other common viral infections like the common cold, often cause irritation of the throat as part of the entire constellation of the illness they cause, which also includes nasal congestion, cough, aches and pains, and, sometimes, fever. Use Program #19 for 7 - 9 Hours overnight for a week or more. Use the GelPads on the neck to assist with Cold and Flu conditions by targeting Bacteria and Virus in the Neck Lymph Glands region. Use the Dr Beck Atomizer for Lung infections

The FLU program can be used also with the [VIRUS] program as well as the [BACTERIA] program

[Gelpads] Use with [FLU] and attach Gelpads to the Neck area



Gel Pads on the Lymph nodes for FLU & Body Fever & Night Sweats

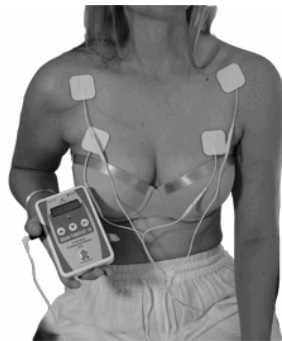

PROGRAM 20 [SILVER] COLLOIDAL SILVER GENERATOR

The [SILVER] Program runs the Alpha Frequency bandwidth on the Dr Beck Machine. Creating Silver in the ALPHA Frequency produces the most bio-active Colloidal Silver available on the market. All features are integral to the Dr Bob Beck Protocol, and you can do this at home on your own machine with a simple press of a button

Use the Colloidal Silver Program as the second most important contribution to your therapy. Use the Colloidal Silver in conjunction with all Dr Beck Protocols and programs. The Colloidal Silver Apparatus is included with your Bob Beck Professional Machine, and is in the Bag. You can also Purchase a separate device, The Dr Beck Atomizer (accessory) for Lung Specific treatments with the Silver you make on the Beck Machine.



Clip the Colloidal Silver Maker on the Cup of water filled to 3/4 full. Plug the Cable into the Dr Beck Machine. Select Program # 20 [SILVER]. Press (▶||) Play, and you will see a small cloud appear in the water, indicating the production of silver at 10 ppm at .0010 microns. The Program will turn off once the Silver is produced. If you wish to double the strength, then Press (▶||) Play one more time. The more times you Press (▶||) Play, the greater the silver PPM. PROGRAM 21 [CANCER] CANCER TREATMENT PROGRAM

The [CANCER] Program is the 20HZ - 780Hz Frequency bandwidth on the Dr Beck Machine. The CANCER Program will run a Set of 10 Frequency programs within the 20Hz - 780Hz Frequencies continuously. Use Program #21 for 7 - 9 Hours overnight for one week to one month. Extended periods are required.

"Bioelectrification" or "blood electrification", a simple electronic therapy that was discovered to assist in the treatment and remission of all Cancers.

GELPADS TO USE WITH PROGRAM #21 [CANCER]


- Use the Blood Cleanse methods for the treatment of Blood type Cancers. Place Gel Pads over the arteries of the wrist & foot
- Use Gel Pads for treatment of Lymph Cancer and Breast Cancer. Place Gel Pads on the Lymph Nodes and Breast
- Place Gel Pads over the location of any organ cancer
- Place Gel Pads on the temples for Brain type Cancers

OTHER ACCESSORIES TO USE WITH PROGRAM #21 [CANCER]

[Earclips] for Ear-Lobes



[Flex] For Hands or Feet

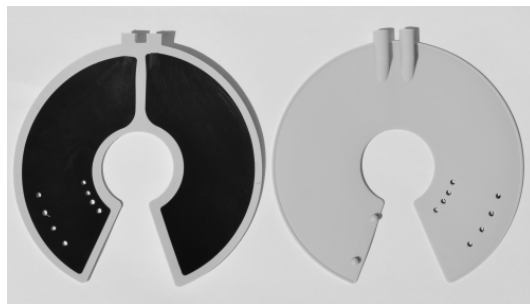

Use Breast Pads if required

PROGRAM 22 [ALLERGY] TREATMENT PROGRAM

The [ALLERGY] Program is the 40HZ - 250Hz Frequency bandwidth on the Dr Beck Machine. The ALLERGY Program will run a Set of 10 Frequency programs within the 40HZ - 250Hz Frequencies continuously. Use Program #22 for 7 - 9 Hours overnight for one week to one month. Use for extended periods if required.

The Dr Beck Allergy frequencies target and stimulate points along the spine from the neck to the sacrum. Lymph nodes are located throughout the body but the largest groupings are found in the neck, which will be the location of your GEL-PADS for Allergy treatment.

GELPADS TO USE WITH PROGRAM #22 [ALLERGY]

- Place Gelpads on the upper neck, on the point of the Lymph Glands under the neck
- Secondary Gelpads position on the face, on each side of the nose.

OTHER ACCESSORIES TO USE WITH PROGRAM #22 [ALLERGY]


[Earclips] for Ear-Lobes



[Flex] For Hands or Feet



Make Colloidal Silver and breath deeply; the Atomized Silver into the Lungs, and then exhale out through the nose.

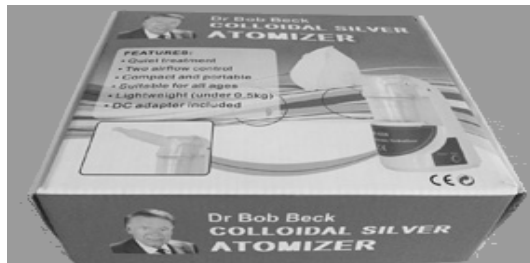

FAQ. [Frequently Asked Questions]

I cant find a frequency for my illness, what other frequency can I use on this machine?

ANSWER: The Original Dr Bob Beck Machine only used 4 Programs for virtually all treatments, so we have totally upgraded the machine and provided 22 Programs including all of the original Dr Beck BT Series. However if you feel that your disease or illness does not appear to match the Programs of Dr Beck, then Dr Beck recommended the BT-6 as the "Magic Machine" as was described in this early literature. He Recommended using BT-5 and BT-6 Together as a means of healing hundreds of different diseases. So In an overnight treatment 7-9 hours, one could use BT-5 for 4 1/2 hours, and then BT-6 for 4 1/2 hours. Also consume Colloidal Silver daily with the program.

How Much Silver Should I drink? And How long can I use the Colloidal Silver Atomizer Machine?

ANSWER: Dr Beck Recommended the drinking of about 4-5 regular size cups of colloidal silver at 10ppm strength per day to assist with a disease treatment program. The Dr Beck Atomizer can be used for a period of 10-15 minutes, for 3 or 4 times per day. Should the Lung, or Throat Bacterial Infection be severe, Colloidal Atomizer treatments can be increased.

What is the Frequency to heal ADHD condition?

ANSWER: Use the BETA wave Frequency for ADHD. You can also incorporate BT-5 for 1 hour in the treatment program. Use the machine overnight, with flex straps on the feet for 7-9 hours, you can use the GelPads on the temples during the day if necessary Beta Waves are associated with normal waking consciousness and a heightened state of alertness, logic and critical reasoning. Having the right amount of Beta allows us to focus to our task whether at school, work or sports performance. [BETA] Program resolves ADHD, daydreaming, depression and poor cognition. Beta waves will permit Conscious focus, memory, problem solving

What's the frequency to heal Anxiety, Depression and Mood Disorder?

ANSWER: For Anxiety Depression and Mood Disorders Treatment use Program # 5 SCHUMMAN EARTH FREQUENCY [SCHUM] Overnight for 7-9 hours. The Schumann Resonance in the body and helps to recover autonomic nerve balance, remove stress. During the day hours, use the Program # 6 [Alpha] Wave Frequency. Alpha Waves promotes a state of meditation and relaxation. [ALPHA] Can also be incorporated into the night program if required.

What's the Frequency for Better Learning and memory retention?

ANSWER: Alpha waves are the most conducive to learning. The Alpha state (about 8-10Hz) occurs when brain activity slows just below the normal waking state of Beta (11-25hz). In Alpha, the mind and body are relaxed but a level of focus is easily maintained. While in this deeper brainwave state, information can be processed consciously without as much "mental activity" to interfere with it. There is also a greater link between the conscious and subconscious mind in Alpha, meaning that while one is consciously learning, their brain is also unconsciously processing what one is learning. Modern brain science has shown that only 1 of 6 of our brain's processing methods happen on the conscious level, so there is certainly a benefit in engaging the deeper parts of the brain, while taking in information. Theta waves can be added to this Learning Program. Theta is known for Super Learning and Long Term Memory

What's the best frequency for General Meditation?

ANSWER: For Normal or General Meditation, Theta waves are the ones that help you to sleep and deeply relax. Those who practice general meditation forms find themselves in a Theta State. When you are in Theta meditation the low frequency theta waves in your brain can help produce intuition, creativity, strong emotional connections as well as lower your stress and anxiety. As you can imagine there are so many fantastic health benefits that can be experienced if you are able to effectively achieve theta meditation and experience these very low brain frequencies.

What are the best frequencies for Deep Meditation?

ANSWER: Deep Meditation or Transcendental Meditation works best with Delta waves, which are the deepest level of meditation. A delta wave is a brain wave with a frequency of oscillation between 0.5–4 Hz. some people who practice deep meditation can remain conscious while in delta-sleep. the delta state is associated with “no thinking” removing thoughts clearing the mind. Delta waves are linked to the unconscious part of our mind. The Best suggestion for Meditation is the Program DELTA.

Can Gamma Waves assist with Schizophrenia and Bi-Polar?

ANSWER: Gamma Waves are known to assist with Schizophrenia They assist with increased mental activity, including perception, fear, problem solving and consciousness, stabilizing Bi-Polar issues. Stabilizing mood swings, and balancing certain brain chemicals (neurotransmitters) that control emotional states and behavior. Gamma Waves can help to treat mania and to prevent the return of both manic and depressive episodes in bipolar disorder. They may also help for treat the mood problems associated with schizophrenia, such as depression.

Why is the Charge is not lasting as long anymore on the Beck Machine? I have had the machine for 3 months

ANSWER: Try to Recharge the machine for about 8 hours; However Rechargeable batteries only have so many recharge cycles. Do not believe the Recharge Claims on the Battery Packages, some manufactures claim "2000" recharges etc. Unfortunately the truth is that Rechargeable batteries do loose charge over time, and may need to be replaced. Purchase a second set of batteries if you are a "Heavy" Beck Machine user. Department stores now sell packets of 4pcs AA 2000+ mAh with a wall charger.. Its great having a portable battery powered Beck Machine, but not great if batteries go flat, when you really need a treatment, so have a second pack handy!

Is Your Beck Professional as powerful as my Rife Digital Professional for Cancer Treatments?

ANSWER: *The Rife Digital Professional uses specific Dr Rife Machine frequencies on the MOR Carrier Wave Sideband, and therefore requires Mains Power to run, and delivers an intensely powerful range of frequencies for healing treatments. It is an excellent tool for Cancer Treatment. The Dr Beck Professional Machine is very powerful in its own right' Producing a much lower HZ frequency bandwidth inherent to the Dr Beck Protocol, and is also proven very effective in Cancer treatments. If you are currently using the Rife Digital Professional, you can use the Dr Beck Machine as the portable treatment device when you are out and about. Remember the Dr Beck Machine produces Colloidal Silver and can be used with the Dr Beck Atomizer machine, and the Silver Program is a highly important composition of the Dr Beck Cancer Treatment Protocol.*

I want to use the Beck Machine at my workplace but I am a little embarrassed to have straps dragging around my workstation, what is the best method for portable treatments?

ANSWER: *Gelpads will assist you here. In a work, or school or even on a train, bus or walking down the street, place the Dr Beck Professional Machine in an inside pocket, and bring the wires and Gelpads through a sleeve onto the skin, or even have the Gelpads attached to your body, out of sight.*

I have a serious chronic disease for about 7 years, how long will It be that I see positive results and remission?

ANSWER: *The Longer the Disease condition, the longer the treatment time will be. Do expect to see positive results in a week or two, however do not stop treatments, but continue them for months. Remember to do the extended "overnight" treatments of 7-9 hours, with the straps attached to the feet, or the Gel Pads on the arterial area of the inner leg each night. The more the treatment time, the faster the positive results. Also add the Colloidal Silver to your program each day, as prescribed by the Dr Bob Beck Protocol.*

I saw this other machine called the SOTA on the Web which says it does the Beck Protocol, but it only has 1 or 2 frequency and no

Colloidal Silver, but it looks like its the most well known machine on the Web for Beck Protocol, Is this an older machine or something?

ANSWER: Yes a "SOTA" machine has been around for 8 years and has not really been updated as far as I see, and does not Provide the Dr Beck Protocol Silver Program on the actual machine. I believe the "SOTA" does not provide the full range of the BT Series of programs as prescribed by Dr Beck, perhaps only 2 or 3 frequencies. The Bob Beck Professional machine however, has 22 Programs built in, which adhere precisely to the Dr Beck Protocol, Including all Brain Tuner Programs and 5 Brain Wave frequencies plus Schumann Earth Frequency, Blood Pulsar, Blood Cleanser and disease frequency sets; and yes we do include the Silver Program on the machine at Program position #20. Our Company also provides a huge number of accessories to make the Dr Beck Protocol work, such as GelPads to be placed at exactly the body points where the programs are required to be transmitted for healing. We even developed special Breast Pads to assist women with Cancer and Breast Lymph treatments. Our company developed the first Colloidal Silver Program which produces silver at the Alpha Frequency bandwidth on the Dr Beck Machine. Colloidal Silver in the ALPHA Frequency produces the most bio-active Colloidal Silver available on the market. All features are integral to the Dr Bob Beck Protocol, and you can do this at home on your own machine with a simple press of a button

Where Can I find out more about Dr Beck and his protocol?, and also where can I buy the Dr Beck Colloidal Silver Atomizer machine?

ANSWER: Further information on Dr Beck can be found on the Official Dr Bob Beck Organization website: drbeck.net. You will likely be able to buy the Dr Beck Atomizer machine there also. Check the CART page on drbeck.net, as we are constantly developing new products and accessories for the Dr Beck Professional Machine.

Does the Dr Beck machine have Japanese Language?

ANSWER: We are working on a multiple language edition of the Dr Beck Machine for a future release, meanwhile use Google Translate to Translate the English menu to another language

Will the Beck Machine Charger work in UK and Australia?

ANSWER: Yes, the Charger for Beck Machine is supplied with Socket Plugs for every country. The Charger is a Switching power supply so will work in USA and UK (110v - 220v - 240v)

Can I Charge the Device and Use it at the same time?

ANSWER: No, you can only use the machine after it has been charged up with the battery charger. If you are a heavy user of Beck Machine, then purchase a second pack of 4 x AA 2000+ mAh Batteries with an external charger, so you never "run out of juice"

How Long Do I need to use the Machine?

ANSWER: This depends upon the type of illness or Disease being treated. The more Chronic condition or serious disease require months of treatment with the Flex Straps attached to the feet overnight for 7-9 hours, or perhaps with the Gel-Pads attached to the area of the body where the disease condition originated. If you are only using the machine for a Detox or Cleanse, or to resolve a short term Allergy problem, then a one week treatment is sufficient. Remember for General health Dr Beck recommended everyone to drink colloidal silver every week. You can make this on the machine with program #20 [SILVER]

Is Colloidal Silver a new invention?

ANSWER: Silver has been used for thousands of years as a healing and preventive health product. Treating wounds with silver was common from the 1800s to the mid-1900s when the use of antibiotics took precedence in the medical field. However, silver is still utilized in many medical circumstances. Hospitals use a silver treatment to dress scars and wounds especially wounds that resist healing. Catheters are lined with silver to prevent infection and water is purified with silver. Now, silver is coming to the at-home health care market

How do you use Colloidal Silver?

ANSWER: *Colloidal silver can be used as a dietary supplement as well as for the treatment of a wide range of ailments. Colloidal silver is used to treat infections due to yeast; bacteria (tuberculosis, Lyme disease, bubonic plague, pneumonia, leprosy, gonorrhea, syphilis, scarlet fever, stomach ulcers, cholera); parasites (ringworm, malaria); and viruses (HIV/AIDS, pneumonia, herpes, shingles, warts).*

Colloidal silver is also used for lung conditions including emphysema and bronchitis; skin conditions including rosacea, cradle cap (atopic dermatitis), eczema, impetigo, and psoriasis; and inflammation (sometimes due to infection) of the bladder (cystitis), prostate (prostatitis), colon (colitis), nose (rhinitis), stomach (gastritis), tonsils (tonsillitis), appendix (appendicitis), and sinuses (sinusitis).

Other uses include treatment of cancer, diabetes, arthritis, lupus, chronic fatigue syndrome, leukemia, hay fever and other allergies, trench foot, and gum disease. Colloidal silver is also used to prevent flu, H1N1 (swine) flu, and the common cold. Some women take colloidal silver during pregnancy to aid the baby's growth and health as well as the mother's delivery and recovery.

Colloidal silver is applied directly to the skin for acne, burns, eye infections, fungal infections, throat infections, skin infections, and Staphylococcus infections.

How Does Colloidal Silver work?

ANSWER: *Colloidal silver works as a catalyst. A catalyst is best described as a substance that brings about (causes) a reaction or occurrence, without itself participating or being consumed. Silver, in its atomic state, has the capacity to absorb oxygen and act as a catalyst to bring about oxidation. Employing a simple catalytic reduction/oxidation reaction, colloidal silver will react with any negative charge presented by the organism's transport or membrane proteins and deactivate them. Colloidal silver can kill certain germs by binding to and destroying bacteria proteins.*

Can I mix Essential Oils with the Colloidal Silver into the Dr Beck Atomizer machine?

ANSWER: The Dr Beck Atomizer uses a smaller nozzle to create a finer mist, and uses a different pressure mechanism to suspend the tiny silver particles (or colloids) more completely in the mist droplets. The Dr Beck Atomizer was developed for Colloidal Silver only, so we suggest not to mix essential oils with the colloidal silver water to make vapor, as the oils could block the tiny vapor outlets. Use the Dr Beck Atomizer for treatment of Cold and Flu, Chest Infections, Lung Infections or for treatment of serious Lung diseases. The Dr Beck Atomizer is an accessory, so is purchased separately to the Dr Beck Professional machine package.

Can I use the Flex Straps and Velcro Straps on the wrists?

ANSWER: YES, In the initial testing stages of the machine, people used the FLEX and VELCRO on the wrists and had no sensitivity issues. You can use the straps on the wrists with the Beck Machine


What is the Warranty Period on the Dr Beck Machine

ANSWER: The Dr Beck machine is Guaranteed for 12 months. Warranty will not cover Broken LCD Screen caused by dropping or other damage, or damage to the actual body cause by dropping. The Batteries are not under Guarantee, however we do use Good Brand Batteries such as SONY for extended life and recharge. We do not Guarantee damage to machine caused by leaking batteries so make sure your AA battery replacements are a good type and only use rechargeable. Accessories such as GEL-PADS are not under warranty as they are multiple use and throw away items. Leads and Wires have a 30 day warranty.

Notes

INTERNATIONAL SALES AND INQUIRIES
USA / CANADA / and World

Sales and Service
healthproducts2@gmail.com
www.drbeck.net
USA ☎ 206-708-6408
International ☎ 0011-1-206-708-6408

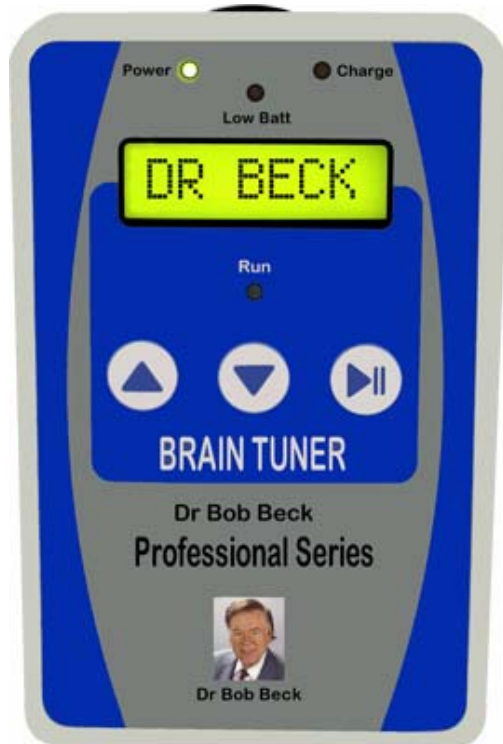

© February/ 2017, Dr Bob Beck Organization USA
© February/ 2017, Electro Therapy Devices Berlin Germany
© February/ 2017, ETDFL Berlin Germany