




## Culique Infrared Facial Mask

Advancing Skin Health and Beauty with Innovative Infrared Technology

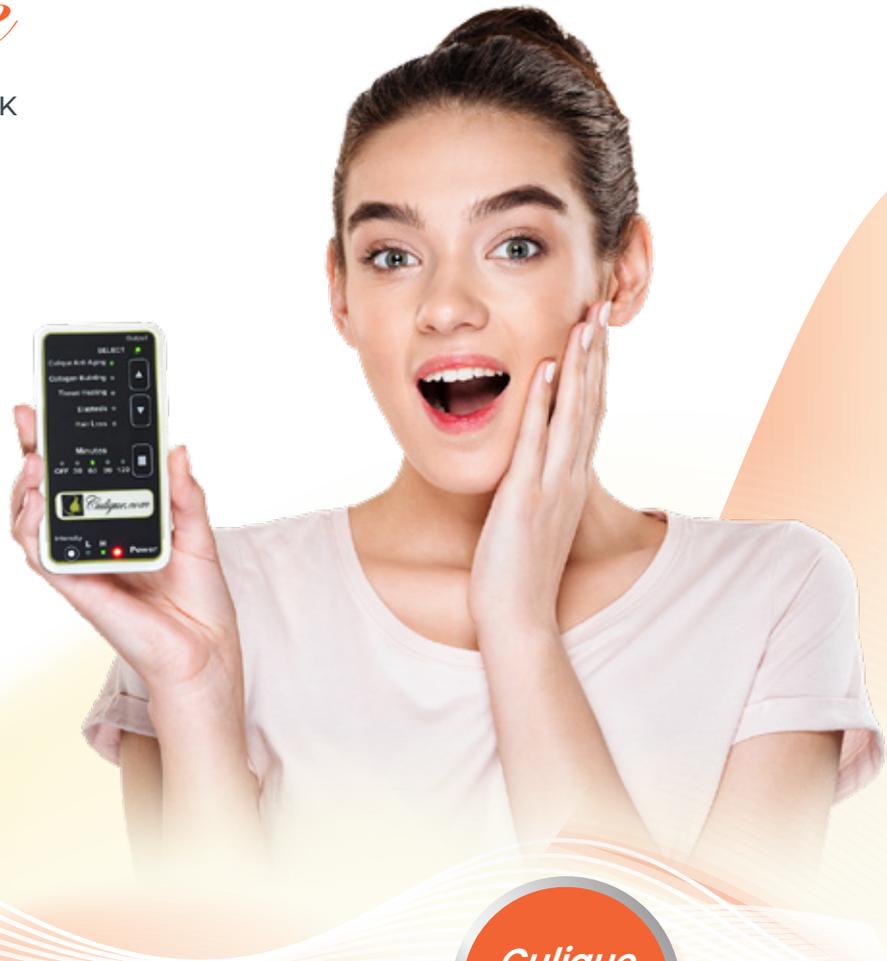

*Culique Mask*

*Culique Controller*

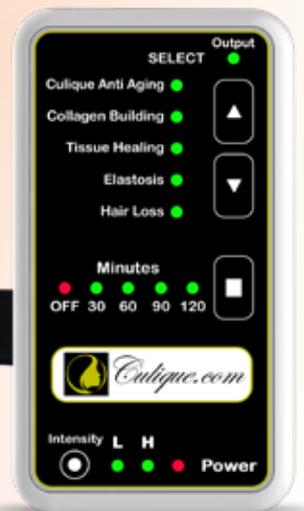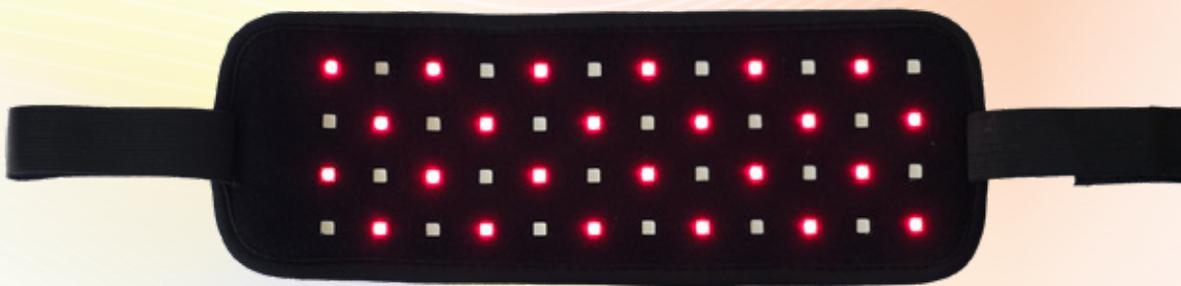


**INSTRUCTION BOOK** ENGLISH




# Index

<b>About</b>	
The Culique Skin Treatment Device .....	03
<b>What is</b>	
Infrared Light Therapy .....	04
<b>Product Description</b>	
Benefits; Components .....	05
<b>Technical Description</b>	
The Culique Infrared Skin Mask .....	06
<b>How To Connect</b>	
The Culique Infrared Skin Mask .....	07
<b>How to Wear</b>	
The Culique Infrared Skin Mask .....	07
<b>User Techniques</b>	
Wearing the Culique Infrared Skin Mask .....	08
<b>Suggested Timing</b>	
For Treatments .....	08
<b>Suggested Intensity Setting</b>	
For Treatments .....	08
<b>Controller Functions</b>	
The Culique Controller .....	09
<b>Connecting 2 People</b>	
to the Culique Controller .....	10
<b>Advanced machine connections</b>	
Connecting 4 people to the RDPV4 .....	11
<b>How Does It Work?</b>	
Skin Concern .....	12
<b>Anti Aging</b> .....	13
<b>Collagen Building</b> .....	14
<b>Tissue Healing</b> .....	15
<b>Elastosis</b> .....	16
<b>Hair Loss</b> .....	17
<b>FAQ</b>	
Frequently Asked Questions .....	18
<b>Warranty &amp; Disclaimer</b> .....	20
<b>ABOUT ETDFL</b> .....	21

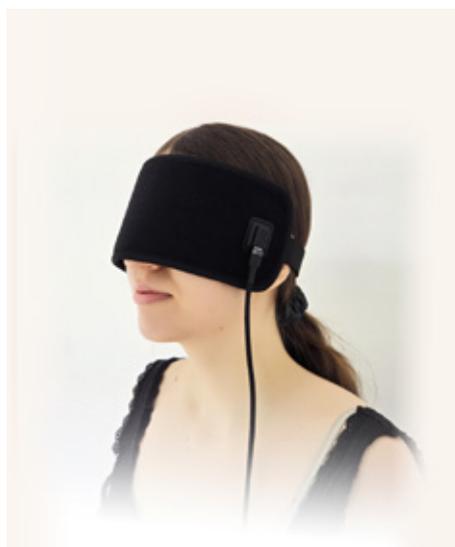

# About The Culique Skin Treatment Device

Dr. Royal Rife is a monumental figure in medical history, celebrated for his pioneering inventions and revolutionary approaches to healing. His groundbreaking work in harnessing the power of light and sound frequencies has left an enduring impact on medical treatment, inspiring technologies that bridge the gap between innovation and wellness.

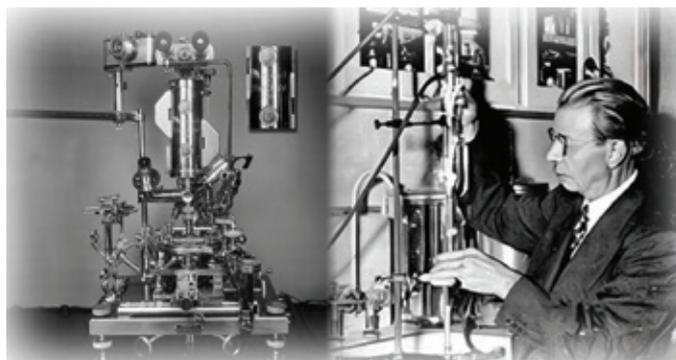

At the core of Dr. Rife's achievements is the invention of the Rife Machine, which adeptly produced sound frequencies and transmitted them to light sources. This innovative fusion of light and sound formed the foundation of Dr. Rife's explorations into therapeutic modalities, revealing the healing potential of frequency therapy.

Dr. Rife's relentless dedication to advancing medical science involved meticulous engineering and experimentation. He sought to unlock the therapeutic benefits of light therapy by integrating it with sound frequencies, pioneering a new era of treatment for various medical conditions. His groundbreaking experiments laid the groundwork for advanced technologies that continue to influence modern medicine.

Among Dr. Rife's notable inventions was the BeamRay, an early prototype that effectively combined light and frequency technologies. This revolutionary device not only showcased the potential of light therapy but also served as a blueprint for the development of contemporary medical devices.

Today, the Culique Infrared Skin Treatment Mask embodies Dr. Rife's principles by seamlessly integrating these foundational technologies. Although Dr. Rife did not invent the Culique Mask, its design utilizes his groundbreaking methods of converting sound waves into light waves, delivering targeted skincare solutions. By employing specific frequencies of light, the Culique Mask promotes cellular rejuvenation, enhances skin texture, and supports overall skin health, echoing Dr. Rife's vision of healing through frequency.

The Culique Mask represents the culmination of decades of research and innovation inspired by Dr. Rife's legacy. It stands as a testament to the transformative power of combining light and sound in therapeutic applications, reminding us of the endless possibilities that await in the realm of skincare and medical science.


Through the integration of Dr. Rife's pioneering technologies, the Culique Infrared Skin Treatment Mask offers a modern approach to skincare, drawing from a rich history of innovation to provide effective, rejuvenating treatments. In honoring Dr. Rife's legacy, the Culique Mask continues to push the boundaries of what's possible in skincare, reflecting a commitment to wellness rooted in scientific advancement.

# WHAT IS INFRARED LIGHT THERAPY?

The Culique Mask employs focused, deep-penetrating wavelengths across five distinct frequency groups, generating over 50 unique frequencies in kilohertz (kHz). These kHz sound frequencies are converted into light frequencies through the infrared LED array. This advanced therapy is designed to rejuvenate and enhance skin health by promoting cellular repair and improving circulation, effectively mimicking the benefits of natural sunlight while eliminating harmful UV rays.

## Culique Controller

The Culique Controller is an intuitive and user-friendly device designed to effortlessly manage your Culique Infrared Facial Mask. The controller transmits frequencies across specific bandwidths with output to each visible and invisible Light on the Mask.

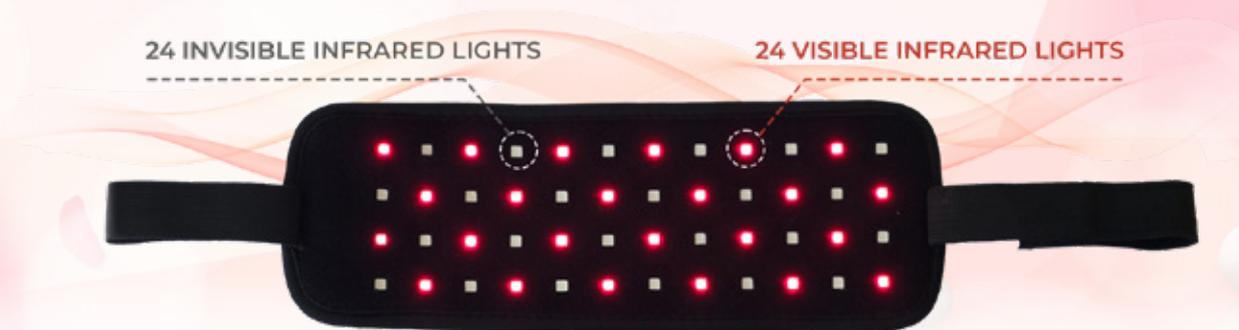

### 5 Treatment Modes

Anti-aging, collagen building, skin tissue healing, Elastosis and Hair Loss.

### Hands-Free

Easy to use, comfortable to wear

### FDA Cleared

This device is FDA-cleared, offering professional-grade infrared light therapy right in the comfort of your own home

### Comes with Adjustable Strap

The Culique Mask includes an adjustable velcro strap, providing a universal fit that ensures the device sits comfortably and securely against the skin for all users

## Culique Infrared Facial Mask

*Experience professional-grade infrared light therapy in the comfort of your home for real, visible results*



*The Culique Controller*

### Automatic Timer

Automatic timer permits 30,60,90 and 120 minutes

### Adjustable Intensity

Adjustable Intensity  
Switch from low intensity to High intensity with the press of a button

### 48 Programmable Infrared LEDs

Each of the 48 LEDs on the Culique Mask operates on a specialized frequency program tailored to address specific skin conditions and facilitate effective treatments for those issues.

## PRODUCT DESCRIPTION

This powerful FDA-cleared light therapy device designed for the full face features a combination of 24 Visible Infrared LED Lights emitting diodes in Visible mode and 24 Invisible Infrared Lights emitting diodes in Invisible mode that work together to help boost collagen production, improve skin density, smooth wrinkles, diminish discoloration, and clear acne for a clear, younger-looking complexion. For all skin types.



## Benefits

- **Smooths lines & wrinkles**

Visible Infrared RED Light penetrates deep within skin's dermis to boost collagen production and smooth fine lines & wrinkles. 97% showed an overall improvement in fine lines, wrinkles, and skin tone after 10 weeks.

- **Clears complexion**

Invisible Infrared Light destroys acne-causing bacteria to help clear mild to moderate acne. Visible Infrared RED Light helps even skin tone & texture to prevent flare-ups. 97% showed overall improvement in acne after 2 weeks.

- **Firms skin & increases skin density**

Visible Infrared RED Light penetrates into skin to support natural collagen and visibly firms skin.

- **Treats signs of sun damage**

Visible Infrared RED Light works to reduce sun spots, age spots, and discoloration.

- **Hair Loss**

Using the special Hair Loss program, the mask can also be used for hair regrowth.

## Components

### Key Components

- 24 Visible Infrared RED Lights emitting diodes penetrate deep within skin's dermis to boost collagen production. Each Led is powered by a specific Khz frequency bandwidth according to the selected program
- 24 Invisible Infrared Lights emitting diodes destroy P. acnes bacteria to reduce redness, inflammation, and prevent future breakouts. Each Led is powered by a specific Khz frequency bandwidth according to the selected program

# Technical Description

## The Culique Infrared Skin Mask

### Frequency range

**Tested** 0:01 Hz – 50 MHz

### Frequency Wave and Sideband

Transmits all ETDFL frequencies Using Multiple Infrared Light Spectrums when connected to RDPV4 machine. All Dr Rife MOR Carrier Wave and RF frequencies are reproduced through the Positive and Grounded areas of the Mask with +99%† efficiency.

### Brightness Settings when Using the Culique Mask with the Culique Controller

Set the Culique Controller to the LOW illumination level to start the treatment. After a while, you can switch to HIGH for more effective Skin treatment protocols. Be sure to keep your Skin open while using the Culique Mask, as the Multiple Mode Infrared Plus frequencies need to reach both Skin. If the light feels too intense, adjust the controller back to the LOW setting.

### Brightness Settings when Using the Culique Mask with the RDPV4 Controller

- Configure the RDPV4 to a voltage level of 12.5 V DC if the illumination level is excessively intense. 12.5 V DC is the lowest illumination level. The middle setting is 14.5 V, while the default maximum level is 17.5 V, which is recommended for Skin treatments.
- Be sure to keep your Skin open while using the Culique Mask, as the Multiple Mode Infrared Plus frequencies need to reach both Skin. If the light feels too intense, adjust the controller back to the LOW setting.

### Clinic Use, and Cleaning

If using the Mask for multiple people, or for a clinic, wipe down the surface area after use with a cloth. Do not use alcohol directly on the surface. You cannot wash the Mask in a washing machine, as this will damage the electronics.



## How To Connect

### The Controller to the Culique Infrared Skin Mask



Connect the Controller to power using the power cable (1), connect the Controller to the Culique Infrared Skin Mask using the USB-C cable (2)


## How To Wear

### The Culique Infrared Skin Mask




24 INVISIBLE INFRARED LIGHTS

24 VISIBLE INFRARED LIGHTS

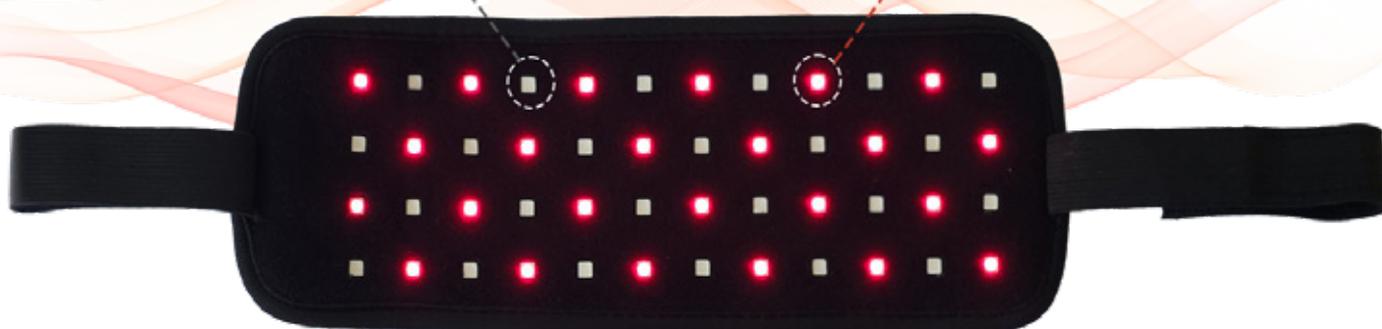

## User Techniques

### Wearing the Culique Infrared Skin Mask

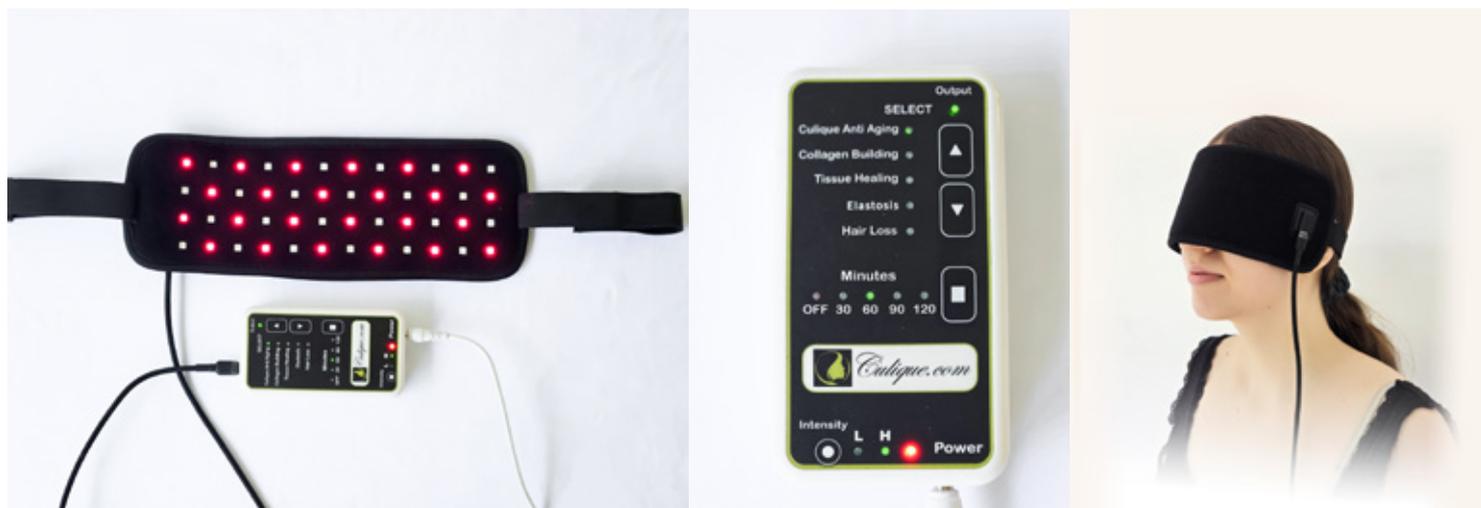

- Connect the Controller to power and the Culique Infrared Skin Mask.
- Check the LED power indicator on the Controller
- Check the infrared light on the mask.
- Select the desired mode and timing.
- Wear the Culique Infrared Skin Mask (see above).
- Remember: The 24 Lights facing the Skin are INVISIBLE. You will not see them, but they are on, and transmitting light at a frequency invisible to the human Skin.

## Suggested Timing

### For Treatments

- Choice: 30, 60, 90, 120 minutes
- 30 minutes, minimum treatment.
- 120 minutes maximum treatment time.


## Suggested Intensity Setting

### For Treatments

- Start with low intensity settings.
- High intensity settings for more serious Skin disorders, however revert back to low intensity settings in the event of discomfort
- Skin must be open during treatment as the treatment is focused on actual Skin treatment

# Controller Functions

## The Culique Controller

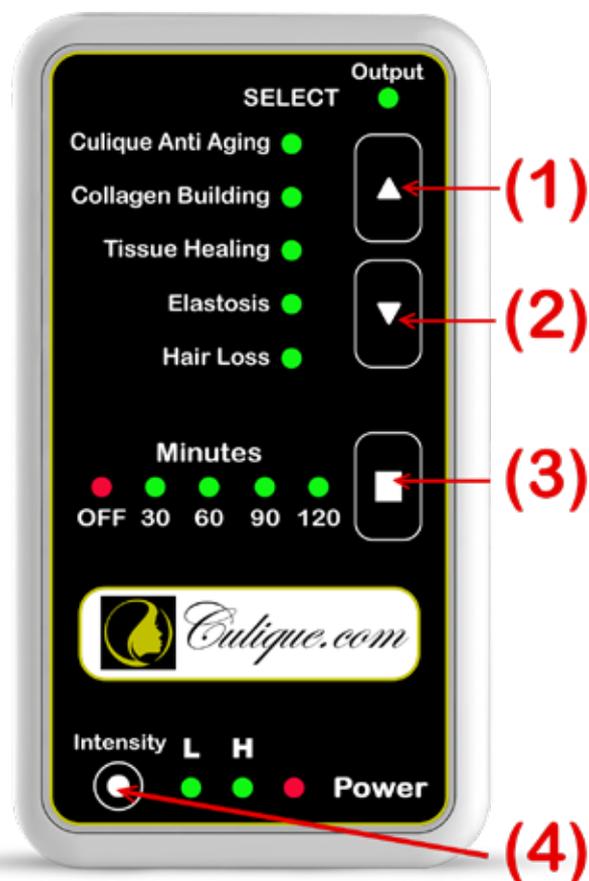

The Culique Controller

### (1) BUTTON 1: UP

Use UP arrow to move UP Skin disorder Categories 1- 5

### (2) BUTTON 2: DOWN

Use DOWN arrow to move down Skin disorder Categories 1- 5

### (3) BUTTON 3: TIME ADJUSTMENT

USE To Switch treatment time from OFF = No Output to Mask 30 Min, 60 Min, 90 Min, 120 Min

### (4) BUTTON 4: POWER INTENSITY

L = Low Power Infrared output  
H = High Power Infrared output

#### Output

**Output LED ON:** The Controller is sending frequency signal output to the Culique Infrared Skin Mat

**Power LED ON:** The Controller is connected to Power Source (Power Adapter)

### \* CONTROLLER PORT SIDE VIEW



Controller Port



The Culique Mask

- (1) USB Port for one Culique Infrared Skin Mask
- (2) USB Port for a Second Culique Infrared Skin Mask
- (3) POWER PORT (Power Input)





USB - C

# Connecting 2 People to the Retyne controller




## Using 2 Culique Infrared Mask for 2 People with the Controller

- Connect the controller to power.
- Connect the controller to 2 Culique Infrared Masks using the 2 USB-C output ports.
- Check the LED power indicator on the controller and the infrared lights on both masks.
- Select the desired mode and timing (both Culique Infrared Masks will use the same mode and timing settings on the controller).
- Wear the Culique Infrared Skin Mask.

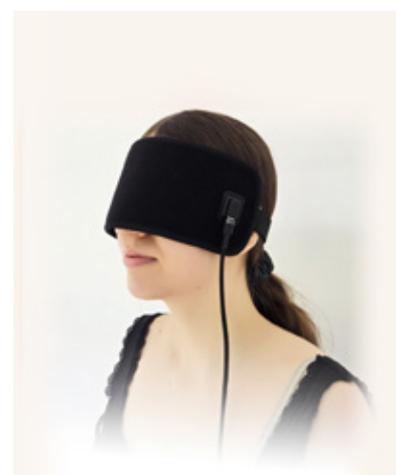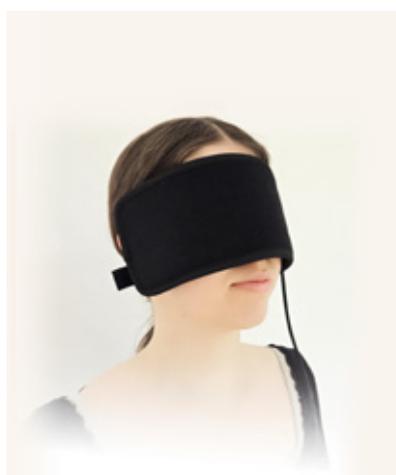

# Advanced machine connections

## Connecting 4 people to the (optional) RDPV4

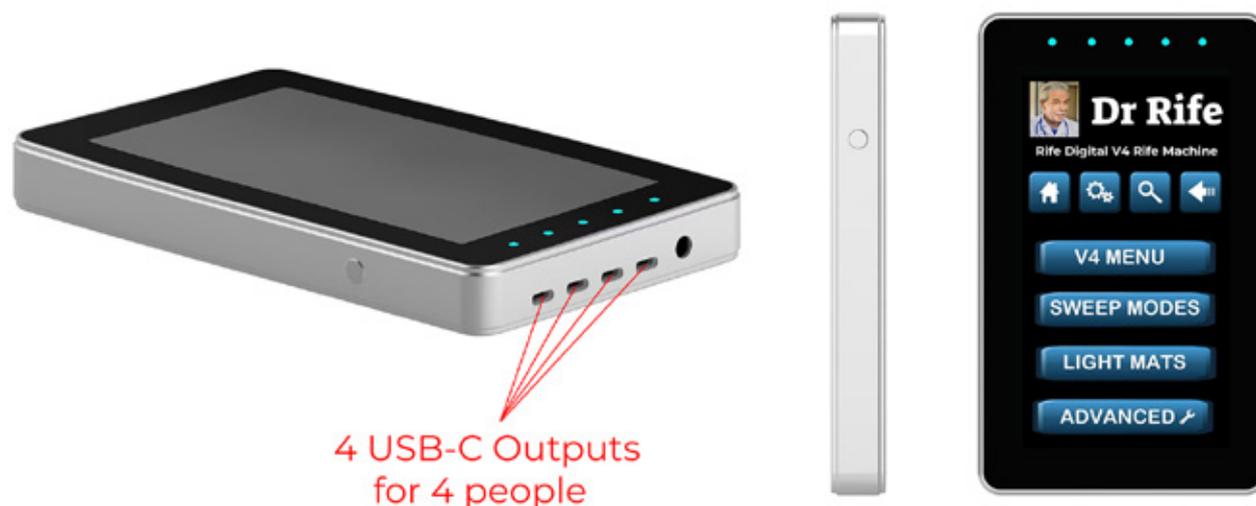

*Rife Digital Professional V4 (RDPV4)*

The Culique controller provided offers ample General Skin Treatment frequencies for various Skin disorders. However, owners of the Rife Machine RDPV4 can access more specific groups by selecting the "Advanced ETDFL Program Menu"

### 1. Accessing Specific Groups:

Navigate to the MENU setting or use the search bar to find the specific "Group" number corresponding to the Skin disorder you wish to treat.

### 2. Connecting the Devices:

Insert the Culique USBC cable into the top of the RDPV4 machine. Ensure the Culique Infrared Skin Mask is securely connected to the machine.

### 3. Using the RDPV4 CULIQUE MENU:

Select Culique Light Mask Menu on the RDPV4 and follow instructions on screen

### 4. Using RDPV4 GROUP MENU:

Select any Program number from the ETDFL Group menu, and sweep that Group while the Culique Mask is connected.

### 5. Adjusting Output Settings:

Set your RDPV4 to the Highest voltage setting when it's connected to the Culique Infrared Skin Mask to optimize treatment effectiveness.

These steps enable advanced connections between the Culique Infrared Skin Mask and Rife Machine RDPV4, facilitating precise and enhanced treatment for various Skin conditions.

# HOW DOES IT WORK?

## 1. Target Aging

The 10 specific wavelengths of light used in these masks target various skin concerns, promoting collagen production and enhancing skin elasticity, which helps to smooth out fine lines and wrinkles. The therapy also stimulates cellular turnover, resulting in a brighter and more even skin tone. Regular use can lead to a more radiant, rejuvenated complexion, making it a convenient and effective tool in the fight against aging.

## 2. Collagen Building


The targeted infrared light therapy stimulates collagen production, enhancing skin elasticity and reducing the appearance of fine lines and wrinkles. By penetrating the skin at various depths, Culique mask promotes cellular regeneration and improve overall skin tone and texture. Regular use can result in a more youthful, radiant complexion, making it a convenient and effective addition to any skincare routine.

## 3. Skin Tissue Healing and Acne

Using an LED mask for skin tissue healing offers a range of transformative benefits, particularly for those recovering from injuries, acne, or other skin concerns. The therapeutic infrared light penetrates deep into the skin, promoting cellular regeneration and accelerating the healing process. By enhancing blood circulation, LED masks facilitate the delivery of oxygen and nutrients to damaged tissues, reducing inflammation and speeding up recovery time. Additionally, the anti-inflammatory properties of certain wavelengths can help soothe irritated skin, making it an excellent option for post-procedure care.

## Pre-Programmed Treatment Time

This FDA cleared LED treatment device utilizes a combination of Visible & Invisible Infrared Light to clear and prevent breakouts and treat aging skin from 30 minutes a day.

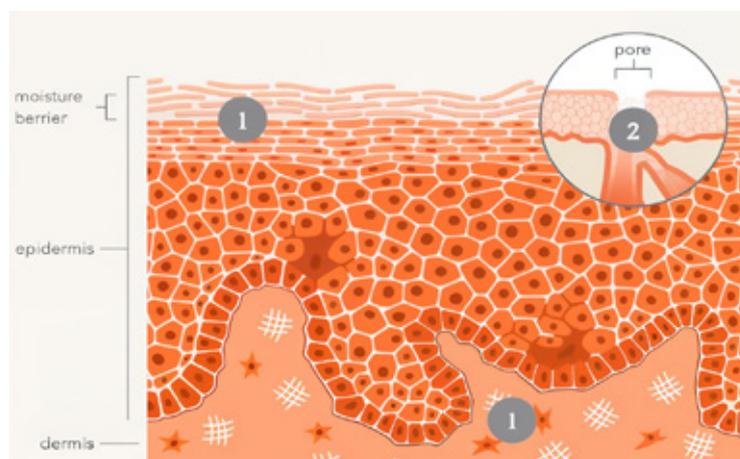

## 4. Elastosis

Elastosis, characterized by the degradation of elastic fibers in the skin, can lead to sagging and a loss of firmness. The Culique mask harnesses specific wavelengths of light that stimulate collagen and elastin production, promoting skin resilience and restoring a more youthful appearance. Additionally, The Culique mask improves blood circulation, aiding in the delivery of essential nutrients while facilitating detoxification. The result is not only improved skin elasticity but also a smoother, more even skin tone.

## 5. Hair Loss

Using the Culique LED mask for hair loss offers a promising solution for those seeking to revitalize thinning hair. Readjust the mask to be worn over the scalp, particularly target the thinning areas, the Culique mask utilizes specific wavelengths of light that stimulate hair follicles and enhance blood circulation. This process promotes cellular activity and encourages the growth of new hair while also helping to strengthen existing strands. The gentle warmth from the LED lights can improve nutrient absorption, further supporting healthy hair growth. Regular use can lead to noticeable improvements in hair thickness and overall scalp health.

# ANTI AGING

## Understanding the Culique Mask and Its Anti-Aging Benefits

The Culique Mask is a cutting-edge skincare tool designed for anti-aging treatments, featuring an "Anti Aging" setting (#1) on its controller. This setting transmits ten specific frequencies (22.5, 0.2, 132.8, 158.6, 429.3, 533.9, 678.21, 834.63, 888.76, 942.32 kHz) from Group 877, aimed at rejuvenating skin by targeting various aspects of aging.

A standout feature is its adjustable design, allowing users to position the mask for targeted treatment. By tilting it upwards, the mask addresses forehead lines, eye wrinkles, and areas around the nose and cheeks. Tipping it downwards focuses on the lower face, treating nasolabial folds and aging around the lips. This versatility ensures infrared light reaches deeper skin layers for maximum benefit.



Moreover, the design of the Culique Mask allows it to be wrapped around the neck, further extending its anti-aging capabilities. This positioning is particularly beneficial for treating areas that are often overlooked but are crucial for a youthful appearance. The infrared light can penetrate the skin on the neck, promoting collagen production and improving skin elasticity, which can help diminish the appearance of sagging skin and fine lines around the mouth and jawline.

The combination of precise frequency transmission and adjustable positioning makes the Culique Mask a versatile tool in any skincare regimen. Whether focusing on the forehead, chin, or neck, users can tailor their treatments to address their unique concerns, ensuring that each session is as effective as possible.

Ultimately, the Culique Mask represents a holistic approach to anti-aging skincare, leveraging advanced technology and user-friendly design to deliver impressive results. With its ability to transmit targeted infrared light through multiple frequencies, it not only rejuvenates the skin but also empowers users to take control of their skincare journey.



*Culique Anti Aging*

# COLLAGEN BUILDING

## Collagen Building Benefits of the Culique Mask

The Culique Mask is an advanced skincare device that uses targeted frequencies to promote collagen production and improve skin health. Setting #2 on the controller, labeled "Collagen Building," emits ten specific frequencies (0.07, 0.32, 18.23, 37.22, 57.5, 62.73, 86.32, 89.7, 96.5, and 236.09 kHz) selected to stimulate collagen synthesis. This setting is ideal for users looking to enhance skin elasticity and reduce signs of aging.

The mask's adjustable design allows precise targeting of collagen-building treatments. By positioning the mask upward, users can treat forehead lines, eye wrinkles, and the areas around the nose and cheeks. Tilting the mask downward focuses on the mouth area, addressing nasolabial folds and lines around the lips. This flexibility ensures that infrared light reaches the deeper skin layers, promoting a smoother, more youthful appearance.



*Collagen Building*



The controller converts these frequencies into visible and invisible infrared light, which penetrates the skin, triggering photobiomodulation (PBM). In PBM, light absorbed by mitochondria increases ATP production, energizing cells and stimulating fibroblast activity. This leads to enhanced collagen production and skin rejuvenation. The mask can also be wrapped around the neck, treating sagging skin and fine lines in this often-overlooked area to support a youthful contour and smooth texture.

Additionally, infrared therapy offers anti-inflammatory benefits, reducing chronic inflammation and creating a favorable environment for collagen synthesis. By promoting the remodeling of collagen fibers, the mask improves skin texture and firmness.

In summary, the Culique Mask's precise frequency transmission and versatile design make it a powerful tool in anti-aging skincare, enhancing collagen production, reducing inflammation, and supporting overall skin resilience for a more youthful complexion.

# TISSUE HEALING

## Tissue Healing Benefits of the Culique Mask

The Culique Mask's "Tissue Healing" setting (#3) offers an advanced solution for skin and scar tissue recovery. It transmits ten specific frequencies (0.12, 0.37, 7.25, 45.67, 45.75, 57.5, 92.5, 93.5, 96.5, and 325 kHz) designed to promote healing across various tissue types. This versatile mask can be used on multiple body areas, including the face, neck, chest, and limbs, making it ideal for treating post-surgical scars, injuries, acne scars, and other skin concerns.

Through photobiomodulation (PBM), the mask enhances mitochondrial function, increasing ATP production, which fuels the body's natural healing processes. The controller converts these frequencies into visible and invisible infrared light, which penetrates the skin at varying depths to support cellular repair and tissue regeneration.



*Tissue Healing*



The "Tissue Healing" setting reduces inflammation, which can otherwise hinder recovery and prolong discomfort. By mitigating inflammation, the mask fosters a conducive healing environment, enabling faster and more effective tissue repair.

Photobiomodulation occurs when the infrared light is absorbed by chromophores, particularly in mitochondria. This stimulation increases ATP, which boosts fibroblast activity—essential for collagen and elastin

Additionally, the mask promotes vasodilation, enhancing blood circulation to increase oxygen and nutrient delivery. This accelerates healing by facilitating the removal of metabolic waste and inflammatory mediators. Its anti-inflammatory effects also reduce the risk of hypertrophic or keloid scars by creating a favorable healing environment.


The Culique Mask aids in all wound healing stages: it reduces inflammation in the initial phases, supports cell proliferation and angiogenesis in the proliferation phase, and assists in collagen remodeling during the final stage. This multifaceted approach makes the mask an invaluable tool for comprehensive skin healing and recovery.

# ELASTOSIS

## Elastosis Treatment with the Culique Mask

The Culique Mask's "Elastosis" setting (#4) is specifically tailored to combat skin sagging and improve elasticity. This setting transmits ten targeted frequencies (0.05, 0.12, 0.87, 5, 27.5, 62.5, 193, 322.53, 475.17, and 527 kHz) designed to address the underlying causes of elastosis, restoring a youthful appearance by enhancing collagen and elastin production.

The mask utilizes photobiomodulation (PBM), where infrared light penetrates the skin and interacts with mitochondria, increasing ATP production. This energy boost fuels cellular processes, promoting the synthesis of collagen and elastin—two essential proteins for skin structure and elasticity. Enhanced fibroblast activity leads to the production of the extracellular matrix, reinforcing the skin's firmness and resilience.

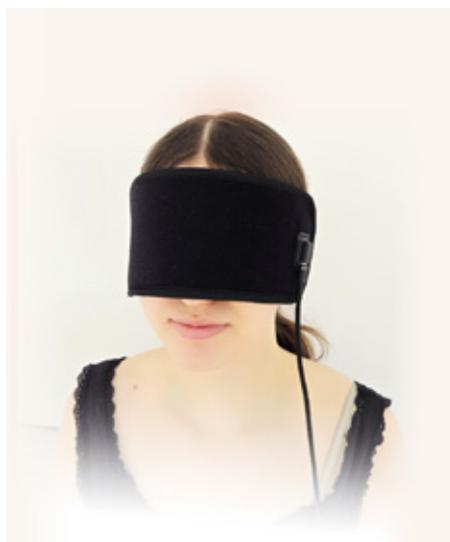

*Elastosis*

Additionally, the mask improves blood circulation, delivering oxygen and nutrients to the skin and aiding in the removal of waste. This optimal environment supports cellular repair, enhances skin texture, and facilitates recovery from sagging.

By stimulating fibroblast activity, the Culique Mask encourages the synthesis and remodeling of collagen and elastin, which are crucial for maintaining firm, elastic skin. Infrared light also reduces inflammation, preserving the skin's structural integrity and counteracting the collagen degradation associated with aging.

Consistent use of the "Elastosis" setting can lead to noticeable long-term improvements in skin texture, firmness, and elasticity, making it a valuable tool for addressing signs of aging and promoting a youthful, resilient complexion.

Regular use of the Culique Mask on the "Elastosis" setting also provides cumulative benefits that enhance skin resilience over time. As collagen and elastin synthesis are gradually stimulated, users can experience a progressive firming of the skin, with visible reductions in sagging and a smoother, more toned complexion.


# HAIR LOSS

## The Culique Mask for Hair Loss Treatment

The Culique Mask's "Hair Loss" program (#5) offers a comprehensive solution for addressing hair thinning and conditions like alopecia areata. Through a set of ten specific frequencies (0.06, 5.07, 95, 127.67, 275.07, 455.82, 515.16, 684.81, 712.23, and 993.41 kHz), the mask stimulates hair follicles and promotes scalp health, leading to improved hair density and vitality.

At the core of its effectiveness is photobiomodulation (PBM). Infrared light penetrates the scalp, where it is absorbed by mitochondria within hair follicle cells, boosting ATP production. This energy enhancement drives cellular functions, increasing hair follicle activity and encouraging growth. Dermal papilla cells, which regulate the hair cycle, are particularly responsive, transitioning from the resting (telogen) phase to the growth (anagen) phase more readily, resulting in thicker hair and reduced shedding.

The mask also improves blood circulation in the scalp, delivering essential nutrients to hair roots, which strengthens dormant or weakened follicles. Additionally, the infrared light's anti-inflammatory properties help to mitigate chronic scalp inflammation, creating a more favorable environment for hair growth. This reduction in inflammation helps alleviate conditions like androgenetic alopecia, where inflammation can disrupt the normal hair cycle.

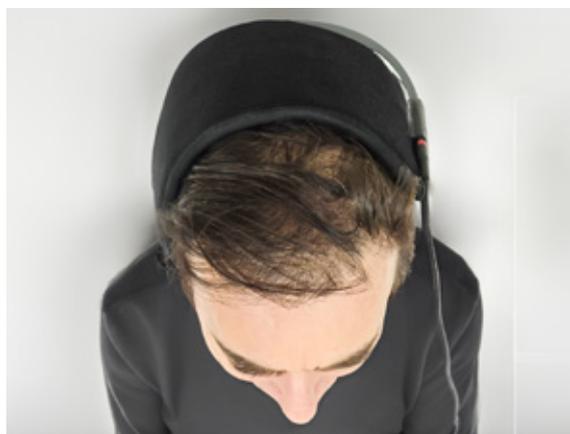

Furthermore, the Culique Mask supports scalp health by promoting collagen and elastin production. These proteins contribute to the scalp's structural integrity, providing a healthy foundation for hair follicles to thrive. The infrared light also activates stem cells within the follicles, which are essential for regenerating damaged hair structures and encouraging new growth.

By leveraging the mask's advanced technology, users can address hair loss holistically, benefiting from improved follicle function, reduced inflammation, enhanced circulation, and a fortified scalp, all of which contribute to revitalized, fuller hair.



*Hair Loss*

## FAQ

### Frequently Asked Questions

#### ■ **Eyes open or closed when using the Culique IR mat on the skin?**

*For optimal results during facial skin treatments, it is essential to keep your eyes closed. This not only protects your eyes but also allows the infrared lights to effectively treat the eyelids and surrounding areas. If you find the light intensity uncomfortable, consider using a blackout mask to cover your eyes while using the Culique mat, or adjust the controller to the "Low" setting for a more comfortable experience.*

#### ■ **What setting on the Rife Professional V4 do I use with the Culique Infrared Mask?**

*We suggest you look up the exact frequency of the Skin-related disorder you want to treat. Often, you will find the disorder listed within the Presets of the machine. If not, use the ETDFL book (Free Download at [ETDFL.COM](http://ETDFL.COM)) to look up the frequencies for Skin-related treatments.*

#### ■ **Where can I buy the most up-to-date Rife Machine to use with my Infrared Mask accessory?**

*You can purchase the Rife Professional V4 directly from the [ETDFL.com](http://ETDFL.com) Store shopping cart. Although availability may not be until 2025*

#### ■ **I feel sensitivity when using the Culique Infrared Mask on my Skin:**

*We suggest you start with the "Low" setting on the controller and ensure that the mask fits comfortably without pressing too hard against your Skin. If any sensitivity is felt, ensure the controller is on the LOW Intensity setting. The mask should be a looser fit around the Skin, the mask should not be tight against the face.*


#### ■ **The mask doesn't fit my face properly, and there's a gap around my Skin. Help!**

*Adjust the Velcro straps of the Culique Infrared Mask to ensure a snug and comfortable fit around your Skin. Ensure that the mask makes contact with the skin around your Skin but isn't too tight, allowing the frequencies to be delivered effectively.*

#### ■ **Can I use these accessories for my Dog or Cat, Animals?**

*No, the Culique Infrared Mask is designed specifically for human Skin treatment and should not be used on animals.*

- What is the purpose of the specific array of 24 visible lights and 24 invisible lights, and why are they placed in this order?

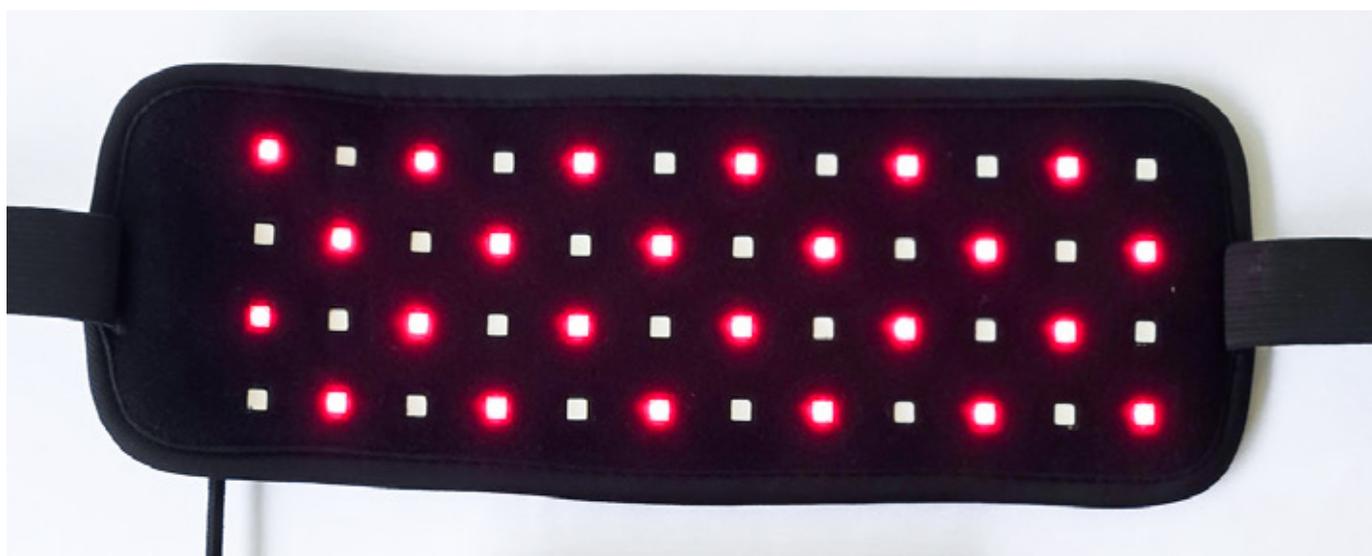

*The 24 invisible infrared lights are engineered to emit specific frequencies across designated bandwidths, directly impacting the Skin. The invisible infrared lights are designed to produce highly harmonious light bandwidths, optimizing compatibility for close Skin contact. These invisible wavelengths generate minimal light stress on the Skin while ensuring effective healing transmissions. Utilizing light as the medium of transmission, the invisible infrared lights convey targeted frequencies to various structures of the Skin*

*The 24 visible red LEDs serve a crucial role in skin treatment by emitting light that is readily absorbed by the skin's surface. This visible light stimulates cellular activity and promotes various beneficial processes, such as collagen synthesis and increased circulation. By enhancing blood flow, the red LEDs help deliver essential nutrients and oxygen to skin cells, facilitating the repair and regeneration of tissues. This process not only improves skin texture and tone but also aids in the reduction of fine lines and wrinkles, contributing to a more youthful appearance.*

*The arrangement of the 24 visible red LEDs in conjunction with the 24 invisible infrared LEDs is deliberate and strategic. This configuration ensures a synergistic effect, allowing both types of light to work together to maximize the therapeutic benefits. The visible red light targets the epidermal layers, stimulating surface-level healing and rejuvenation, while the invisible infrared light penetrates deeper into the dermis, addressing underlying issues such as collagen degradation and inflammation. This harmonious interaction between the visible and invisible light sources creates a comprehensive treatment approach that addresses multiple layers of skin concerns effectively.*


# Warranty & Disclaimer

## Warranty

- DOA Warranty: If your Culique Infrared Mask arrives with any damage caused by transport/courier, it will be replaced immediately.
- 3-Month Replacement Warranty: If your Culique Infrared Mask breaks under normal use, we will replace the item within 90 days. The cost of shipping will be charged for the return of a new item to your address.
- The Culique mat is equipped with 44 Infrared and Invisible Infrared LEDs, 24 of which are designed to be invisible to the naked eye. These invisible LEDs can only be detected with a phone camera, which confirms the mat's functionality and effectiveness during use, even if the 24 invisible LEDs aren't visible. Regarding warranty conditions, a single LED malfunction, such as a dead pixel, is often categorized as "normal wear and tear" rather than a defect. Typically, an LED array is considered defective only if there are multiple dead 'LED' pixels in a specific area or if they form a cluster.
- Due to sanitary reasons, the Culique Infrared Mask will not be repaired by our service center. We will always issue a new product under warranty.

## Disclaimer

- The E.T.D.F.L. (ETDFL.com) Labs holds no guarantee that the Culique Infrared Mask will cure any specific Skin disorders. The Infrared Mask technology has been tested for over 3 years in the EU by researchers who have had successful results using the mask in conjunction with the RDPV4 and the Culique controller.
- The Culique Infrared Mask Technology is grounded in extensive research, however, the collaborating researchers do not guarantee that their findings or technologies will result in a cure for Skin disorders.
- The Culique Infrared Mask is a unique product available for individual or clinical purchase.

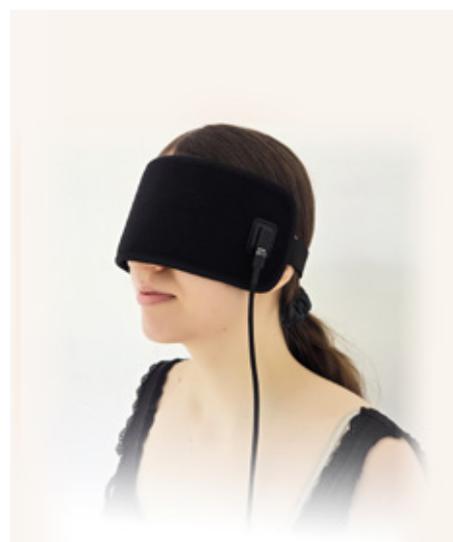



ETDFL.com



E.T.D.F.L is an ongoing project by 12 Clinics across the planet who use the Quantum SCIO Bioresonance machine to record Disease frequencies of actual clients, and then produce the 10 most accurate frequencies which are deemed suitable for treatment with the Royal Rife Machine. This project began in 2006 when the original CAFL frequency list was discontinued. ETDFL became the new list, now used by more than 100,000 Royal Rife Machine users across the planet. The accuracy of the frequencies is seen in the positive results of the clients who use this Frequency List in clinics and homes in every country.

ETDFL Berlin Germany, develop the latest accessories and equipment for Electro Frequency Therapy. The Culique Mask is a unique skin treatment product developed in conjunction with Culique Labs, Berlin, Germany and ETDFL Labs Germany.

**ETDFL** is updated once a year. Latest PDF versions available at: [ETDFL.com](http://ETDFL.com)

For additional information on skin diseases, skin treatments, and the Culique Mask, including skin concern groups not covered in this book, please visit the website

**Culique.com**

PDF versions of ETDFL at **ETDFL.com** allow you to search the ETDFL Database for Skin Concern frequencies which may be more up to date than the List provided within this book.

The Culique.com Skin Treatment Mask is licensed to use all ETDFL Royal Rife Frequency Sets © ETDFL Berlin Germany



The Culique Store prioritizes specialized skincare solutions by utilizing advanced infrared technology, designed to address a wide range of skin concerns, from acne to signs of aging

## Contact Us

 [www.culique.com](http://www.culique.com)

 206-708-6408

 [etdfl.web@gmail.com](mailto:etdfl.web@gmail.com)

 **Culique Store USA**  
1712 Pioneer Avenue,  
Cheyenne, WY 82001, USA

 **Culique Labs**  
Internationales Handelszentrum,  
Berlin Friedrichstrasse 95,  
Germany