




# Beamer PEMF

Pulsed Electro Frequency Magnetic Field Mat

[www.beamer.tech](http://www.beamer.tech)



*Beamer PEMF Mat*

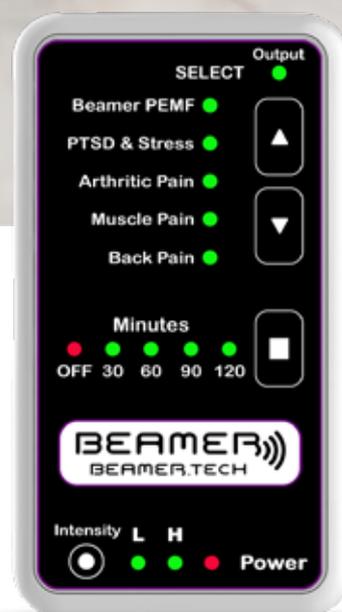

*Beamer Controller*

# INSTRUCTION BOOK

ENGLISH



# Index

<b>About</b>		
Legacy of Beamer PEMF Technology .....	03	
<b>What is Beamer PEMF Mat?</b> .....	04	
<b>Product Description</b> .....	05,06	
<b>How PEMF Works?</b> .....	07	
<b>Technical Description</b>		
The Beamer PEMF Therapy Device .....	08	
<b>How To Connect</b> .....	09	
<b>Controller Functions</b>		
The Beamer Controller .....	10,11	
<b>Timing Selection</b>		
For Treatments .....	11	
<b>Suggested Intensity Setting</b>		
For Treatments .....	11	
<b>Advanced Machine</b>		
Connections .....	12	
<b>BEAMER PEMF</b>		
PROGRAM #1 .....	13	
<b>PTSD &amp; Stress</b>		
PROGRAM #2 .....	14	
<b>Arthritic Pain</b>		
PROGRAM #3 .....	15	
<b>Muscle Pain</b>		
PROGRAM #4 .....	16	
<b>Back Pain</b>		
PROGRAM #5 .....	17	
<b>Advanced Treatment</b> .....	18	
<b>CONDITIONS</b> .....	19,20	
<b>Integumentary System</b> .....	21	
<b>Skeletal System</b> .....	22	
<b>Muscular System</b> .....	23	
<b>Lymphatic System</b> .....	24	
<b>Respiratory System</b> .....	25	
<b>Digestive System</b> .....	26	
<b>Nervous System</b> .....	27	
<b>Endocrine System</b> .....	28	
<b>Circulatory System</b> .....	29	
<b>Urinary System</b> .....	30	
<b>Reproductive Systems</b> .....	31	
<b>Other Categories</b> .....	32	
<b>FAQ</b>		
Frequently Asked Questions .....	33,34	
<b>Warranty &amp; Disclaimer</b> .....	34	
<b>ABOUT ETDFL</b> .....	35	

## About

# Legacy of Beamer PEMF Technology

Dr. Royal Rife was a pioneering figure in frequency-based wellness technology, known for his groundbreaking research on bioelectrical frequencies and their impact on the body's natural healing processes. His work laid the foundation for modern frequency-based therapy, inspiring advancements in non-invasive wellness solutions.

Building on this legacy, Beamer-Tech has revolutionized frequency therapy through Pulsed Electromagnetic Field (PEMF) technology. Unlike light-based devices, Beamer-Tech delivers targeted electromagnetic pulses to stimulate cellular repair, improve circulation, and enhance overall recovery.



Dr. Rife's research emphasized the importance of precise frequency application, a principle that remains central to Beamer-Tech's PEMF therapy programs. By transmitting optimized PEMF frequencies, Beamer-Tech provides users with a powerful, drug-free approach to wellness—supporting pain relief, muscle recovery, and relaxation.

Today, the advancement of PEMF technology continues to shape the future of holistic health and wellness. Beamer-Tech honors this legacy by integrating scientifically validated frequency therapy into a modern, user-friendly PEMF system, making high-quality, non-invasive recovery solutions accessible to all.

The Beamer PEMF technology uses Pulsed frequencies, transmitted through a large mat to balance the human body organ system as a whole body healing modality. The PEMF method of treatment has been in use successfully for over 60 years, and continues to be a primary method of achieving homeostasis and a state of balance within the body.

# WHAT IS BEAMER PEMF MAT?

The Beamer technology utilizes targeted, deep-penetrating Pulsed Electromagnetic Fields (PEMF) across multiple frequency ranges, delivering optimized energy waves to the body. These carefully calibrated PEMF signals are designed to enhance cellular function, improve circulation, and support faster recovery by promoting natural self-healing processes.



Beamer PEMF Mat

The Beamer Controller is an intuitive and user-friendly device designed to seamlessly manage your Beamer PEMF Mat.



Beamer Controller

It transmits targeted pulsed electromagnetic frequencies across specific bandwidths, delivering optimized output to the Beamer Mat for enhanced therapeutic benefits.

## 5 Treatment Modes

- Beamer PEMF
- PTSD & Stress
- Arthritic Pain
- Muscle Pain
- Back Pain

## Full-Body Comfort and Versatility

The Beamer Mat is designed to provide a comfortable and adaptable therapy experience, allowing users to lie down and receive full-body PEMF treatment. Its flexible and lightweight structure ensures that it can be easily placed on a bed, recliner, or treatment table, making it ideal for various therapy settings.

## Automatic Timer

Automatic timer permits 30,60,90 and 120 minutes

# Beamer PEMF

*Pulsed Electro Frequency Magnetic Field Mat*



## Beamer PEMF Mat

Beamer Mat – A full-body PEMF mat designed to deliver targeted pulsed electromagnetic fields for enhanced recovery, relaxation, and pain relief.

## Adjustable Intensity

Adjustable Intensity  
Switch from low intensity to High intensity with the press of a button

## Optimized PEMF Frequency Programs

The Beamer Mat operates on specialized PEMF frequency programs, designed to target key areas of the body. These carefully calibrated electromagnetic pulses help enhance circulation, support cellular recovery, and promote overall well-being.

# PRODUCT DESCRIPTION

The Beamer Mat is a powerful PEMF therapy device designed to deliver targeted pulsed electromagnetic fields for optimal recovery and well-being. Using scientifically calibrated PEMF frequency programs, the Beamer Mat helps enhance circulation, support muscle recovery, and promote overall body balance.

With a versatile and user-friendly design, the Beamer Mat provides a non-invasive, drug-free approach to managing stress, pain relief, and rehabilitation. Whether used for daily wellness or targeted therapy, the Beamer Mat ensures effective and comfortable treatment sessions.

## Our Promise

### Advancing Therapeutic Innovation: Inspired by Pioneers

Dr. Royal Rife was a visionary in medical technology, renowned for his groundbreaking research in frequency-based therapy. His pioneering approach utilized specific frequency transmissions to support healing and wellness, laying the foundation for modern therapeutic advancements.

Building on these principles, Beamer-Tech continues this legacy by employing Pulsed Electromagnetic Field (PEMF) therapy—a scientifically backed technology designed to enhance circulation, support cellular regeneration, and promote overall well-being.

The Beamer-Tech Mat represents the next evolution in non-invasive therapy, using targeted PEMF frequencies to optimize recovery, pain management, and relaxation. This advancement is transforming the way individuals experience holistic wellness and therapeutic care, offering a powerful, drug-free solution for modern health needs.

## About [beamer.tech](https://beamer.tech)

### GLOBAL PRESENCE

Beamer-Tech has a strong international presence, with its PEMF therapy solutions available in key markets, including Germany, the UK, USA, Hong Kong, and Japan. Dedicated to advancing non-invasive wellness technology, Beamer-Tech provides individuals with access to cutting-edge PEMF therapy, supporting recovery, relaxation, and overall well-being.

Beamer-Tech remains committed to continuous innovation, delivering state-of-the-art PEMF technology to individuals seeking a drug-free, non-invasive approach to enhanced wellness and recovery.

### INNOVATION AND DEVELOPMENT

The Beamer-Tech PEMF Mat was developed through years of research and engineering excellence, incorporating the latest advancements in Pulsed Electromagnetic Field (PEMF) therapy. Designed for both home and clinical use, this innovative mat has been scientifically calibrated to optimize circulation, pain relief, and cellular regeneration, ensuring users receive a proven, effective therapy solution.

# PRODUCT DESCRIPTION

## Key Components

The Beamer.Tech system includes a versatile 5-program controller, designed to manage the Beamer PEMF Mat efficiently. The controller allows users to select from multiple therapy modes, adjusting both treatment duration and intensity levels for a customized experience.

### The system is equipped with:

- A single high-performance Beamer Mat optimized for PEMF therapy
- A dedicated Beamer Mat Controller with an intuitive interface
- Durable connection cables for secure and efficient power transmission
- A universal power adapter with interchangeable plugs for international compatibility
- A convenient storage and travel bag for easy portability

## Advanced PEMF Therapy for Holistic Wellness

Unlike traditional light therapy devices, the Beamer Mat utilizes Pulsed Electromagnetic Field (PEMF) technology to deliver targeted frequencies that help:

- Enhance circulation and oxygenation
- Support muscle recovery and pain relief
- Promote stress reduction and relaxation

Designed for both home and clinical use, the Beamer Mat provides a non-invasive, drug-free approach to recovery and overall well-being.

## Benefits

### • Enhancing Cellular Function

The Beamer-Tech Mat is designed to support the body's natural healing processes by delivering targeted Pulsed Electromagnetic Field (PEMF) frequencies. These optimized energy waves help enhance cellular regeneration, promote tissue repair, and improve overall well-being.

### • Advanced PEMF Technology for Recovery

By emitting precise PEMF frequencies, the Beamer-Tech Mat helps stimulate key biological processes, including:

- Reducing inflammation
- Alleviating pain
- Enhancing the body's immune response
- Accelerating recovery times

### • A Non-Invasive Solution for Wellness

Beamer-Tech's PEMF technology is a drug-free, non-invasive approach to supporting health and recovery. It helps optimize cellular function, promote relaxation, and aid in pain management, making it an ideal tool for daily wellness and therapeutic use.

### • Full-Body PEMF Therapy

The Beamer-Tech Mat is equipped with a versatile, full-body therapy design, allowing users to experience deep-penetrating PEMF treatment for improved circulation, muscle relaxation, and overall health benefits.

# How PEMF Works?

Built on safe, trusted, and certified PEMF technology, Beamer PEMF delivers proven benefits for relaxation, sleep improvement, meditation support, pain management, and overall well-being. It is the only wearable device that works directly at the cellular level, offering unique advantages unmatched by any other product

## PEMF: A Natural and Innovative Approach to Wellness

Discover how PEMF technology can enhance your well-being by promoting healing, relaxation, and vitality.

- **Stimulates the Vagus Nerve:** Helps regulate the nervous system by activating the vagus nerve, fostering relaxation, reducing stress, and improving overall well-being.
- **Promotes Relaxation:** Functions like a cellular spa, easing stress and tension while enhancing sleep, focus, and meditation.
- **Relieves Pain:** Effectively alleviates pain by reducing inflammation and improving circulation in targeted areas.
- **Balances Cellular Functions:** Optimizes cellular activity to support overall health and equilibrium.
- **Enhances Blood Flow:** Boosts circulation, ensuring efficient delivery of essential nutrients and oxygen.

## Different Frequencies, Different Benefits

Frequencies play a vital role in influencing our well-being and optimizing body functions in a way that feels natural and familiar.

- **Delta Waves:** Promote deep relaxation and better sleep quality.
- **Theta Waves:** Encourage relaxation, reduce stress, and enhance memory.
- **Alpha Waves:** Create calmness and support learning and mental clarity.
- **Beta Waves:** Boost focus, quick thinking, and mental performance.
- **Gamma Waves:** Enhance clarity and complex information processing.

## WHAT WILL YOU FEEL AFTER USING OUR PEMF DEVICE?



### Deep Relaxation

Experience a sense of calm and stress relief.



### Enhanced Mood

Feel more balanced, positive, and refreshed.



### Increased Energy

Enjoy a natural boost in vitality and focus.

- **Natural & Non-Intrusive:** PEMF enhances well-being naturally, without invasive methods.
- **Time-Efficient:** Quick sessions deliver benefits fast—ideal for a clear, energized start to your day.
- **Full-Body Biohacking:** PEMF promotes both mental and physical wellness for total-body health.
- **Science-Backed:** Research confirms PEMF's role in relaxation, recovery, and overall wellness.
- **Easy to Use:** NeoRhythm is simple, accessible, and user-friendly for everyone.

# Technical Description

## The Beamer-Tech PEMF Therapy Device

### Frequency range

**Tested** 0:01 Hz - 50 MHz

Transmits optimized PEMF frequencies to support circulation, pain relief, and recovery

### PEMF Frequency Wave and Efficiency

Delivers all PEMF therapy frequencies through targeted electromagnetic pulses for maximum therapeutic benefit

Generates optimized waveforms to ensure high efficiency and effectiveness in full-body treatment

### Intensity Settings for Beamer-Tech PEMF Mat

To begin treatment, set the Beamer-Tech Controller to the LOW intensity level. Gradually increase to HIGH for enhanced therapy sessions. Users can customize their treatment duration with selectable 30, 60, 90, or 120-minute programs.

- Ensure the targeted area is fully in contact with the PEMF mat for optimal therapy.
- If higher intensity is not required, revert to the LOW setting for comfort.
- For sensitive users, a thin cotton sheet can be placed over the mat for added comfort.

### Clinic Use and Cleaning

- Designed for both home and clinical use, making it an ideal non-invasive wellness solution.
- For multi-user environments, wipe the surface with a damp cloth after each use.
- Avoid using alcohol-based cleaners or submerging the mat in water, as this will damage the electronics.

# Beamer PEMF

Pulsed Electro Frequency Magnetic Field Mat


[www.beamer.tech](http://www.beamer.tech)

# How To Connect

## Connecting Power to the Beamer Mat Controller

The Power adapter supplied with the Beamer Mat Controller can be used on a variation of world voltages (110v – 250v). The power adapter is supplied with multiple country plugs for international use.

Connect the power adapter to the mains, and connect the power plug to the Mat Controller power port. Once you have connected the power plug to the Mat Controller, the Red power indicator lights up

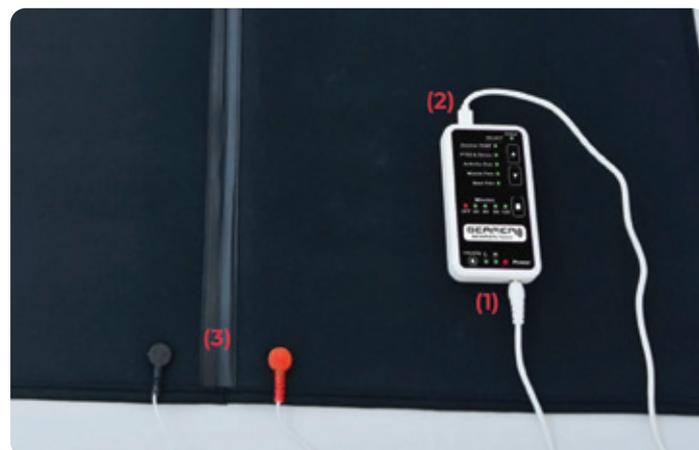

## Connecting the Beamer Mat to the Mat Controller

- First, Connect the Red and Black Buttons into the Beamer PEMF Mat.
- Located on the top of the Beamer Mat are a red and black spheres
- This will indicate which color plug should be connected to the mat. The corresponding plug connection points are located on the bottom of the mat. Plug each connector in firmly.




## How to Use Beamer Mat

- Lay the mat on a flat surface – It can be placed on a bed, couch, or treatment table.
- Connect the Beamer Mat Controller – Plug in the red and black cables to the corresponding ports on the mat.
- Select your therapy mode – Choose from five targeted programs:
  - Beamer PEMF – General wellness
  - PTSD & Stress – Relaxation and mental clarity
  - Arthritic Pain – Joint pain relief
  - Muscle Pain – Muscle recovery
  - Back Pain – Spinal and lower back support
- Set the treatment duration – Options include 30, 60, 90, or 120 minutes.
- Adjust intensity – Choose between Low (L) and High (H) settings.
- Relax and enjoy the session – Users can lay directly on the mat or use a thin sheet for added comfort.

# Controller Functions

## The Beamer Controller

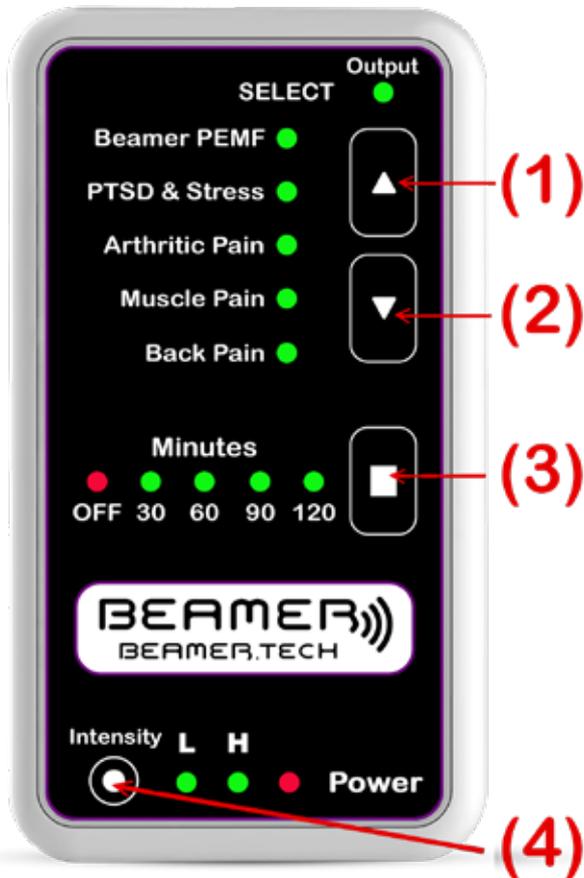

The Beamer Controller

### (1) BUTTON 1: UP

Use the UP arrow to navigate through Beamer PEMF Categories 1-5.

### (2) BUTTON 2: DOWN

Use the DOWN arrow to navigate through Beamer PEMF Categories 1-5.

### (3) BUTTON 3: TIME ADJUSTMENT

Use to switch treatment time from OFF = No Output to Device: 30 Min, 60 Min, 90 Min, 120 Min

### (4) BUTTON 4: POWER INTENSITY

L = Low Power PEMF output

H = High Power PEMF output



**Output LED ON:** The Controller is actively transmitting PEMF frequency signals to the Beamer PEMF Mat, delivering targeted electromagnetic pulses for optimized recovery and wellness.

**Power LED ON:** The Controller is connected to Power Source (Power Adapter)

### \* CONTROLLER PORT SIDE VIEW

(1) POWER PORT (Power Input)

(2) USB-C to

CS0126 Socket Cable (3)



### Device Compatibility

The Beamer-Tech Mat is fully compatible with:

- RDPV3
- RDPV4
- Beamer-Tech Controller

With high-quality connectors and a universal design, the Beamer-Tech Mat ensures reliable operation across multiple platforms, making it a versatile and effective PEMF therapy solution.

## Beamer-Tech Mat Connectivity & Compatibility

The Beamer-Tech Mat is designed for seamless connectivity with its controller and other compatible devices.

### Connection Cables Specifications

1. RDPV4 & Beamer Controller Interface Cable
  - Type: USB-C to Dual CS0126 Socket (Black/Red)
  - Function: Provides standardized connectivity between the RDPV4 unit & Beamer PEMF Mat.
  - Function: Provides standardized connectivity between the Beamer Controller & Beamer PEMF Mat.
2. RDPV3 Legacy Interface Cable
  - Type: RCA to Dual CS0126 Socket (Black/Red)
  - Function: Enables backward compatibility with RDPV3 systems.



USB-C



CS0126 Socket Cable



Black+Red CS0126 Connector Cable

## Timing Selection For Treatments

### Available Treatment Time Options:

- 30 minutes – Minimum treatment time, ideal for quick relaxation and circulation support.
- 60 minutes – Standard session for muscle recovery and pain relief.
- 90 minutes – Extended therapy for deep relaxation and targeted healing.
- 120 minutes – Maximum session time for comprehensive full-body PEMF benefits.

## Suggested Intensity Setting For Treatments

### Recommended Settings:

- Start with Low Intensity (L) – Ideal for first-time users or those with sensitive areas.
- Increase to High Intensity (H) – Recommended for deeper muscle relaxation, enhanced circulation, and recovery support.
- Adjust as Needed – If any discomfort occurs, revert to Low Intensity for a more comfortable experience.

# Advanced machine connections

Connecting the Beamer Mat to the Rife Digital Professional V4



4 USB-C Outputs  
for 4 People or 4 Mats



*Rife Digital Professional V4 (RDPV4)*



If the practitioner or personal user desires to select specific PEMF frequency groups; then we suggest to use the Rife Digital Professional V4 Rife Machine in conjunction with the Beamer Mat, and select a preset program group that is the closest match to your requirement.

Connect the White cable with the Red and Black plugs to the Beamer PEMF Mat.

The Beamer Mat has a Red and Black indicator on the mat to show you where you which color plug to connect to mat

Connect the Beamer Mat directly to the Rife Machine, and select "Pulse Sweep" on the Rife Digital Professional V4 Menu. Pulse Sweep is the PEMF "SWEEP" command. Pulse Sweep is the recommended output command to be used with Beamer PEMF mat, and Beamer PEMF accessories.




When selecting an individual Group number (1-30) on Rife Digital Professional V4, then select "PULSE SWEEP" on the menu as the recommended output command


## BEAMER PEMF (PROGRAM #1)

### Beamer PEMF Mat – Program 1 (Pain Relief)

The Beamer PEMF (Pulsed Electromagnetic Field) device is a full-body therapeutic system consisting of a laydown mat and a five-program controller. Program 1 is designed specifically to relieve general body pain. It features ten targeted frequencies, sourced from ETDFL.com, ranging from 0.16 kHz to 750 kHz. Each frequency runs for 10 minutes, working in synergy to promote healing, reduce pain, and support overall well-being.

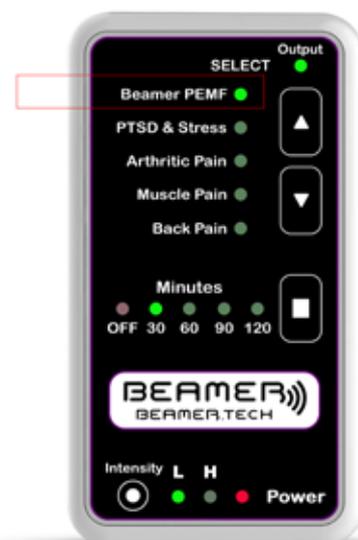

PEMF therapy emits pulsed electromagnetic waves that stimulate the body's natural electrical and cellular functions. In Program 1, these low-frequency pulses penetrate the body, helping to enhance blood flow, improve cellular metabolism, reduce inflammation, and repair tissues—key elements in pain reduction.

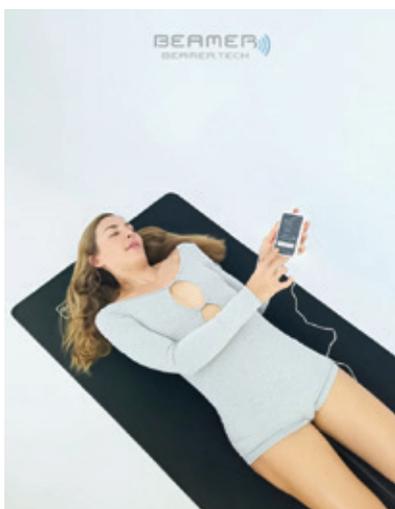

Lower frequencies such as 0.16 kHz and 0.35 kHz target the nervous system and penetrate deeply into tissues to reduce the brain's perception of pain. Mid-range frequencies like 5.26 kHz and 27.5 kHz support muscle relaxation and promote the repair of strained or overworked muscles. Higher frequencies such as 225.47 kHz and 750 kHz improve oxygenation, detoxification, and microcirculation—further enhancing pain relief and tissue recovery.

Users can select session durations of 30, 60, 90, or 120 minutes. For chronic pain conditions such as arthritis, fibromyalgia, or muscle stiffness, 90 to 120 minutes is recommended for best results. Extended sessions stimulate ATP (cellular energy) production, improve nutrient delivery, and reduce stiffness and inflammation.

The full-body Beamer Mat offers flexibility—place it on a bed, sofa, or recliner for comfortable use. Its non-invasive design makes it a safe, natural alternative or complement to medication and physical therapy in pain management.

## PTSD & Stress (PROGRAM #2)

### Beamer PEMF Mat – Program 2 (Stress & PTSD Relief)

The Beamer.Tech PEMF Mat features Program 2, designed to alleviate stress and PTSD symptoms using pulsed electromagnetic field (PEMF) therapy. It delivers ten specific frequencies ranging from 0.14 kHz to 912.33 kHz, each running for 10 minutes. Sessions can be set to 30, 60, 90, or 120 minutes, with 90–120 minutes recommended for deep nervous system engagement and long-lasting benefits.



PEMF works by emitting low-frequency pulses that interact with the body's bioelectrical system, enhancing circulation, cellular repair, and neural balance. Program 2 targets the limbic and autonomic nervous systems—key areas linked to emotional regulation and the fight-or-flight response. Lower frequencies like 0.68 kHz and 2.5 kHz activate the parasympathetic system to promote calmness, while higher frequencies like 655.2 kHz and 912.33 kHz support neural synchronization and reduce inflammation.



This program also promotes brainwave entrainment. Frequencies in the 0.14–60 kHz range stimulate delta and theta waves for deep relaxation, while mid-range frequencies (122.53–300 kHz) support alpha waves for calm alertness. These patterns can reduce cortisol levels, improve sleep, and ease physical symptoms like muscle tension and fatigue.

Notably, 496.01 kHz and 750 kHz may stimulate neuroplasticity—helping the brain form healthier responses to traumatic triggers. The mat's full-body design allows users to recline comfortably, making it ideal for those with mobility issues or high anxiety. Wearing light clothing during sessions ensures optimal field absorption.

Clinical feedback indicates that regular use of Program 2 may reduce anxiety, improve mood, and enhance resilience. As a drug-free, non-invasive option, it offers a powerful complement to traditional therapies for stress and trauma recovery.

## Arthritic Pain (PROGRAM #3)

### Beamer PEMF Mat – Program 3 (Arthritis & Joint Relief)

The Beamer.Tech PEMF Mat offers Program 3, a targeted sequence for arthritis, osteoarthritis, and rheumatism. It delivers ten carefully selected frequencies (from 0.05 kHz to 642.91 kHz), each pulsed for 10 minutes over sessions lasting 30 to 120 minutes. A 90–120 minute duration is recommended for deeper penetration into joints and connective tissues, enhancing relief and healing.



PEMF therapy works by stimulating cellular repair and reducing inflammation. Program 3's frequencies, sourced from ETDFL.com, are chosen to address cartilage regeneration, synovial fluid production, and pain modulation. Lower frequencies (0.05–0.9 kHz) target nerve pain and joint stiffness. Mid-range frequencies (9–55.33 kHz) enhance ATP production in affected areas. Higher frequencies (398.4–642.91 kHz) may suppress pro-inflammatory cytokines and support tissue regeneration.

Notably, frequencies like 11.09 kHz and 225.47 kHz may promote collagen production and slow cartilage breakdown. Frequencies such as 425.71 kHz and 522.53 kHz are believed to support synovial fluid balance, improving joint lubrication. This frequency blend supports both wear-and-tear and autoimmune joint conditions.

Pain relief is achieved through neuromodulation, reducing nerve signal transmission and stimulating natural endorphin release. Frequencies like 642.91 kHz are particularly useful in decreasing pain sensitivity and promoting relaxation.



The mat's full-body design allows users to recline comfortably while targeting specific joints. Placing painful areas directly on the mat maximizes local benefits, while systemic coverage addresses widespread inflammation.

Program 3 is non-invasive and drug-free, with research supporting its use in slowing arthritis progression. Regular 3–5 weekly sessions may reduce reliance on NSAIDs and improve joint mobility when paired with lifestyle care such as movement and anti-inflammatory diets.

## Muscle Pain (PROGRAM #4)

### Beamer PEMF Mat – Program 4 (Muscle & Nerve Pain Relief)


Program 4 of the Beamer.Tech PEMF Mat is designed to relieve muscle tension and nerve-related pain using ten targeted frequencies (0.03 kHz to 568.43 kHz), each pulsed for 10 minutes. Sessions range from 30 to 120 minutes, with 90–120 minutes recommended for optimal results. The full-body mat allows users to lie or sit comfortably, offering both localized and systemic relief.



PEMF therapy addresses pain through cellular repair, reduced inflammation, and neural modulation. Lower frequencies (0.03–0.9 kHz) relax deep muscle fibers and calm peripheral nerves. Mid-range frequencies (7.5–20 kHz) boost ATP production in stressed tissues, aiding recovery. Higher frequencies (322.06–568.43 kHz) are effective in disrupting pain signals in acute or chronic nerve pain.

Program 4 also targets nerve function by interacting with voltage-gated ion channels. Frequencies like 0.78 kHz and 377.3 kHz help stabilize overactive nerves, common in sciatica and neuropathy. Meanwhile, 425.71 kHz may reduce substance P, a chemical linked to heightened pain perception.

Muscle relaxation is enhanced through calcium ion modulation. Frequencies such as 0.37 kHz and 7.5 kHz assist in releasing contracted myofibrils, similar to the effects of massage therapy. Users often describe a melting sensation in tight areas and improved flexibility post-session.

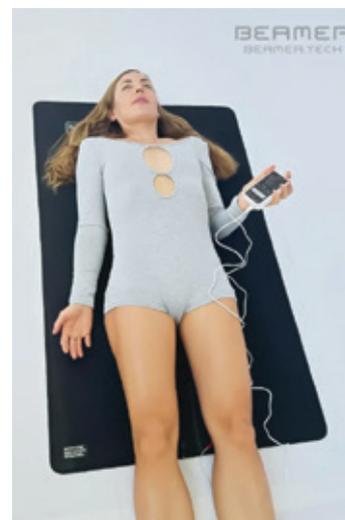

The mat's broad coverage addresses both direct pain points (e.g., stiff shoulders) and secondary tension zones (e.g., back pain from muscle guarding). Light clothing is advised to ensure effective field transmission during longer sessions.

Clinical data supports Program 4's effectiveness, showing up to 50% pain reduction after repeated use. For lasting results, 5+ sessions per week are recommended, paired with hydration and gentle movement to help detoxify and repair treated tissues.

## Back Pain (PROGRAM #5)

### Beamer PEMF Mat – Program 5 (Back Pain Relief)

Program 5 is designed to target multiple causes of back pain using ten specialized frequencies (0.14–932 kHz) delivered in 10-minute pulses. Recommended sessions range from 90 to 120 minutes for optimal results. The program addresses inflammation, muscle tightness, nerve irritation, and degenerative changes—common sources of both acute and chronic back pain.



Lower frequencies (0.14–7.5 kHz) release tension in deep muscles like the erector spinae and QL, which often contribute to postural and mechanical pain. Mid-range frequencies (55–96.6 kHz) enhance tissue repair in spinal discs, ligaments, and surrounding connective structures. Higher frequencies (376.29–932 kHz) target spinal nerves to reduce pain signal transmission and promote neural balance.

For nerve-related back pain, frequencies such as 425.09 kHz and 571 kHz help calm hyperactive nerve membranes and lower levels of substance P, a pain-related neurotransmitter. The 833 kHz frequency has shown promise in reducing radiating pain from lumbar compression, while 932 kHz may help break muscle spasm cycles in the paravertebral region.



Structural pain, including osteoarthritis in facet joints, also benefits from frequencies like 96.6 kHz, which may slow cartilage breakdown and improve joint fluid production. Fibroblast activity is supported by 55 kHz and 376.29 kHz, promoting soft tissue repair.

The Beamer mat's full-body layout allows users to treat specific areas like lumbar or thoracic regions while also addressing compensatory tension and systemic inflammation. Light clothing is advised during sessions.

Clinical studies suggest consistent use may reduce pain scores by 40–60%, improve mobility, and reduce dependence on medication. For best results, 4–5 weekly sessions paired with posture care and light movement are recommended.



# Advanced Treatments

## Beamer-Tech PEMF System and RDPV4 Compatibility

The Beamer-Tech PEMF System offers a comprehensive solution for recovery, circulation, and pain relief using Pulsed Electromagnetic Field (PEMF) therapy. Its mat can be used independently or paired with the Rife Digital Professional V4 (RDPV4) to boost frequency-based benefits.

## Expanding Frequency Applications with RDPV4

A key feature of RDPV4 is its precise frequency transmission for targeted wellness. When integrated with the Beamer-Tech Mat, users can access a wide range of preset programs to optimize their therapy sessions.

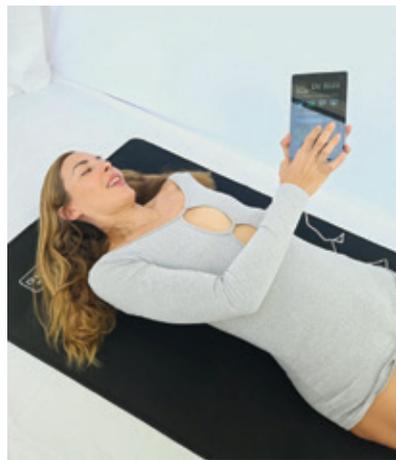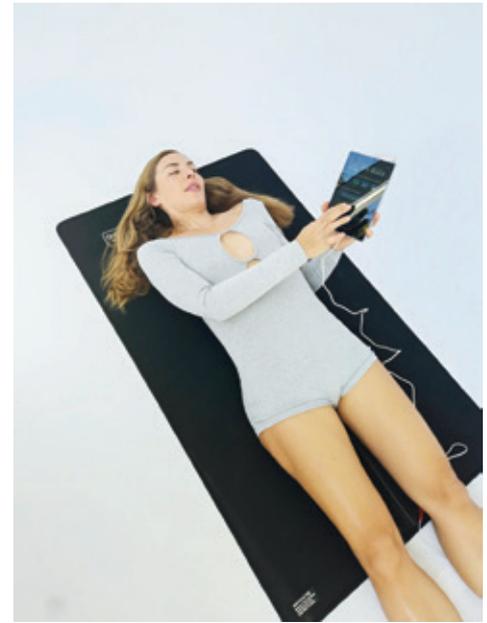

Together, the mat and RDPV4 deliver targeted PEMF pulses to support:

- Enhanced circulation and tissue oxygenation.
- Reduced inflammation and muscle tension.
- Deep relaxation and stress relief.
- Recovery from muscle and joint strain.

With RDPV4's dual output, users can target specific areas while enjoying full-body PEMF therapy—offering more flexibility in treatment.

## Advanced Beamer-Tech Controller

The Beamer-Tech controller includes:

- Five preset PEMF programs for relaxation, recovery, and circulation.
- Low/High intensity options for individual comfort.
- Timer settings for 30, 60, 90, or 120 minutes.
- Secure ports for stable PEMF transmission.



## A Non-Invasive, Drug-Free Wellness Solution

By combining PEMF with RDPV4 frequencies, the Beamer-Tech system delivers a powerful, customizable, and non-invasive therapy—ideal for both home and professional wellness use.



Our bodies are made up of a total of eleven organ systems that manage all the essential body functions. While interdependent, each system can become out of tune. When one is not functioning properly, the others will attempt to correct the problem. All body systems will work together to try to create a state of balance within the body.

The 11 organ systems include the integumentary system, skeletal system, muscular system, lymphatic system, respiratory system, digestive system, nervous system, endocrine system, cardiovascular system, urinary system, and reproductive systems.

The Beamer PEMF technology uses Pulsed frequencies, transmitted through a large mat to balance the human body organ system as a whole body healing modality. The PEMF method of treatment has been in use successfully for over 60 years, and continues to be a primary method of achieving homeostasis and a state of balance within the body.

The following chapters will explain how the Beamer Pulsed frequency technology will assist with the Human body and organ systems



1. Integumentary System  
(Skin, Hair, Nails)



2. Skeletal System  
(Bones, Joints)



3. Muscular System (Cardiac,  
Smooth, and Skeletal Muscles)



4. Lymphatic &  
Immune System



5. Respiratory System  
(Lungs, Airways)



6. Digestive System  
(Stomach, Liver, Intestines)



7. Nervous System  
(Brain, Spinal Cord, Nerves)



8. Endocrine System  
(Hormonal Balance)



9. Cardiovascular (Circulatory)  
System (Heart, Blood Vessels)



10. Urinary System  
(Kidney, Bladder, Urethra)




11. Reproductive System  
(Male & Female Health)



12. Other Categories

# INTEGUMENTARY SYSTEM

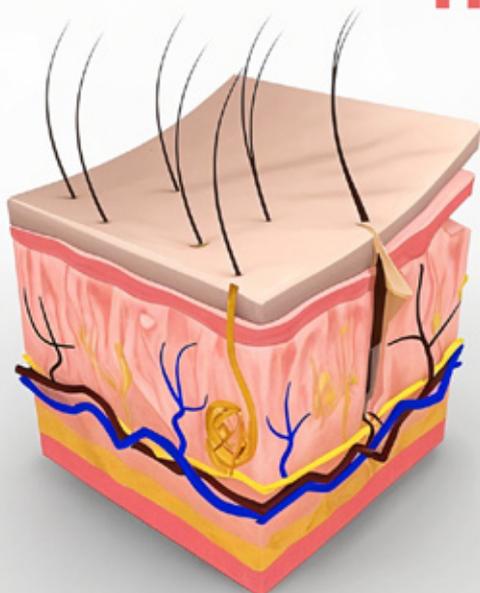

Equaling 15-20% of our total body mass, this system includes the **skin, hair, and nails.**

## Integumentary System (Skin, Hair, Nails)

The integumentary system is the body's largest organ, making up 15–20% of body mass. It serves as a protective barrier against physical, chemical, and biological threats. The skin regulates perspiration and body temperature, while also transmitting sensations like touch, pain, and pleasure. Hair helps regulate scalp temperature and secretes pheromones. Nails protect the fingertips, essential for tasks like self-defense, food gathering, and sensing. When the skin is wounded, it becomes vulnerable to infection.

## Beamer PEMF Targets

Beamer PEMF targets various skin-related conditions including acne, chloracne, dermatitis, eczema, hives, lupus, psoriasis, autoimmune disorders, keratin and pore issues, and skin cancers like melanoma and basal cell carcinoma. Using pulsed electromagnetic frequencies, the Beamer mat promotes whole-body healing.

## Suggested Treatment Method

Overnight treatments are most effective for skin-related conditions. Daytime use is also beneficial—simply lay the mat on a recliner, sofa, or bed. You can lie directly on the mat or place a sheet over it. Connect the red and black cables to the mat's base near the feet, then follow the setup guide to link with the controller or Rife device.

# SKELETAL SYSTEM

---

The skeletal system, containing bones and joints, supports and protects the body's internal organs.



## Skeletal System (Bones, Joints)


The skeletal system protects vital organs and supports body structure. Ribs shield abdominal organs, while the skull protects the brain, which controls all bodily functions. Bones provide shape, support movement through muscle attachment, store minerals like calcium, and produce blood cells in the marrow.

## Beamer PEMF Targets

Beamer PEMF addresses skeletal conditions such as arthritis, cancer, scoliosis, osteoporosis, gout, bursitis, fractures, and amputations. Through pulsed electromagnetic frequencies, it supports whole-body healing and bone health.

## Suggested Treatment Method

Overnight use is recommended for skeletal issues. For daytime sessions, place the mat on a recliner, sofa, or bed. You may lie directly on the mat or use a bedsheet. Connect the red and black cables at the base (near the feet) and follow the setup guide to link with the controller or Rife Machine. This setup promotes consistent therapy for bone and joint recovery.

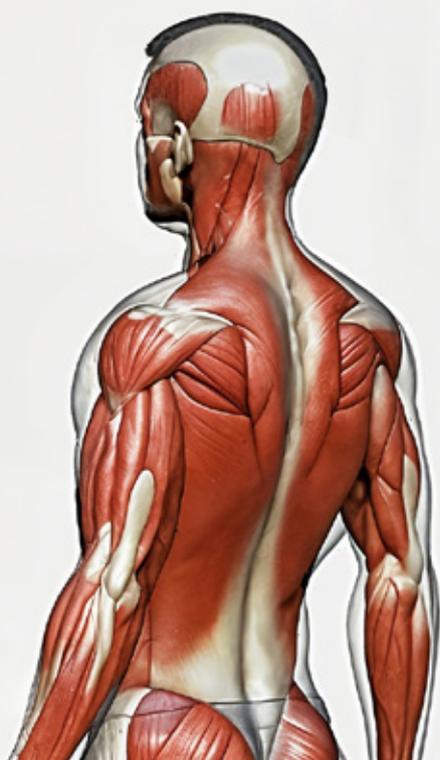

# MUSCULAR SYSTEM

---

Voluntary, involuntary and cardiac muscles compose the muscular system

## **Muscular System (Cardiac, Smooth, and Skeletal Muscles)**

Cardiac muscles, found in the heart, power blood flow throughout the body. Smooth (involuntary) muscles line organs and control internal movements like digestion.

Skeletal (voluntary) muscles respond to brain signals via the nervous system, enabling movement, maintaining posture, and producing body heat. A lack of oxygen can cause muscle cramps or tears, and disuse can lead to atrophy.

## **Beamer PEMF Targets**

Beamer PEMF supports muscular health by targeting conditions such as muscular dystrophy, fibromyalgia, tendinitis, multiple sclerosis, muscle strains, sprains, and hernias. Its pulsed electromagnetic frequencies promote whole-body healing by improving muscle function and recovery.

## **Suggested Treatment Method**

Overnight sessions are ideal for muscular disorders. For daytime use, lay the Beamer mat on a bed, recliner, sofa, or lounge. A sheet may be placed over the mat, or users can lie directly on it. Connect the red and black cables near the mat's foot end, then follow the setup guide to link with the controller or Rife Machine. This setup allows effective, targeted support for cardiac, smooth, and skeletal muscles.

# LYMPHATIC SYSTEM

---

Red bone marrow, thymus, lymphatic vessels, thoracic duct, spleen, and lymph nodes make up this system responsible for fluid transport.



## **Lymphatic System (Bone Marrow, Thymus, Lymphatic Vessels, Thoracic Duct, Spleen, Lymph Nodes)**


The lymphatic system helps maintain fluid balance by returning clean fluids to the bloodstream and draining excess fluids and waste from tissues. It also plays a key role in immunity by housing lymphocytes (white blood cells) that defend the body against infections.

### **Beamer PEMF Targets**

Beamer PEMF supports lymphatic health by targeting conditions like anemia, leukemia, tuberculosis of the lymph nodes, Hodgkin's disease, and other blood-related disorders. Using pulsed electromagnetic frequencies through a full-body mat, the technology promotes healing and balance across the body's systems.

### **Suggested Treatment Method**

Overnight sessions are ideal for lymphatic-related conditions. For daytime treatments, place the mat on a recliner, bed, or sofa. You may lie directly on it or use a bedsheet. Attach the red and black connectors at the base of the mat (near the feet), and follow the setup instructions to connect with the controller or Rife Machine. This method provides consistent and non-invasive support for the lymphatic system's cleansing and immune functions.

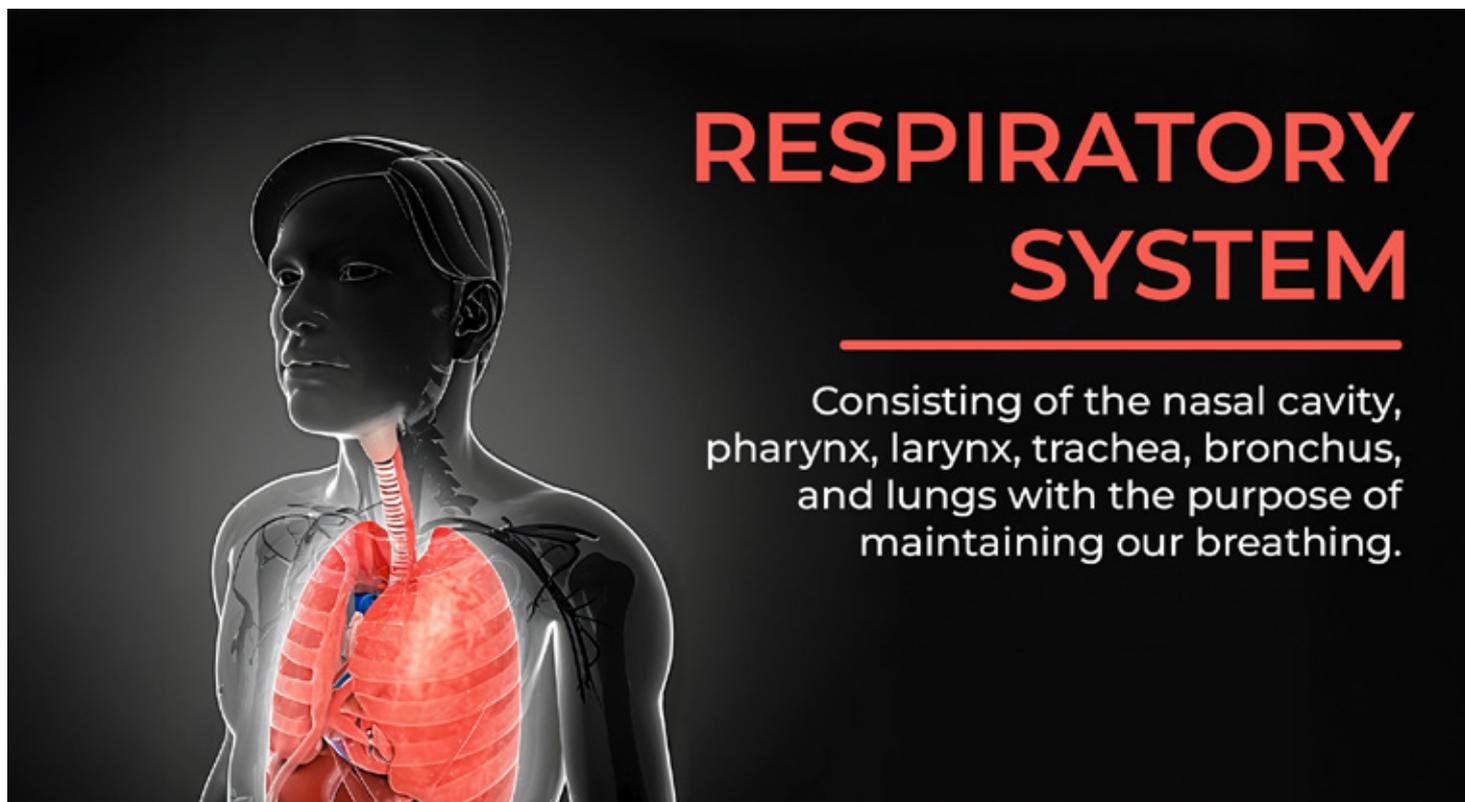

# RESPIRATORY SYSTEM

Consisting of the nasal cavity, pharynx, larynx, trachea, bronchus, and lungs with the purpose of maintaining our breathing.

## Respiratory System (Nasal Cavity, Pharynx, Larynx, Trachea, Bronchus, Lungs)

The respiratory system enables breathing by delivering oxygen to the lungs for cellular energy and removing carbon dioxide as a waste product. It also supports speech and the sense of smell.

### Beamer PEMF Targets

Beamer PEMF addresses various respiratory conditions including allergies, rhinitis, sinusitis, laryngitis, COPD, pleurisy, bronchitis, emphysema, asthma, sarcoidosis, fibrosis, asbestosis, pulmonary vascular disease, lung infections, sleep apnea, tuberculosis, and cancers of the lungs or throat. Using pulsed electromagnetic fields, the Beamer mat helps balance and heal the body through full-body frequency therapy.


### Suggested Treatment Method

Overnight sessions are recommended for respiratory issues. For daytime use, place the Beamer mat on a bed, recliner, or sofa. You can lie directly on it or with a bedsheet. Connect the red and black cables to the mat's base (near the feet), then follow the setup guide to connect with the controller or Rife Machine. This setup offers a gentle, non-invasive approach to respiratory health support.

# DIGESTIVE SYSTEM

---

For absorption of nutrients, the digestive system includes the oral cavity, esophagus, liver, stomach, small intestine, large intestine, rectum, and anus.

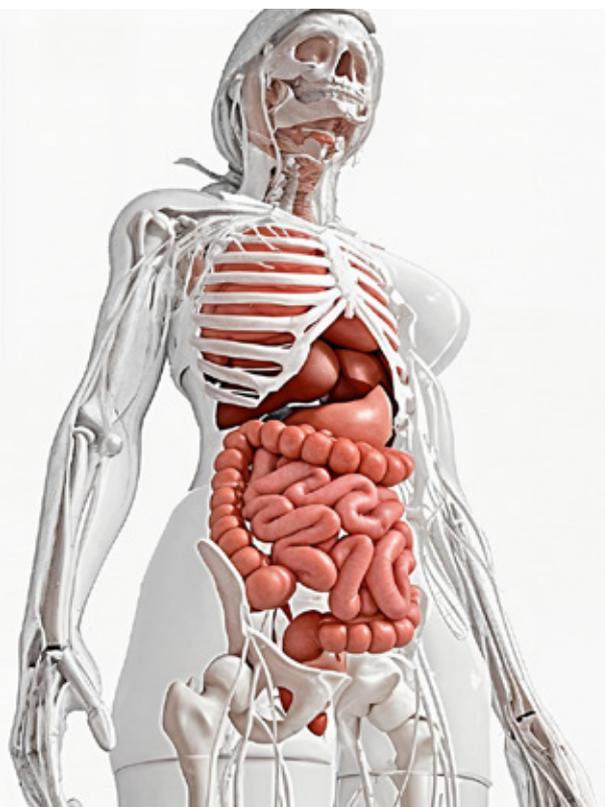

## **Digestive System (Oral Cavity, Esophagus, Liver, Stomach, Small & Large Intestines, Rectum, Anus)**


The digestive system begins in the mouth and is responsible for breaking down food, absorbing nutrients, and eliminating waste. It determines which vitamins, minerals, and nutrients can be used or stored by the body and removes the rest.

### **Beamer PEMF Targets**

Beamer PEMF supports digestive health by targeting conditions such as diverticulitis, gastritis, pancreatitis, cholecystitis, cirrhosis, hepatitis, liver cancer, IBS, Crohn's disease, and hemorrhoids. Pulsed electromagnetic frequencies help restore balance and support healing across the digestive tract.

### **Suggested Treatment Method**

Overnight therapy is recommended for digestive issues. For daytime sessions, place the mat on a bed, recliner, or sofa. You may lie directly on it or with a bedsheet. Attach the red and black connectors at the base of the mat (near the feet) and follow the setup guide to link it to the controller or Rife Machine. This setup delivers gentle, whole-body PEMF support for digestive system function.

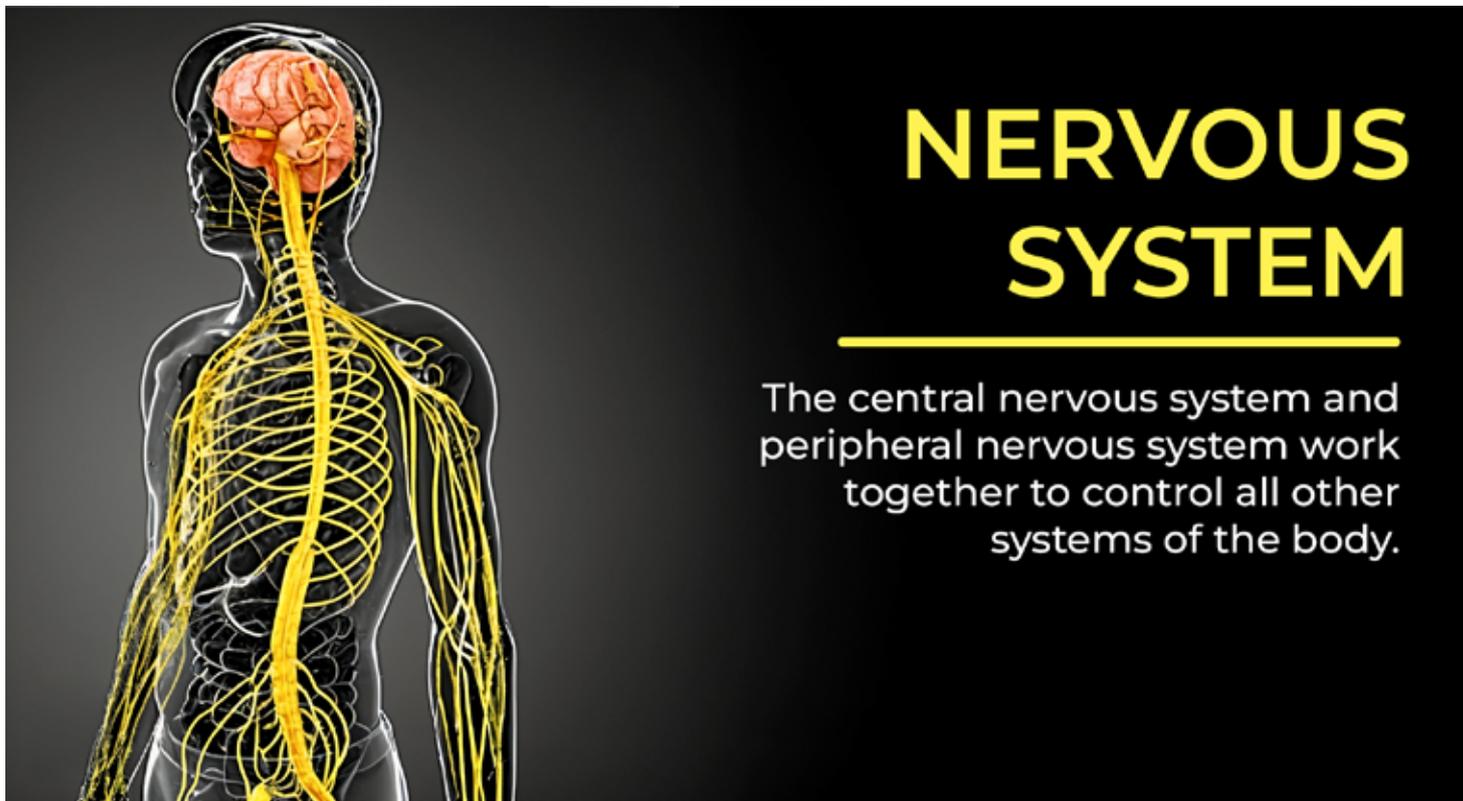

# NERVOUS SYSTEM

---

The central nervous system and peripheral nervous system work together to control all other systems of the body.

## **Nervous System (Brain, Spinal Cord, Nerves)**

The nervous system consists of the central nervous system (CNS)—brain and spinal cord—and the peripheral nervous system, which includes nerves connecting the CNS to the rest of the body. It controls all other body systems, including digestion and heartbeat, and responds to both internal and external stimuli, transmitting sensations such as pain and touch.

## **Beamer PEMF Targets**

Beamer PEMF targets a wide range of nervous system conditions, including paralysis, Parkinson's disease, palsy, embolisms, thrombosis, arteriosclerosis, polio, ALS, meningitis, multiple sclerosis, muscular tics, Huntington's disease, brain or spinal cancers, epilepsy, seizures, narcolepsy, migraines, neuropathy, and traumatic brain injury. Using pulsed electromagnetic fields, Beamer PEMF supports neurological balance and healing.

## **Suggested Treatment Method**

Overnight use is recommended for nervous system issues. For daytime sessions, lay the mat on a bed, recliner, or sofa. You may use a sheet or lie directly on the mat. Connect the red and black cables near the feet of the mat, then follow the setup guide to connect with the controller or Rife Machine. This approach provides non-invasive, full-body nervous system support.

# ENDOCRINE SYSTEM

---

Responsible for releasing hormones, this system includes the pineal gland, pituitary gland, thyroid, thymus, adrenal gland, pancreas, ovary, and testis.



## Endocrine System (Pineal, Pituitary, Thyroid, Thymus, Adrenal, Pancreas, Ovary, Testis)


The endocrine system regulates essential body functions by releasing hormones that control growth, reproduction, metabolism, and blood glucose levels.

### Beamer PEMF Targets

Beamer PEMF supports endocrine health by targeting conditions like Type 1 and Type 2 diabetes, hypoglycemia, Addison's disease, Cushing's syndrome, hyper- and hypothyroidism, thyroid cancer, and other cancers of endocrine glands. Pulsed electromagnetic frequencies transmitted through the mat help restore balance across hormone-regulating organs.

### Suggested Treatment Method

Overnight sessions are ideal for endocrine-related conditions. For daytime treatments, place the mat on a bed, recliner, or sofa. You can lie directly on the mat or use a sheet. Connect the red and black connectors near the base (at the foot of the mat), and follow the setup guide to link with the controller or Rife Machine. This method provides a safe, non-invasive way to support endocrine balance and hormonal wellness.

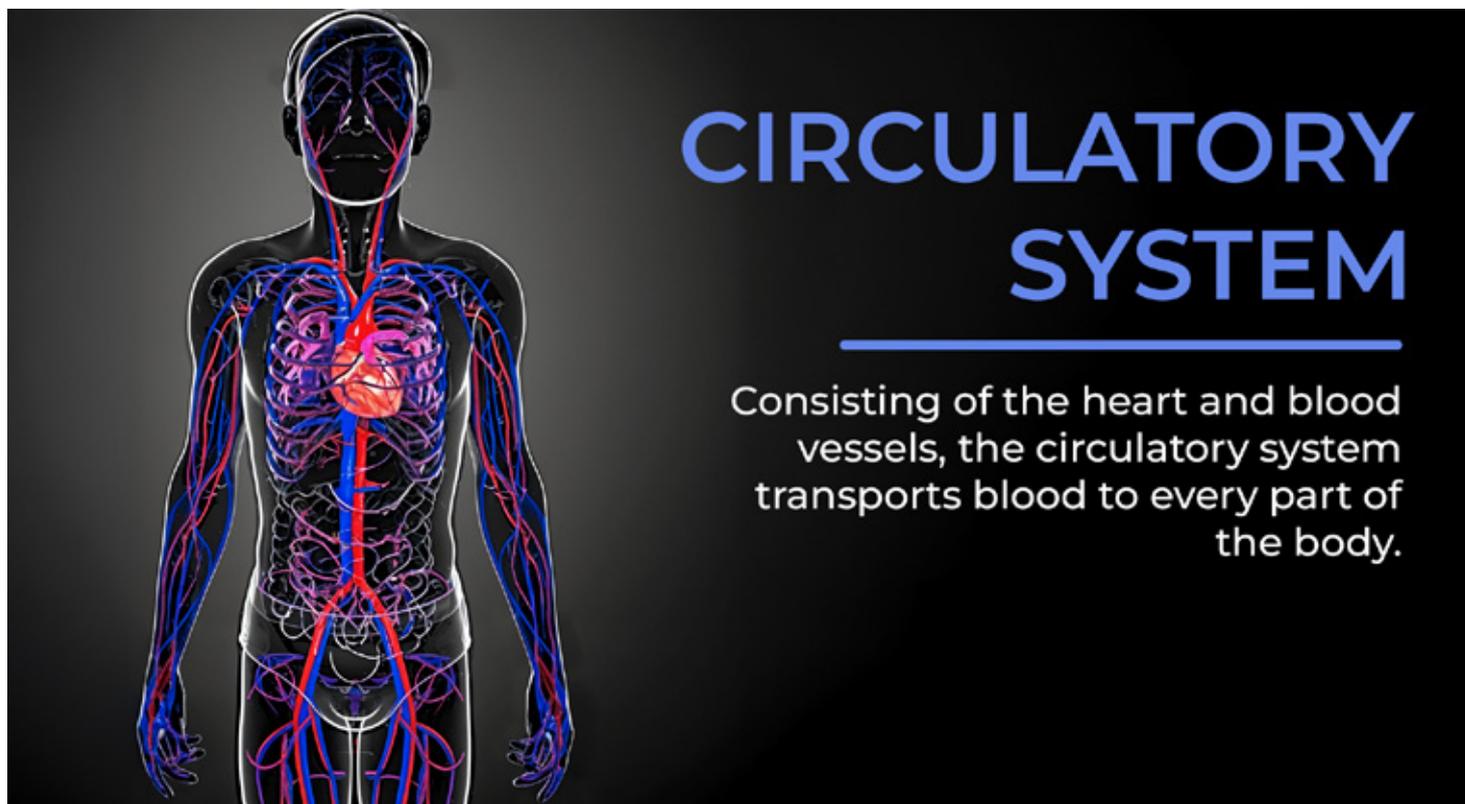

# CIRCULATORY SYSTEM

Consisting of the heart and blood vessels, the circulatory system transports blood to every part of the body.

## Cardiovascular (Circulatory) System (Heart, Blood Vessels)

The cardiovascular system circulates blood throughout the body. The heart pumps blood, while arteries and veins transport oxygen, carbon dioxide, nutrients, and waste, delivering nourishment to organs and muscles.

### Beamer PEMF Targets

Beamer PEMF supports circulatory health by targeting conditions such as myocardial infarction, arrhythmias, valve disorders, pacemaker-related conditions, transplants, ischemic and hypertensive heart disease, arteriosclerosis, aneurysms, fistulas, anaphylactic shock, varicose veins, cold weather injuries, and blood vessel cancers like sarcomas. Pulsed electromagnetic frequencies from the mat help restore vascular balance and function.


### Suggested Treatment Method

Overnight sessions are ideal for cardiovascular issues. For daytime use, lay the mat on a recliner, sofa, or bed. You may lie directly on it or use a sheet. Attach the red and black cables to the connectors at the foot of the mat, then follow the setup instructions to link with the controller or Rife Machine. This setup provides full-body, non-invasive circulatory support.

# URINARY SYSTEM

---

Responsible for eliminating waste, this system includes the kidneys, ureter, urinary bladder, and urethra.

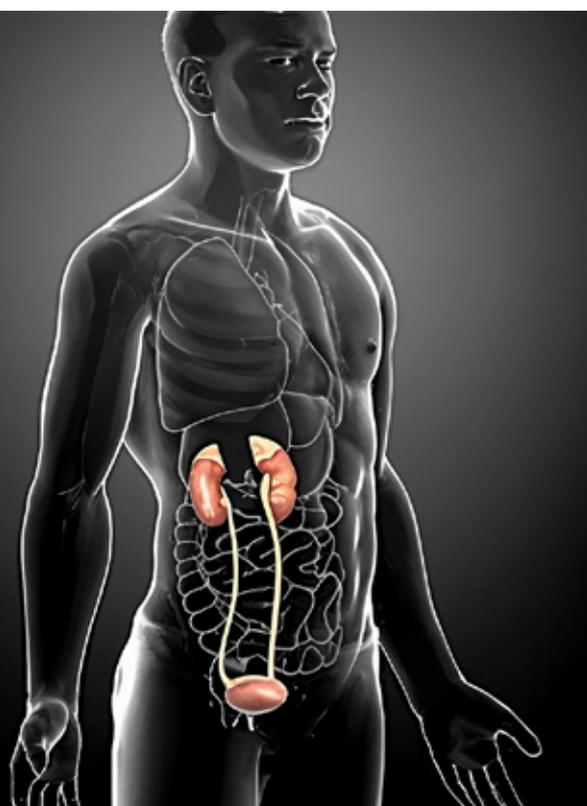

## Urinary System (Kidneys, Ureters, Bladder, Urethra)

The urinary system is responsible for removing metabolic waste and excess substances from the body. It regulates fluid volume, balances electrolytes, and maintains proper blood pH—key functions for overall detoxification and internal balance.

### Beamer PEMF Targets


Beamer PEMF targets a wide range of urinary system conditions, including nephrosis, bladder cancer, urethritis, bedwetting (enuresis), kidney and urinary stones, infections, renal failure, incontinence, blood in the urine, and interstitial cystitis. Using pulsed electromagnetic frequencies, the mat supports the healing and rebalancing of urinary tract functions.

### Suggested Treatment Method

Overnight sessions are highly recommended for treating urinary conditions, offering deep, uninterrupted therapy. Daytime treatments are also effective—place the mat on a recliner, bed, or sofa. You can lie directly on the mat or use a sheet if preferred. Connect the red and black cables near the foot of the mat and follow the setup instructions to link it with the controller or Rife Machine. This full-body PEMF approach offers a non-invasive and convenient way to support urinary health and waste elimination.

# REPRODUCTIVE SYSTEMS

This system is responsible for creating human life, and has varied organs between men and women.

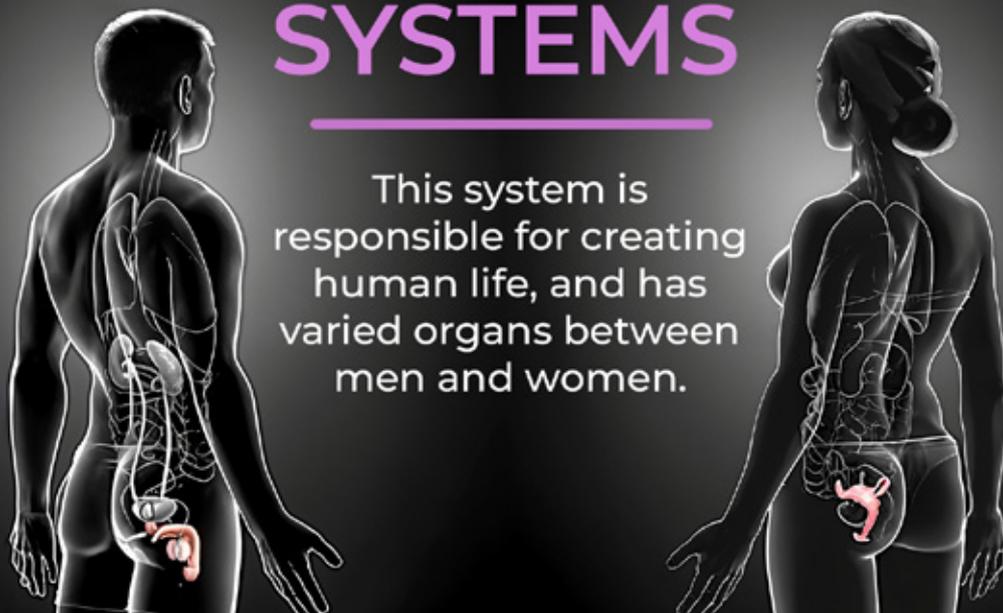

## Reproductive Systems

The reproductive system consists of organs involved in reproduction. In males: testes, penis, seminal vesicles, prostate, and urethra. In females: ovaries, fallopian tubes, uterus, vagina, and vulva. These systems play a key role in fertility, hormone production, and reproductive health.

## Beamer PEMF Targets

Beamer PEMF supports reproductive health by targeting conditions in females such as breast cancer, post-mastectomy or lumpectomy recovery, hysterectomy, endometriosis, pregnancy complications, and cancers or injuries of the ovaries, cervix, uterus, or vulva.

In males, it addresses prostate cancer, testicular cancer, enlarged prostate (BPH), prostatitis, infertility, and testosterone deficiency. The Beamer mat uses pulsed electromagnetic frequencies to restore balance and promote healing in the reproductive system.

## Suggested Treatment Method

Overnight treatments are ideal for reproductive conditions, offering sustained support. For daytime sessions, place the mat on a recliner, bed, or sofa. Lie directly on it or use a bedsheet. Connect the red and black cables near the base of the mat (at the feet) and follow the setup guide to connect with the controller or Rife Machine. This whole-body PEMF therapy provides gentle, non-invasive care for reproductive wellness.

## OTHER CATEGORIES

### **Organs of Special Sense (Eyes)**

The eyes are responsible for vision and are highly sensitive to structural and neurological conditions. Beamer PEMF supports disorders such as blurry vision, visual field loss, scotoma, muscle dysfunction, intraocular hemorrhage, retinal detachment, retinopathy, conjunctivitis, corneal damage and transplants, aphakia, cataracts, glaucoma, optic neuropathy, tear gland issues, and cancers of the eye or related structures.

### **Impairment of Auditory Acuity (Ears)**

The ears regulate hearing and balance. Beamer PEMF may support conditions like hearing loss, tinnitus, chronic ear infections, balance disorders, and ear-related cancers by stimulating tissue repair and neural balance.

### **Mental Disorders (Psychological Issues)**

Beamer PEMF may help regulate brain activity and support recovery in mental health conditions such as anxiety, PTSD, OCD, panic attacks, adjustment disorder, bipolar disorder, depression, schizophrenia, neuro-cognitive disorders, eating disorders, and phobias. If your mental health rating predates 2014, consult a VSO or attorney to review potential updates in classification and benefits.

### **Infectious Diseases, Immune Disorders & Nutritional Deficiencies**

Beamer PEMF targets conditions such as malaria, cholera, trench fever, leprosy, syphilis, HIV, COVID/SARS, lupus, Lyme disease, typhoid, tuberculosis (non-respiratory), and chronic fatigue syndrome. It may help modulate immune response and support cellular healing.

### **Usage with Rife Digital**

Beamer PEMF can be combined with the Rife Digital Professional series for more specific frequency targeting. Users can reference the ETDFL Frequency Manual from ETDFL.com to select the appropriate program, then run it through the Beamer Mat for enhanced results.

## FAQ Frequently Asked Questions

### ■ What Setting on the Rife Professional V4 Should I Use with the Beamer-Tech PEMF Mat?

When using the Beamer-Tech PEMF Mat with the Rife Professional V4, it is recommended to select a preset frequency program that aligns with your wellness goals. The Rife Professional V4 includes multiple frequency sets designed to support pain relief, circulation enhancement, muscle recovery, and relaxation. If not, refer to the ETDFL book (available for free download at [ETDFL.COM](http://ETDFL.COM)) to find the appropriate frequencies for additional wellness-based frequency recommendations

### ■ Where Can I Buy the Latest V4 Rife Machine to Use with My Beamer Accessory?

To access the 3,500 groups and advanced features of the controller, you can purchase the Rife Professional V4 through the Beamer.Tech shopping cart

### I Feel Sensitivity When Using the Beamer-Tech PEMF Mat

If you experience sensitivity while using the Beamer-Tech PEMF Mat, we recommend starting with the "Low" intensity setting on the controller. Ensure that the mat is positioned comfortably on a flat surface and does not apply excessive pressure to the body.

If any discomfort persists, keep the controller at the LOW intensity setting and allow your body to gradually adapt to the therapy. The mat should provide gentle and even coverage without tightness, ensuring a relaxing and effective PEMF session.

### ■ The Beamer-Tech PEMF Mat Doesn't Fit Properly. Help!

The Beamer-Tech PEMF Mat is designed for full-body coverage and should be placed on a flat and comfortable surface for optimal use. Unlike wearable devices, the PEMF mat does not require direct skin contact to be effective.

To ensure the best results:

- Lay the mat flat on a bed, recliner, or treatment table for full-body exposure.
- Position yourself comfortably on the mat, ensuring even distribution of PEMF therapy.
- If needed, use a thin cotton sheet for added comfort without interfering with the PEMF transmission.

\* No straps or adjustments are necessary, as the PEMF waves penetrate the body naturally, providing effective therapy without requiring a tight fit.

### ■ Can I Use the Beamer-Tech PEMF Mat for My Dog, Cat, or Other Animals?

Yes, the Beamer-Tech PEMF Mat can be safely used for dogs, cats, and other animals to support relaxation, circulation, and recovery. Many pet owners use PEMF therapy for animals experiencing joint discomfort, mobility issues, or general wellness support.

However, please be aware that animals may chew or tug on the mat or cables, potentially damaging the equipment. The Beamer-Tech PEMF system operates at a safe low voltage, ensuring no harm to pets in case of accidental contact.

\* Important Note: Any physical damage to the mat or controller caused by pets (such as chewing or tearing) is not covered under warranty. To prevent damage, monitor your pet during use and keep cables out of reach when not in use.

## ■ How Does the Beamer-Tech Controller Actually Work?

*The Beamer-Tech Controller is designed to precisely generate and transmit Pulsed Electromagnetic Field (PEMF) frequencies to the Beamer-Tech PEMF Mat, delivering targeted therapy for circulation, pain relief, and recovery support.*

*Unlike light-based devices, the Beamer-Tech Controller operates by:*

- *Generating scientifically calibrated PEMF signals that interact with the body's natural bioelectrical system.*
- *Delivering pulsed electromagnetic waves through the mat to enhance cellular function and circulation.*
- *Offering adjustable intensity levels and multiple preset frequency programs to cater to individual wellness needs.*

*By optimizing PEMF output, the Beamer-Tech system supports muscle relaxation, inflammation reduction, and overall well-being in a non-invasive and natural way.*

## Warranty & Disclaimer

### Warranty

- **DOA Warranty:** If your Beamer™ arrives with any damage caused by transport/courier, it will be replaced immediately
- **3 Month Replacement Warranty:** If your Beamer™ becomes broken with normal use, we will replace the item within 90 days. Cost of shipping will be charged for the return of a new item to your address.
- Due to sanitary reasons the Beamer™ mat and controller will not be repaired by our service center. We will always issue a new product, under warranty.

### Disclaimer

- Beamer™ holds no guarantee, that the BEAMER™ Mat and controller device will cure any specific diseases. The Beamer™ technologies have been tested for over 3 years in EU by researchers, who have had successful results using the Beamer™ in conjunction with the Rife Digital Series only. The Beamer™ has not been tested with any other brand Rife Machine, so there are no guarantees should you alter the equipment for other uses. The Beamer™ mat itself is a Rubber Conductive Conduit, and in itself does not perform any process without a machine attached.



ETDFL.com



## About E.T.D.F.L and Frequency-Based Therapy with Beamer-Tech

**E.T.D.F.L** (Electro Therapy Device Frequency List) is a global initiative led by 12 clinics around the world using Quantum SCIO Bioresonance technology to collect and refine high-accuracy frequency sets. First launched in 2006 to replace the original CAFL list, ETDFL has since become the standard for over 100,000 Rife Machine users worldwide.

Each frequency set within ETDFL is carefully developed from real clinical data and includes the 10 most effective frequencies for each condition. These are designed for use with devices such as the Rife Digital Professional V4 (RDPV4) and can be transmitted to compatible PEMF systems like the Beamer-Tech PEMF Mat, offering a safe, non-invasive option for promoting recovery, reducing inflammation, and improving circulation.

The Beamer-Tech PEMF System is designed to work seamlessly with the RDPV4, allowing users to benefit from the full range of frequency sets offered by ETDFL. This integration supports targeted wellness goals, including relaxation, muscle recovery, pain relief, and overall cellular balance.

The ETDFL list is updated annually, and the latest searchable PDF version is available for free at: **ETDFL.com**

For up-to-date frequency programs and broader wellness applications not listed in printed materials, users are encouraged to explore the digital database at ETDFL.com.

# Pulsed Electro Frequency Magnetic Field Mat



Beamer PEFM

**BEAMER**   
BEAMER.TECH

The Beamer-Tech Store focuses on advanced wellness solutions by leveraging cutting-edge PEMF technology to support recovery, enhance circulation, and promote overall well-being.

## Contact Us

 [www.beamer.tech](http://www.beamer.tech)

 206-708-6408

 [etdfl.web@gmail.com](mailto:etdfl.web@gmail.com)

### **Beamer Tech Solutions**

1712 Pioneer Avenue, Cheyenne, WY  
82001, USA E Parcel & Mail  
Distribution office, no Pickups